



---

# Brain Imaging Data Structure Specification

---

v1.2.2

2020-02-12

# Contents

|          |  |           |
|----------|--|-----------|
| <b>1</b> | <b>The Brain Imaging Data Structure</b>                            | <b>4</b>  |
| <b>2</b> | <b>Introduction</b>  | <b>5</b>  |
|          | Motivation . . . . .   | 5         |
|          | Extensions . . . . .   | 5         |
|          | Citing BIDS . . . . .  | 5         |
| <b>3</b> | <b>Common principles</b>   | <b>7</b>  |
|          | Definitions . . . . .  | 7         |
|          | Compulsory, optional, and additional data and metadata . . . . .   | 7         |
|          | File name structure . . . . .                                      | 8         |
|          | Source vs. raw vs. derived data . . . . .                          | 8         |
|          | The Inheritance Principle . . . . .                                | 9         |
|          | Good practice recommendations . . . . .                            | 10        |
|          | File Formation specification . . . . .                             | 10        |
|          | Imaging files . . . . .  | 10        |
|          | Tabular files . . . . .  | 10        |
|          | Key/value files (dictionaries) . . . . .                           | 11        |
|          | Participant names and other labels . . . . .                       | 11        |
|          | Units . . . . .  | 11        |
|          | Directory structure . . . . .                                      | 12        |
|          | Single session example . . . . .                                   | 12        |
|          | Unspecified data . . . . .   | 13        |
| <b>4</b> | <b>Modality-agnostic files</b>                                     | <b>14</b> |
|          | Dataset description . . . . .                                      | 14        |
|          | dataset_description.json . . . . .                                 | 14        |
|          | README . . . . .   | 15        |
|          | CHANGES . . . . .  | 15        |
|          | Participants file . . . . .  | 15        |
|          | Phenotypic and assessment data . . . . .                           | 15        |
|          | Scans file . . . . .   | 16        |
|          | Code . . . . .   | 17        |
| <b>5</b> | <b>Magnetic Resonance Imaging data</b>                             | <b>18</b> |
|          | Common metadata fields . . . . .                                   | 18        |
|          | Scanner Hardware . . . . .   | 18        |
|          | Anatomy imaging data . . . . .                                     | 22        |
|          | Task (including resting state) imaging data . . . . .              | 23        |
|          | Diffusion imaging data . . . . .                                   | 27        |
|          | Fieldmap data . . . . .  | 28        |
| <b>6</b> | <b>Magnetoencephalography (MEG)</b>                                | <b>30</b> |
|          | MEG recording data . . . . .                                       | 30        |
|          | Sidecar JSON (*_meg.json) . . . . .                                | 30        |
|          | Channels description (*_channels.tsv) . . . . .                    | 33        |
|          | Coordinate System JSON (*_coordsystem.json) . . . . .              | 35        |
|          | Landmark photos (*_photo.jpg) . . . . .                            | 37        |
|          | Head shape and electrode description (*_headshape.<ext>) . . . . . | 38        |

|  |           |
|--|-----------|
| Empty-room MEG recordings . . . . .                              | 38        |
| <b>7 Electroencephalography (EEG)</b>                            | <b>40</b> |
| EEG recording data . . . . .                                     | 40        |
| Sidecar JSON (*_eeg.json) . . . . .                              | 41        |
| Channels description (*_channels.tsv) . . . . .                  | 44        |
| Electrodes description (*_electrodes.tsv) . . . . .              | 45        |
| Coordinate System JSON (*_coordsystem.json) . . . . .            | 46        |
| Landmark photos (*_photo.jpg) . . . . .                          | 48        |
| <b>8 intracranial Electroencephalography (iEEG)</b>              | <b>49</b> |
| iEEG recording data . . . . .                                    | 49        |
| Terminology: Electrodes vs. Channels . . . . .                   | 50        |
| Sidecar JSON (*_ieeg.json) . . . . .                             | 50        |
| Channels description (*_channels.tsv) . . . . .                  | 53        |
| Electrode description (*_electrodes.tsv) . . . . .               | 55        |
| Coordinate System JSON (*_coordsystem.json) . . . . .            | 56        |
| Recommended 3D coordinate systems . . . . .                      | 57        |
| Allowed 2D coordinate systems . . . . .                          | 57        |
| Multiple coordinate systems . . . . .                            | 57        |
| Photos of the electrode positions (*_photo.jpg) . . . . .        | 57        |
| Electrical stimulation . . . . .                                 | 59        |
| <b>9 Task events</b>   | <b>60</b> |
| <b>10 Physiological and other continuous recordings</b>          | <b>62</b> |
| <b>11 Behavioral experiments (with no MRI)</b>                   | <b>64</b> |
| <b>12 Longitudinal studies with multiple sessions (visits)</b>   | <b>65</b> |
| Sessions file . . . . .  | 66        |
| Multi-site or multi-center studies . . . . .                     | 66        |
| Treat each site/center as a separate dataset . . . . .           | 66        |
| Option 2: Combining sites/centers into one dataset . . . . .     | 66        |
| <b>13 BIDS Extension Proposals</b>                               | <b>67</b> |
| <b>14 Appendix I: Contributors</b>                               | <b>68</b> |
| <b>15 Appendix II: Licenses</b>                                  | <b>71</b> |
| <b>16 Appendix III: Hierarchical Event Descriptor (HED) Tags</b> | <b>72</b> |
| <b>17 Appendix IV: Entity table</b>                              | <b>74</b> |
| <b>18 Appendix V: Units</b>                                      | <b>75</b> |
| Prefixes . . . . .   | 75        |
| Multiples . . . . .  | 75        |
| Submultiples . . . . .   | 76        |
| <b>19 Appendix VI: MEG file formats</b>                          | <b>77</b> |
| CTF . . . . .  | 77        |
| Neuromag/Elekta/MEGIN . . . . .                                  | 78        |
| BTi/4D neuroimaging . . . . .                                    | 78        |
| KIT/Yokogawa/Ricoh . . . . .                                     | 79        |
| KRISS . . . . .  | 79        |
| ITAB . . . . .   | 80        |
| Aalto MEG-MRI . . . . .  | 80        |
| <b>20 Appendix VII: preferred names of MEG systems</b>           | <b>81</b> |
| <b>21 Appendix VIII: Preferred Names of Coordinate Systems</b>   | <b>82</b> |
| Coordinate Systems applicable to MEG, EEG, and iEEG . . . . .    | 82        |

|   |           |
|---|-----------|
| MEG Specific Coordinate Systems . . . . .   | 82        |
| EEG Specific Coordinate Systems . . . . .   | 83        |
| iEEG Specific Coordinate Systems . . . . .  | 83        |
| Template Based Coordinate Systems . . . . . | 83        |
| Standard template identifiers . . . . .     | 83        |
| <b>22 Changelog</b>                         | <b>85</b> |
| v1.2.2 (2020-02-12) . . . . .               | 85        |
| v1.2.1 (2019-08-14) . . . . .               | 85        |
| v1.2.0 (2019-03-04) . . . . .               | 86        |
| v1.1.2 (2019-01-10) . . . . .               | 87        |
| 1.1.1 . . . . .                             | 88        |
| 1.1.0 . . . . .                             | 88        |
| 1.0.2 . . . . .                             | 88        |
| 1.0.1 . . . . .                             | 88        |
| 1.0.1-rc1 . . . . .                         | 88        |
| 1.0.0 . . . . .                             | 89        |
| 1.0.0-rc4 . . . . .                         | 89        |
| 1.0.0-rc3 . . . . .                         | 89        |
| 1.0.0-rc2 . . . . .                         | 89        |
| 1.0.0-rc1 . . . . .                         | 90        |
| 1.1.1 . . . . .                             | 90        |
| 1.1.0 . . . . .                             | 90        |
| 1.0.2 . . . . .                             | 90        |
| 1.0.1 . . . . .                             | 90        |
| 1.0.1-rc1 . . . . .                         | 90        |
| 1.0.0 . . . . .                             | 91        |
| 1.0.0-rc4 . . . . .                         | 91        |
| 1.0.0-rc3 . . . . .                         | 91        |
| 1.0.0-rc2 . . . . .                         | 91        |
| 1.0.0-rc1 . . . . .                         | 92        |

# Chapter 1

## The Brain Imaging Data Structure

This site serves as an online resource to see the current state of the Brain Imaging Data Structure (BIDS) specification. It contains information about the core specification, as well as many modality-specific extensions.

To get started, check out the introduction. If you'd like more information on how to adapt your own datasets to match the BIDS specification, we recommend exploring the [bids-specification starter kit](#).

For an overview of the BIDS ecosystem, visit the [BIDS homepage](#).

# Chapter 2

## Introduction

### Motivation

Neuroimaging experiments result in complicated data that can be arranged in many different ways. So far there is no consensus how to organize and share data obtained in neuroimaging experiments. Even two researchers working in the same lab can opt to arrange their data in a different way. Lack of consensus (or a standard) leads to misunderstandings and time wasted on rearranging data or rewriting scripts expecting certain structure. Here we describe a simple and easy-to-adopt way of organising neuroimaging and behavioral data. By using this standard you will benefit in the following ways:

- It will be easy for another researcher to work on your data. To understand the organisation of the files and their format you will only need to refer them to this document. This is especially important if you are running your own lab and anticipate more than one person working on the same data over time. By using BIDS you will save time trying to understand and reuse data acquired by a graduate student or postdoc that has already left the lab.
- There are a growing number of data analysis software packages that can understand data organised according to BIDS (see <http://bids.neuroimaging.io> for the most up to date list).
- Databases such as OpenNeuro.org accept datasets organised according to BIDS. If you ever plan to share your data publicly (nowadays some journals require this) you can minimize the additional time and energy spent on publication, and speed up the curation process by using BIDS to structure and describe your data right after acquisition.
- There are [validation tools](#) that can check your dataset integrity and let you easily spot missing values.

BIDS is heavily inspired by the format used internally by OpenfMRI.org and has been supported by the International Neuroinformatics Coordinating Facility and the Neuroimaging Data Sharing Task Force. While working on BIDS we consulted many neuroscientists to make sure it covers most common experiments, but at the same time is intuitive and easy to adopt. The specification is intentionally based on simple file formats and folder structures to reflect current lab practices and make it accessible to a wide range of scientists coming from different backgrounds.

### Extensions

The BIDS specification can be extended in a backwards compatible way and will evolve over time. This is accomplished through community-driven BIDS Extension Proposals (BEPs). For more information about the BEP process, see [Extending the BIDS specification](#).

### Citing BIDS

When referring to BIDS in context of academic literature please cite:

Gorgolewski, K.J., Auer, T., Calhoun, V.D., Craddock, R.C., Das, S., Duff, E.P., Flandin, G., Ghosh, S.S., Glatard, T., Halchenko, Y.O., Handwerker, D.A., Hanke, M., Keator, D., Li, X., Michael, Z., Maumet, C., Nichols, B.N., Nichols, T.E., Pellman, J., Poline, J.-B., Rokem, A., Schaefer, G., Sochat, V., Triplett, W., Turner, J.A., Varoquaux, G., Poldrack, R.A., 2016. [The brain imaging](#)

[data structure](#), a format for organizing and describing outputs of neuroimaging experiments. Sci Data 3, 160044.

as well as other papers describing specific BIDS extensions (see below).

BIDS has also a [Research Resource Identifier \(RRID\)](#) - `RRID:SCR_016124` - which you can also include in your manuscript in addition to citing the paper.

# Chapter 3

## Common principles

### Definitions

The keywords "MUST", "MUST NOT", "REQUIRED", "SHALL", "SHALL NOT", "SHOULD", "SHOULD NOT", "RECOMMENDED", "MAY", and "OPTIONAL" in this document are to be interpreted as described in [RFC2119].

Throughout this specification we use a list of terms. To avoid misunderstanding we clarify them here.

1. Dataset - a set of neuroimaging and behavioral data acquired for a purpose of a particular study. A dataset consists of data acquired from one or more subjects, possibly from multiple sessions.
2. Subject - a person or animal participating in the study.
3. Session - a logical grouping of neuroimaging and behavioral data consistent across subjects. Session can (but doesn't have to) be synonymous to a visit in a longitudinal study. In general, subjects will stay in the scanner during one session. However, for example, if a subject has to leave the scanner room and then be re-positioned on the scanner bed, the set of MRI acquisitions will still be considered as a session and match sessions acquired in other subjects. Similarly, in situations where different data types are obtained over several visits (for example fMRI on one day followed by DWI the day after) those can be grouped in one session. Defining multiple sessions is appropriate when several identical or similar data acquisitions are planned and performed on all -or most- subjects, often in the case of some intervention between sessions (e.g., training).
4. Data acquisition - a continuous uninterrupted block of time during which a brain scanning instrument was acquiring data according to particular scanning sequence/protocol.
5. Data type - a functional group of different types of data. In BIDS we define eight data types: func (task based and resting state functional MRI), dwi (diffusion weighted imaging), fmap (field inhomogeneity mapping data such as field maps), anat (structural imaging such as T1, T2, etc.), meg (magnetoencephalography), eeg (electroencephalography), ieeg (intracranial electroencephalography), beh (behavioral).
6. Task - a set of structured activities performed by the participant. Tasks are usually accompanied by stimuli and responses, and can greatly vary in complexity. For the purpose of this specification we consider the so-called "resting state" a task. In the context of brain scanning, a task is always tied to one data acquisition. Therefore, even if during one acquisition the subject performed multiple conceptually different behaviors (with different sets of instructions) they will be considered one (combined) task.
7. Event - a stimulus or subject response recorded during a task. Each event has an onset time and duration. Note that not all tasks will have recorded events (e.g., resting state).
8. Run - an uninterrupted repetition of data acquisition that has the same acquisition parameters and task (however events can change from run to run due to different subject response or randomized nature of the stimuli). Run is a synonym of a data acquisition.

### Compulsory, optional, and additional data and metadata

The following standard describes a way of arranging data and writing down metadata for a subset of neuroimaging experiments. Some aspects of the standard are compulsory. For example a particular file name format is



required when storing structural scans. Some aspects are regulated but optional. For example a T2 volume does not need to be included, but when it is available it should be saved under a particular file name specified in the standard. This standard aspires to describe a majority of datasets, but acknowledges that there will be cases that do not fit. In such cases one can include additional files and subfolders to the existing folder structure following common sense. For example one may want to include eye tracking data in a vendor specific format that is not covered by this standard. The most sensible place to put it is next to the continuous recording file with the same naming scheme but different extensions. The solutions will change from case to case and publicly available datasets will be reviewed to include common data types in the future releases of the BIDS specification.

## File name structure

A file name consists of a chain of *entities*, or key-value pairs, a *suffix* and an *extension*. Two prominent examples of entities are **subject** and **session**.

For a data file that was collected in a given **session** from a given **subject**, the file name will begin with the string `sub-<label>_ses-<label>`.

Note that `sub-<label>` corresponds to the **subject** entity because it has the `sub-` "key" and `<label>` "value", where `<label>` would in a real data file correspond to a unique identifier of that subject, such as 01. The same holds for the **session** entity with its `ses-` key and its `<label>` value.

A chain of entities, followed by a suffix, connected by underscores (`_`) produces a human readable file name, such as `sub-01_task-rest_eeg.edf`. It is evident from the file name alone that the file contains resting state data from subject 01. The suffix `eeg` and the extension `.edf` depend on the imaging modality and the data format and indicate further details of the file's contents.

A summary of all entities in BIDS and the order in which they **MUST** be specified is available in the entity table in the appendix.

## Source vs. raw vs. derived data

BIDS in its current form is designed to harmonize and describe raw (unprocessed or minimally processed due to file format conversion) data. During analysis such data will be transformed and partial as well as final results will be saved. Derivatives of the raw data (other than products of DICOM to NIfTI conversion) **MUST** be kept separate from the raw data. This way one can protect the raw data from accidental changes by file permissions. In addition it is easy to distinguish partial results from the raw data and share the latter. Similar rules apply to source data which is defined as data before harmonization and/or file format conversion (for example E-Prime event logs or DICOM files).

This specification currently does not go into details of recommending a particular naming scheme for including different types of source data (raw event logs, parameter files, etc. before conversion to BIDS) and data derivatives (correlation maps, brain masks, contrasts maps, etc.). However, in the case that these data are to be included:

1. These data **MUST** be kept in separate **sourcedata** and **derivatives** folders each with a similar folder structure as presented below for the BIDS-managed data. For example: `derivatives/fmriprep/sub-01/ses-pre/sub-01_ses-pre_mask.nii.gz` or `sourcedata/sub-01/ses-pre/func/sub-01_ses-pre_task-rest_bold.dicom.tgz` or `sourcedata/sub-01/ses-pre/func/MyEvent.sce`.
2. A README file **SHOULD** be found at the root of the **sourcedata** or the **derivatives** folder (or both). This file should describe the nature of the raw data or the derived data. In the case of the existence of a **derivatives** folder, we **RECOMMEND** including details about the software stack and settings used to generate the results. Inclusion of non-imaging objects that improve reproducibility are encouraged (scripts, settings files, etc.).
3. We **RECOMMEND** including the PDF print-out with the actual sequence parameters generated by the scanner in the **sourcedata** folder.

Alternatively one can organize their data in the following way

```
1 my_dataset/  
2   sourcedata/  
3   ...  
4   rawdata/
```

```

5   dataset_description.json
6   participants.tsv
7   sub-01/
8   sub-02/
9   ...
10  derivatives/
11  ...

```

In this example **only the rawdata subfolder needs to be a BIDS compliant dataset**. This specification does not prescribe anything about the contents of `sourcedata` and `derivatives` folders in the above example - nor does it prescribe the `sourcedata`, `derivatives`, or `rawdata` folder names. The above example is just a convention that can be useful for organizing raw, source, and derived data while maintaining BIDS compliancy of the raw data folder.

## The Inheritance Principle

Any metadata file (`.json`, `.bvec`, `.tsv`, etc.) may be defined at any directory level, but no more than one applicable file may be defined at a given level (Example 1). The values from the top level are inherited by all lower levels unless they are overridden by a file at the lower level. For example, `sub-*_task-rest_bold.json` may be specified at the participant level, setting TR to a specific value. If one of the runs has a different TR than the one specified in that file, another `sub-*_task-rest_bold.json` file can be placed within that specific series directory specifying the TR for that specific run. There is no notion of "unsetting" a key/value pair. For example if there is a JSON file corresponding to particular participant/run defining a key/value and there is a JSON file on the root level of the dataset that does not define this key/value it will not be "unset" for all subjects/runs. Files for a particular participant can exist only at participant level directory, i.e. `/dataset/sub-*/[ses-*/sub-*_T1w.json`. Similarly, any file that is not specific to a participant is to be declared only at top level of dataset for eg: `task-sist_bold.json` must be placed under `/dataset/task-sist_bold.json`

Example 1: Two JSON files that are erroneously at the same level

```

1 sub-01/
2   ses-test/
3     sub-01_ses-test_task-overtverbgeneration_bold.json
4     sub-01_ses-test_task-overtverbgeneration_run-2_bold.json
5     anat/
6       sub-01_ses-test_T1w.nii.gz
7     func/
8       sub-01_ses-test_task-overtverbgeneration_run-1_bold.nii.gz
9       sub-01_ses-test_task-overtverbgeneration_run-2_bold.nii.gz

```

In the above example, two JSON files are listed under `sub-01/ses-test/`, which are each applicable to `sub-01_ses-test_task-overtverbgeneration_run-2_bold.nii.gz`, violating the constraint that no more than one file may be defined at a given level of the directory structure. Instead `sub-01_ses-test_task-overtverbgeneration_run-2_bold.json` should have been under `sub-01/ses-test/func/`.

Example 2: Multiple runs and recs with same acquisition (acq) parameters

```

1 sub-01/
2   anat/
3   func/
4     sub-01_task-xyz_acq-test1_run-1_bold.nii.gz
5     sub-01_task-xyz_acq-test1_run-2_bold.nii.gz
6     sub-01_task-xyz_acq-test1_rec-recon1_bold.nii.gz
7     sub-01_task-xyz_acq-test1_rec-recon2_bold.nii.gz
8     sub-01_task-xyz_acq-test1_bold.json

```

For the above example, all NIfTI files are acquired with same scanning parameters (`acq-test1`). Hence a JSON file describing the acq parameters will apply to different runs and rec files. Also if the JSON file (`task-xyz_acq-test1_bold.json`) is defined at dataset top level directory, it will be applicable to all task runs with `test1` acquisition parameter.

Example 3: Multiple JSON files at different levels for same task and acquisition parameters

```

1 task-xyz_acq-test1_bold.json
2 sub-01/
3   anat/
4   func/
5       sub-01_task-xyz_acq-test1_run-1_bold.nii.gz
6       sub-01_task-xyz_acq-test1_rec-recon1_bold.nii.gz
7       sub-01_task-xyz_acq-test1_rec-recon2_bold.nii.gz
8       sub-01_task-xyz_acq-test1_bold.json

```

In the above example, the fields from the `task-xyz_acq-test1_bold.json` file at the top directory will apply to all bold runs. However, if there is a key with different value in the `sub-01/func/sub-01_task-xyz_acq-test1_bold.json` file defined at a deeper level, that value will be applicable for that particular run/task NIfTI file/s. In other words, the `.json` file at the deeper level overrides values that are potentially also defined in the `.json` at a more shallow level. If the `.json` file at the more shallow level contains key-value-pairs that are not present in the `.json` file at the deeper level, these key-value-pairs are inherited by the `.json` file at the deeper level (but NOT vice versa!).

## Good practice recommendations

**Try to avoid excessive amount of overrides.** Do not specify a field value in the upper levels if lower levels have more or less even distribution of multiple possible values. E.g., if a field `X` has one value for all `ses-01/` and another for all `ses-02/` it better not to be defined at all in the `.json` at the upper level.

## File Formation specification

### Imaging files

All imaging data **MUST** be stored using the NIfTI file format. We **RECOMMEND** using compressed NIfTI files (`.nii.gz`), either version 1.0 or 2.0. Imaging data **SHOULD** be converted to the NIfTI format using a tool that provides as much of the NIfTI header information (such as orientation and slice timing information) as possible. Since the NIfTI standard offers limited support for the various image acquisition parameters available in DICOM files, we **RECOMMEND** that users provide additional meta information extracted from DICOM files in a sidecar JSON file (with the same filename as the `.nii[.gz]` file, but with a `.json` extension). Extraction of BIDS compatible metadata can be performed using `dcm2niix` and `dicm2nii` DICOM to NIfTI converters. A provided `validator` will check for conflicts between the JSON file and the data recorded in the NIfTI header.

### Tabular files

Tabular data **MUST** be saved as tab delimited values (`.tsv`) files, i.e., CSV files where commas are replaced by tabs. Tabs **MUST** be true tab characters and **MUST NOT** be a series of space characters. Each TSV file **MUST** start with a header line listing the names of all columns (with the exception of physiological and other continuous acquisition data - see below for details). Names **MUST** be separated with tabs. String values containing tabs **MUST** be escaped using double quotes. Missing and non-applicable values **MUST** be coded as `n/a`. Numerical values **MUST** employ the dot (`.`) as decimal separator and **MAY** be specified in scientific notation, using `e` or `E` to separate the significand from the exponent. TSV files **MUST** be in UTF-8 encoding.

Example:

```

1 onset duration response_time correct stop_trial go_trial
2 200 200 0 n/a n/a n/a

```

Tabular files **MAY** be optionally accompanied by a simple data dictionary in a JSON format (see below). The data dictionaries **MUST** have the same name as their corresponding tabular files but with `.json` extensions. If a JSON file is provided, it **MAY** contain one or more fields describing the columns found in the TSV file (in addition to any other metadata one wishes to include that describe the file as a whole). Note that if a field name included in the JSON sidecar matches a column name in the TSV file, then that field **MUST** contain a description of the corresponding column, using an object containing the following fields:

| Field    |  |
|----------|--|
| name     | Definition                               |
| LongName | Long (unabbreviated) name of the column. |

| Field       |  |
|-------------|--|
| name        | Definition   |
| Description | Description of the column.   |
| Levels      | For categorical variables: a dictionary of possible values (keys) and their descriptions (values).   |
| Units       | Measurement units. [ <b>&lt;prefix symbol&gt;</b> ] <b>&lt;unit symbol&gt;</b> format following the SI standard is RECOMMENDED (see Appendix V). |
| TermURL     | URL pointing to a formal definition of this type of data in an ontology available on the web.  |

Example:

```

1 {
2   "test": {
3     "LongName": "Education level",
4     "Description": "Education level, self-rated by participant",
5     "Levels": {
6       "1": "Finished primary school",
7       "2": "Finished secondary school",
8       "3": "Student at university",
9       "4": "Has degree from university"
10    }
11  },
12  "bmi": {
13    "LongName": "Body mass index",
14    "Units": "kilograms per squared meters",
15    "TermURL": "http://purl.bioontology.org/ontology/SNOMEDCT/60621009"
16  }
17 }
```

## Key/value files (dictionaries)

JavaScript Object Notation (JSON) files MUST be used for storing key/value pairs. JSON files MUST be in UTF-8 encoding. Extensive documentation of the format can be found here: <http://json.org/>. Several editors have built-in support for JSON syntax highlighting that aids manual creation of such files. An online editor for JSON with built-in validation is available at: <http://jsoneditoronline.org>.

Example:

```

1 {
2   "RepetitionTime": 3,
3   "Instruction": "Lie still and keep your eyes open"
4 }
```

## Participant names and other labels

BIDS uses custom user-defined labels in several situations (naming of participants, sessions, acquisition schemes, etc.) Labels are strings and MUST only consist of letters (lower or upper case) and/or numbers. If numbers are used we RECOMMEND zero padding (e.g., 01 instead of 1 if you have more than nine subjects) to make alphabetical sorting more intuitive.

Please note that a given label is distinct from the "prefix" it refers to. For example **sub-01** refers to the **sub** entity (a subject) with the label 01. The **sub-** prefix is not part of the subject label, but must be included in file names (similarly to other key names). In contrast to other labels, **run** and **echo** labels MUST be integers. Those labels MAY include zero padding, but this is NOT RECOMMENDED to maintain their uniqueness.

## Units

All units SHOULD be specified as per International System of Units (abbreviated as SI, from the French *Système international* (d'unités)) and can be SI units or SI derived units. In case there are valid reasons to deviate from SI units or SI derived units, the units MUST be specified in the sidecar JSON file. In case data is expressed in

SI units or SI derived units, the units MAY be specified in the sidecar JSON file. In case prefixes are added to SI or non-SI units (e.g., mm), the prefixed units MUST be specified in the JSON file (see Appendix V: Units). In particular:

- Elapsed time SHOULD be expressed in seconds. Please note that some DICOM parameters have been traditionally expressed in milliseconds. Those need to be converted to seconds.
- Frequency SHOULD be expressed in Hertz.

Describing dates and timestamps:

- Date time information MUST be expressed in the following format `YYYY-MM-DDThh:mm:ss` (one of the [ISO8601](#) date-time formats). For example: `2009-06-15T13:45:30`
- Time stamp information MUST be expressed in the following format: `13:45:30`
- Dates can be shifted by a random number of days for privacy protection reasons. To distinguish real dates from shifted dates always use year 1925 or earlier when including shifted years. For longitudinal studies please remember to shift dates within one subject by the same number of days to maintain the interval information. Example: `1867-06-15T13:45:30`
- Age SHOULD be given as the number of years since birth at the time of scanning (or first scan in case of multi session datasets). Using higher accuracy (weeks) should in general be avoided due to privacy protection, unless when appropriate given the study goals, e.g., when scanning babies.

## Directory structure

### Single session example

This is an example of the folder and file structure. Because there is only one session, the session level is not required by the format. For details on individual files see descriptions in the next section:

```

1 sub-control01/
2   anat/
3     sub-control01_T1w.nii.gz
4     sub-control01_T1w.json
5     sub-control01_T2w.nii.gz
6     sub-control01_T2w.json
7   func/
8     sub-control01_task-nback_bold.nii.gz
9     sub-control01_task-nback_bold.json
10    sub-control01_task-nback_events.tsv
11    sub-control01_task-nback_physio.tsv.gz
12    sub-control01_task-nback_physio.json
13    sub-control01_task-nback_sbref.nii.gz
14   dwi/
15     sub-control01_dwi.nii.gz
16     sub-control01_dwi.bval
17     sub-control01_dwi.bvec
18   fmap/
19     sub-control01_phasediff.nii.gz
20     sub-control01_phasediff.json
21     sub-control01_magnitude1.nii.gz
22     sub-control01_scans.tsv
23 code/
24   deface.py
25 derivatives/
26 README
27 participants.tsv
28 dataset_description.json
29 CHANGES

```

## Unspecified data

Additional files and folders containing raw data MAY be added as needed for special cases. All non-standard file entities SHOULD conform to BIDS-style naming conventions, including alphabetic entities and suffixes and alphanumeric labels/indices. Non-standard suffixes SHOULD reflect the nature of the data, and existing entities SHOULD be used when appropriate. For example, an ASSET calibration scan might be named `sub-01_acq-ASSET_calibration.nii.gz`.

Non-standard files and directories should be named with care. Future BIDS efforts may standardize new entities and suffixes, changing the meaning of file names and setting requirements on their contents or metadata. Validation and parsing tools MAY treat the presence of non-standard files and directories as an error, so consult the details of these tools for mechanisms to suppress warnings or provide interpretations of your file names.

## Chapter 4

# Modality-agnostic files

### Dataset description

Template: `dataset_description.json` README CHANGES

#### `dataset_description.json`

The file `dataset_description.json` is a JSON file describing the dataset. Every dataset MUST include this file with the following fields:

| Field name         | Definition   |
|--------------------|--|
| Name               | REQUIRED. Name of the dataset.   |
| BIDSVersion        | REQUIRED. The version of the BIDS standard that was used.  |
| License            | RECOMMENDED. What license is this dataset distributed under? The use of license name abbreviations is suggested for specifying a license. A list of common licenses with suggested abbreviations can be found in Appendix II.        |
| Authors            | OPTIONAL. List of individuals who contributed to the creation/curation of the dataset.   |
| Acknowledgements   | OPTIONAL. Text acknowledging contributions of individuals or institutions beyond those listed in Authors or Funding.   |
| HowToAcknowledge   | OPTIONAL. Text containing instructions on how researchers using this dataset should acknowledge the original authors. This field can also be used to define a publication that should be cited in publications that use the dataset. |
| Funding            | OPTIONAL. List of sources of funding (grant numbers)   |
| ReferencesAndLinks | OPTIONAL. List of references to publication that contain information on the dataset, or links.   |
| DatasetDOI         | OPTIONAL. The Document Object Identifier of the dataset (not the corresponding paper).   |

Example:

```
1 {
2   "Name": "The mother of all experiments",
3   "BIDSVersion": "1.0.1",
4   "License": "CC0",
5   "Authors": [
6     "Paul Broca",
7     "Carl Wernicke"
8   ],
9   "Acknowledgements": "Special thanks to Korbinian Brodmann for help in formatting this dataset
    in BIDS. We thank Alan Lloyd Hodgkin and Andrew Huxley for helpful comments and
    discussions about the experiment and manuscript; Hermann Ludwig Helmholtz for
    administrative support; and Claudius Galenus for providing data for the medial-to-lateral
    index analysis."
```

```

10  "HowToAcknowledge": "Please cite this paper: https://www.ncbi.nlm.nih.gov/pubmed
    /001012092119281",
11  "Funding": [
12    "National Institute of Neuroscience Grant F378236MFH1",
13    "National Institute of Neuroscience Grant 5RMZ0023106"
14  ],
15  "ReferencesAndLinks": [
16    "https://www.ncbi.nlm.nih.gov/pubmed/001012092119281",
17    "Alzheimer A., & Kraepelin, E. (2015). Neural correlates of presenile dementia in humans.
    Journal of Neuroscientific Data, 2, 234001. http://doi.org/1920.8/jndata.2015.7"
18  ],
19  "DatasetDOI": "10.0.2.3/dfjj.10"
20 }

```

## README

In addition a free form text file (README) describing the dataset in more details SHOULD be provided. The README file MUST be either in ASCII or UTF-8 encoding.

## CHANGES

Version history of the dataset (describing changes, updates and corrections) MAY be provided in the form of a CHANGES text file. This file MUST follow the CPAN Changelog convention: <http://search.cpan.org/~haarg/CPAN-Changes-0.400002/lib/CPAN/Changes/Spec.pod>. The CHANGES file MUST be either in ASCII or UTF-8 encoding.

Example:

```

1 1.0.1 2015-08-27
2 - Fixed slice timing information.
3
4 1.0.0 2015-08-17
5 - Initial release.

```

## Participants file

Template:

```

1 participants.tsv
2 participants.json
3 phenotype/<measurement_tool_name>.tsv
4 phenotype/<measurement_tool_name>.json

```

Optional: Yes

The purpose of this file is to describe properties of participants such as age, handedness, sex, etc. In case of single session studies this file has one compulsory column **participant\_id** that consists of **sub-<label>**, followed by a list of optional columns describing participants. Each participant needs to be described by one and only one row.

participants.tsv example:

```

1 participant_id age sex group
2 sub-control01 34 M control
3 sub-control02 12 F control
4 sub-patient01 33 F patient

```

## Phenotypic and assessment data

If the dataset includes multiple sets of participant level measurements (for example responses from multiple questionnaires) they can be split into individual files separate from **participants.tsv**.



Each of the measurement files **MUST** be kept in a `/phenotype` directory placed at the root of the BIDS dataset and **MUST** end with the `.tsv` extension. File names **SHOULD** be chosen to reflect the contents of the file. For example, the "Adult ADHD Clinical Diagnostic Scale" could be saved in a file called `/phenotype/acds_adult.tsv`.

The files can include an arbitrary set of columns, but one of them **MUST** be `participant_id` and the entries of that column **MUST** correspond to the subjects in the BIDS dataset and `participants.tsv` file.

As with all other tabular data, the additional phenotypic information files **MAY** be accompanied by a JSON file describing the columns in detail (see Tabular files). In addition to the column description, a section describing the measurement tool (as a whole) **MAY** be added under the name `MeasurementToolMetadata`. This section consists of two keys:

- **Description**: A free text description of the measurement tool
- **TermURL**: A link to an entity in an ontology corresponding to this tool.

As an example, consider the contents of a file called `phenotype/acds_adult.json`:

```
1 {
2   "MeasurementToolMetadata": {
3     "Description": "Adult ADHD Clinical Diagnostic Scale V1.2",
4     "TermURL": "http://www.cognitiveatlas.org/task/id/trm_5586ff878155d"
5   },
6   "adhd_b": {
7     "Description": "B. CHILDHOOD ONSET OF ADHD (PRIOR TO AGE 7)",
8     "Levels": {
9       "1": "YES",
10      "2": "NO"
11    }
12  },
13  "adhd_c_dx": {
14    "Description": "As child met A, B, C, D, E and F diagnostic criteria",
15    "Levels": {
16      "1": "YES",
17      "2": "NO"
18    }
19  }
20 }
```

Please note that in this example `MeasurementToolMetadata` includes information about the questionnaire and `adhd_b` and `adhd_c_dx` correspond to individual columns.

In addition to the keys available to describe columns in all tabular files (`LongName`, `Description`, `Levels`, `Units`, and `TermURL`) the `participants.json` file as well as phenotypic files can also include column descriptions with a `Derivative` field that, when set to true, indicates that values in the corresponding column is a transformation of values from other columns (for example a summary score based on a subset of items in a questionnaire).

## Scans file

Template:

```
1 sub-<label>/[ses-<label>/]
2   sub-<label>[_ses-<label>]_scans.tsv
```

Optional: Yes

The purpose of this file is to describe timing and other properties of each imaging acquisition sequence (each run `.nii.gz` file) within one session. Each `.nii.gz` file should be described by at most one row. Relative paths to files should be used under a compulsory `filename` header. If acquisition time is included it should be under `acq_time` header. Datetime should be expressed in the following format 2009-06-15T13:45:30 (year, month, day, hour (24h), minute, second; this is equivalent to the RFC3339 "date-time" format, time zone is always assumed as local time). For anonymization purposes all dates within one subject should be shifted by a randomly chosen (but common across all runs etc.) number of days. This way relative timing would be preserved, but chances of identifying a person based on the date and time of their scan would be decreased.

Dates that are shifted for anonymization purposes should be set to a year 1925 or earlier to clearly distinguish them from unmodified data. Shifting dates is RECOMMENDED, but not required.

Additional fields can include external behavioral measures relevant to the scan. For example vigilance questionnaire score administered after a resting state scan.

Example:

```
1 filename  acq_time
2 func/sub-control01_task-nback_bold.nii.gz 1877-06-15T13:45:30
3 func/sub-control01_task-motor_bold.nii.gz 1877-06-15T13:55:33
```

## Code

Template: `code/*`

Source code of scripts that were used to prepare the dataset (for example if it was anonymized or defaced) MAY be stored here.<sup>1</sup> Extra care should be taken to avoid including original IDs or any identifiable information with the source code. There are no limitations or recommendations on the language and/or code organization of these scripts at the moment.

<sup>1</sup>Storing actual source files with the data is preferred over links to external source repositories to maximize long term preservation (which would suffer if an external repository would not be available anymore).

## Chapter 5

# Magnetic Resonance Imaging data

### Common metadata fields

MR Data described in sections 8.3.x share the following RECOMMENDED metadata fields (stored in sidecar JSON files). MRI acquisition parameters are divided into several categories based on ["A checklist for fMRI acquisition methods reporting in the literature"](#) by Ben Inglis:

#### Scanner Hardware

| Field name                    | Definition  |
|-------------------------------|---|
| Manufacturer                  | RECOMMENDED. Manufacturer of the equipment that produced the composite instances. Corresponds to DICOM Tag 0008, 0070<br><b>Manufacturer</b>  |
| ManufacturersModelName        | RECOMMENDED. Manufacturer's model name of the equipment that produced the composite instances. Corresponds to DICOM Tag 0008, 1090<br><b>Manufacturers Model Name</b>   |
| DeviceSerialNumber            | RECOMMENDED. The serial number of the equipment that produced the composite instances. Corresponds to DICOM Tag 0018, 1000<br><b>DeviceSerialNumber</b> . A pseudonym can also be used to prevent the equipment from being identifiable, so long as each pseudonym is unique within the dataset |
| StationName                   | RECOMMENDED. Institution defined name of the machine that produced the composite instances. Corresponds to DICOM Tag 0008, 1010<br><b>Station Name</b>  |
| SoftwareVersions              | RECOMMENDED. Manufacturer's designation of software version of the equipment that produced the composite instances. Corresponds to DICOM Tag 0018, 1020 <b>Software Versions</b>  |
| HardcopyDeviceSoftwareVersion | (Deprecated) Manufacturer's designation of the software of the device that created this Hardcopy Image (the printer). Corresponds to DICOM Tag 0018, 101A <b>Hardcopy Device Software Version</b>   |
| MagneticFieldStrength         | RECOMMENDED. Nominal field strength of MR magnet in Tesla. Corresponds to DICOM Tag 0018,0087 <b>Magnetic Field Strength</b>  |
| ReceiveCoilName               | RECOMMENDED. Information describing the receiver coil. Corresponds to DICOM Tag 0018, 1250 <b>Receive Coil Name</b> , although not all vendors populate that DICOM Tag, in which case this field can be derived from an appropriate private DICOM field   |

| Field name                | Definition   |
|---------------------------|--|
| ReceiveCoilActiveElements | RECOMMENDED. Information describing the active/selected elements of the receiver coil. This doesn't correspond to a tag in the DICOM ontology. The vendor-defined terminology for active coil elements can go in this field. As an example, for Siemens, coil channels are typically not activated/selected individually, but rather in pre-defined selectable "groups" of individual channels, and the list of the "groups" of elements that are active/selected in any given scan populates the <b>Coil String</b> entry in Siemens' private DICOM fields (e.g., <b>HEA;HEP</b> for the Siemens standard 32 ch coil when both the anterior and posterior groups are activated). This is a flexible field that can be used as most appropriate for a given vendor and coil to define the "active" coil elements. Since individual scans can sometimes not have the intended coil elements selected, it is preferable for this field to be populated directly from the DICOM for each individual scan, so that it can be used as a mechanism for checking that a given scan was collected with the intended coil elements selected |
| GradientSetType           | RECOMMENDED. It should be possible to infer the gradient coil from the scanner model. If not, e.g. because of a custom upgrade or use of a gradient insert set, then the specifications of the actual gradient coil should be reported independently   |
| MRTransmitCoilSequence    | RECOMMENDED. This is a relevant field if a non-standard transmit coil is used. Corresponds to DICOM Tag 0018, 9049<br><b>MR Transmit Coil Sequence</b>   |
| MatrixCoilMode            | RECOMMENDED. (If used) A method for reducing the number of independent channels by combining in analog the signals from multiple coil elements. There are typically different default modes when using un-accelerated or accelerated (e.g. GRAPPA, SENSE) imaging  |
| CoilCombinationMethod     | RECOMMENDED. Almost all fMRI studies using phased-array coils use root-sum-of-squares (rSOS) combination, but other methods exist. The image reconstruction is changed by the coil combination method (as for the matrix coil mode above), so anything non-standard should be reported   |

### Sequence Specifics

| Field name                  | Definition   |
|-----------------------------|--|
| PulseSequenceType           | RECOMMENDED. A general description of the pulse sequence used for the scan (i.e. MPRAGE, Gradient Echo EPI, Spin Echo EPI, Multiband gradient echo EPI).   |
| ScanningSequence            | RECOMMENDED. Description of the type of data acquired. Corresponds to DICOM Tag 0018, 0020 <b>Scanning Sequence</b> .  |
| SequenceVariant             | RECOMMENDED. Variant of the ScanningSequence. Corresponds to DICOM Tag 0018, 0021 <b>Sequence Variant</b> .  |
| ScanOptions                 | RECOMMENDED. Parameters of ScanningSequence. Corresponds to DICOM Tag 0018, 0022 <b>Scan Options</b> .   |
| SequenceName                | RECOMMENDED. Manufacturer's designation of the sequence name. Corresponds to DICOM Tag 0018, 0024 <b>Sequence Name</b> .   |
| PulseSequenceDetails        | RECOMMENDED. Information beyond pulse sequence type that identifies the specific pulse sequence used (i.e. "Standard Siemens Sequence distributed with the VB17 software," "Siemens WIP ### version #.##," or "Sequence written by X using a version compiled on MM/DD/YYYY"). |
| NonlinearGradientCorrection | RECOMMENDED. Boolean stating if the image saved has been corrected for gradient nonlinearities by the scanner sequence.  |

### In-Plane Spatial Encoding

| Field name                     | Definition   |
|--------------------------------|--|
| NumberShots                    | RECOMMENDED. The number of RF excitations need to reconstruct a slice or volume. Please mind that this is not the same as Echo Train Length which denotes the number of lines of k-space collected after an excitation.  |
| ParallelReductionFactorInPlane | RECOMMENDED. The parallel imaging (e.g, GRAPPA) factor. Use the denominator of the fraction of k-space encoded for each slice. For example, 2 means half of k-space is encoded. Corresponds to DICOM Tag 0018, 9069<br><b>Parallel Reduction Factor In-plane.</b>  |
| ParallelAcquisitionTechnique   | RECOMMENDED. The type of parallel imaging used (e.g. GRAPPA, SENSE). Corresponds to DICOM Tag 0018, 9078<br><b>Parallel Acquisition Technique.</b>   |
| PartialFourier                 | RECOMMENDED. The fraction of partial Fourier information collected. Corresponds to DICOM Tag 0018, 9081<br><b>Partial Fourier.</b>   |
| PartialFourierDirection        | RECOMMENDED. The direction where only partial Fourier information was collected. Corresponds to DICOM Tag 0018, 9036<br><b>Partial Fourier Direction.</b>  |
| PhaseEncodingDirection         | RECOMMENDED. Possible values: i, j, k, i-, j-, k-. The letters i, j, k correspond to the first, second and third axis of the data in the NIFTI file. The polarity of the phase encoding is assumed to go from zero index to maximum index unless - sign is present (then the order is reversed - starting from the highest index instead of zero). <b>PhaseEncodingDirection</b> is defined as the direction along which phase is was modulated which may result in visible distortions. Note that this is not the same as the DICOM term <b>InPlanePhaseEncodingDirection</b> which can have ROW or COL values. This parameter is REQUIRED if corresponding fieldmap data is present or when using multiple runs with different phase encoding directions (which can be later used for field inhomogeneity correction). |
| EffectiveEchoSpacing           | RECOMMENDED. The "effective" sampling interval, specified in seconds, between lines in the phase-encoding direction, defined based on the size of the reconstructed image in the phase direction. It is frequently, but incorrectly, referred to as "dwell time" (see <b>DwellTime</b> parameter below for actual dwell time). It is required for unwarping distortions using field maps. Note that beyond just in-plane acceleration, a variety of other manipulations to the phase encoding need to be accounted for properly, including partial fourier, phase oversampling, phase resolution, phase field-of-view and interpolation. <sup>2</sup> This parameter is REQUIRED if corresponding fieldmap data is present.  |
| TotalReadoutTime               | RECOMMENDED. This is actually the "effective" total readout time, defined as the readout duration, specified in seconds, that would have generated data with the given level of distortion. It is NOT the actual, physical duration of the readout train. If <b>EffectiveEchoSpacing</b> has been properly computed, it is just <b>EffectiveEchoSpacing * (ReconMatrixPE - 1)</b> . <sup>3</sup> This parameter is REQUIRED if corresponding "field/distortion" maps acquired with opposing phase encoding directions are present (see 8.9.4).   |

<sup>2</sup>Conveniently, for Siemens' data, this value is easily obtained as  $1/[\text{BWPPPE} * \text{ReconMatrixPE}]$ , where BWPPPE is the **"BandwidthPerPixelPhaseEncode"** in DICOM tag (0019,1028) and ReconMatrixPE is the size of the actual reconstructed data in the phase direction (which is NOT reflected in a single DICOM tag for all possible aforementioned scan manipulations). See [here](#) and [here](#)

<sup>3</sup>We use the "FSL definition", i.e, the time between the center of the first "effective" echo and the center of the last "effective" echo.

**Timing Parameters**

| Field name             | Definition   |
|------------------------|--|
| EchoTime               | RECOMMENDED. The echo time (TE) for the acquisition, specified in seconds. This parameter is REQUIRED if corresponding fieldmap data is present or the data comes from a multi echo sequence. Corresponds to DICOM Tag 0018, 0081 <b>Echo Time</b> (please note that the DICOM term is in milliseconds not seconds).   |
| InversionTime          | RECOMMENDED. The inversion time (TI) for the acquisition, specified in seconds. Inversion time is the time after the middle of inverting RF pulse to middle of excitation pulse to detect the amount of longitudinal magnetization. Corresponds to DICOM Tag 0018, 0082 <b>Inversion Time</b> (please note that the DICOM term is in milliseconds not seconds).  |
| SliceTiming            | RECOMMENDED. The time at which each slice was acquired within each volume (frame) of the acquisition. Slice timing is not slice order -- rather, it is a list of times (in JSON format) containing the time (in seconds) of each slice acquisition in relation to the beginning of volume acquisition. The list goes through the slices along the slice axis in the slice encoding dimension (see below). Note that to ensure the proper interpretation of the <b>SliceTiming</b> field, it is important to check if the OPTIONAL <b>SliceEncodingDirection</b> exists. In particular, if <b>SliceEncodingDirection</b> is negative, the entries in <b>SliceTiming</b> are defined in reverse order with respect to the slice axis (i.e., the final entry in the <b>SliceTiming</b> list is the time of acquisition of slice 0). This parameter is REQUIRED for sparse sequences that do not have the <b>DelayTime</b> field set. In addition without this parameter slice time correction will not be possible. |
| SliceEncodingDirection | RECOMMENDED. Possible values: i, j, k, i-, j-, k- (the axis of the NIfTI data along which slices were acquired, and the direction in which <b>SliceTiming</b> is defined with respect to). i, j, k identifiers correspond to the first, second and third axis of the data in the NIfTI file. A - sign indicates that the contents of <b>SliceTiming</b> are defined in reverse order - that is, the first entry corresponds to the slice with the largest index, and the final entry corresponds to slice index zero. When present, the axis defined by <b>SliceEncodingDirection</b> needs to be consistent with the 'slice_dim' field in the NIfTI header. When absent, the entries in <b>SliceTiming</b> must be in the order of increasing slice index as defined by the NIfTI header.   |
| DwellTime              | RECOMMENDED. Actual dwell time (in seconds) of the receiver per point in the readout direction, including any oversampling. For Siemens, this corresponds to DICOM field (0019,1018) (in ns). This value is necessary for the optional readout distortion correction of anatomicals in the HCP Pipelines. It also usefully provides a handle on the readout bandwidth, which isn't captured in the other metadata tags. Not to be confused with <b>EffectiveEchoSpacing</b> , and the frequent mislabeling of echo spacing (which is spacing in the phase encoding direction) as "dwell time" (which is spacing in the readout direction).   |

**RF & Contrast**

| Field name                  | Definition   |
|-----------------------------|--|
| FlipAngle                   | RECOMMENDED. Flip angle for the acquisition, specified in degrees. Corresponds to: DICOM Tag 0018, 1314 <b>Flip Angle</b> .  |
| MultibandAccelerationFactor | RECOMMENDED. The multiband factor, for multiband acquisitions.   |
| NegativeContrast            | OPTIONAL. Boolean ( <b>true</b> or <b>false</b> ) value specifying whether increasing voxel intensity (within sample voxels) denotes a decreased value with respect to the contrast suffix. This is commonly the case when Cerebral Blood Volume is estimated via usage of a contrast agent in conjunction with a T2* weighted acquisition protocol. |

## Slice Acceleration

| Field name                  | Definition   |
|-----------------------------|--|
| MultibandAccelerationFactor | RECOMMENDED. The multiband factor, for multiband acquisitions. |

## Anatomical landmarks

Useful for multimodal co-registration with MEG, (S)EEG, TMS, etc.

| Field name                    | Definition   |
|-------------------------------|--|
| AnatomicalLandmarkCoordinates | RECOMMENDED. Key:value pairs of any number of additional anatomical landmarks and their coordinates in voxel units (where first voxel has index 0,0,0) relative to the associated anatomical MRI, (e.g. {"AC": [127,119,149], "PC": [128,93,141], "IH": [131,114,206]}, or {"NAS": [127,213,139], "LPA": [52,113,96], "RPA": [202,113,91]}). |

## Institution information

| Field name                  | Definition  |
|-----------------------------|---|
| InstitutionName             | RECOMMENDED. The name of the institution in charge of the equipment that produced the composite instances. Corresponds to DICOM Tag 0008, 0080 <b>InstitutionName</b> .                     |
| InstitutionAddress          | RECOMMENDED. The address of the institution in charge of the equipment that produced the composite instances. Corresponds to DICOM Tag 0008, 0081 <b>InstitutionAddress</b> .               |
| InstitutionalDepartmentName | RECOMMENDED. The department in the institution in charge of the equipment that produced the composite instances. Corresponds to DICOM Tag 0008, 1040 <b>Institutional Department Name</b> . |

When adding additional metadata please use the camelcase version of [DICOM ontology terms](#) whenever possible.

## Anatomy imaging data

Template:

```

1 sub-<label>/[ses-<label>/]
2   anat/
3     sub-<label>[_ses-<label>][_acq-<label>][_ce-<label>][_rec-<label>][_run-<index>]_<
4     modality_label>.nii.gz]
5     sub-<label>[_ses-<label>][_acq-<label>][_ce-<label>][_rec-<label>][_run-<index>][_mod-<
6     label>]_defacemask.nii.gz]
```

Anatomical (structural) data acquired for that participant. Currently supported modalities include:

| Name        | modality_label | Description  |
|-------------|----------------|--|
| T1 weighted | T1w            |  |
| T2 weighted | T2w            |  |
| T1 Rho map  | T1rho          | Quantitative T1rho brain imaging<br><a href="https://www.ncbi.nlm.nih.gov/pubmed/24474423">https://www.ncbi.nlm.nih.gov/pubmed/24474423</a><br><a href="https://www.ncbi.nlm.nih.gov/pmc/articles/PMC4346383/">https://www.ncbi.nlm.nih.gov/pmc/articles/PMC4346383/</a> |
| T1 map      | T1map          | quantitative T1 map  |
| T2 map      | T2map          | quantitative T2 map  |
| T2*         | T2star         | High resolution T2* image  |
| FLAIR       | FLAIR          |  |
| FLASH       | FLASH          |  |

| Name               | modality_label | Description  |
|--------------------|----------------|--|
| Proton density     | PD             |  |
| Proton density map | PDmap          |  |
| Combined PD/T2     | PDT2           |  |
| Inplane T1         | inplaneT1      | T1-weighted anatomical image matched to functional acquisition |
| Inplane T2         | inplaneT2      | T2-weighted anatomical image matched to functional acquisition |
| Angiography        | angio          |  |

### The run entity

If several scans of the same modality are acquired they **MUST** be indexed with a key-value pair: `_run-1`, `_run-2`, `_run-3` etc. (only integers are allowed as run labels). When there is only one scan of a given type the run key **MAY** be omitted. Please note that diffusion imaging data is stored elsewhere (see below).

### The acq entity

The OPTIONAL `acq-<label>` key/value pair corresponds to a custom label the user **MAY** use to distinguish a different set of parameters used for acquiring the same modality. For example this should be used when a study includes two T1w images - one full brain low resolution and one restricted field of view but high resolution. In such case two files could have the following names: `sub-01_acq-highres_T1w.nii.gz` and `sub-01_acq-lowres_T1w.nii.gz`, however the user is free to choose any other label than `highres` and `lowres` as long as they are consistent across subjects and sessions. In case different sequences are used to record the same modality (e.g. RARE and FLASH for T1w) this field can also be used to make that distinction. At what level of detail to make the distinction (e.g. just between RARE and FLASH, or between RARE, FLASH, and FLASHsubsampling) remains at the discretion of the researcher.

### The ce entity

Similarly the OPTIONAL `ce-<label>` key/value can be used to distinguish sequences using different contrast enhanced images. The label is the name of the contrast agent. The key `ContrastBolusIngredient` **MAY** be also be added in the JSON file, with the same label.

### The rec entity

Similarly the OPTIONAL `rec-<label>` key/value can be used to distinguish different reconstruction algorithms (for example ones using motion correction).

If the structural images included in the dataset were defaced (to protect identity of participants) one **CAN** provide the binary mask that was used to remove facial features in the form of `_defacemask` files. In such cases the OPTIONAL `mod-<label>` key/value pair corresponds to modality label for eg: T1w, inplaneT1, referenced by a defacemask image. E.g., `sub-01_mod-T1w_defacemask.nii.gz`.

Some meta information about the acquisition **MAY** be provided in an additional JSON file. See Common metadata fields for a list of terms and their definitions. There are also some OPTIONAL JSON fields specific to anatomical scans:

| Field name              | Definition  |
|-------------------------|---|
| ContrastBolusIngredient | OPTIONAL. Active ingredient of agent. Values <b>MUST</b> be one of: IODINE, GADOLINIUM, CARBON DIOXIDE, BARIUM, XENON Corresponds to DICOM Tag 0018,1048. |

## Task (including resting state) imaging data

Currently supported image contrasts include:



| Name  | contrast_label | Description   |
|-------|----------------|---|
| BOLD  | bold           | Blood-Oxygen-Level Dependent contrast (specialized T2* weighting)                                   |
| CBV   | cbv            | Cerebral Blood Volume contrast (specialized T2* weighting or difference between T1 weighted images) |
| Phase | phase          | Phase information associated with magnitude information stored in BOLD contrast                     |

Template:

```

1 sub-<label>/[ses-<label>/]
2   func/
3     sub-<label>[_ses-<label>]_task-<label>[_acq-<label>][_ce-<label>][_dir-<label>][_rec-<
4     label>][_run-<index>][_echo-<index>]<contrast_label>.nii.gz
      sub-<label>[_ses-<label>]_task-<label>[_acq-<label>][_ce-<label>][_dir-<label>][_rec-<
      label>][_run-<index>][_echo-<index>]_sbref.nii.gz

```

Imaging data acquired during functional imaging (i.e. imaging which supports rapid temporal repetition). This includes but is not limited to task based fMRI as well as resting state fMRI (i.e. rest is treated as another task). For task based fMRI a corresponding task events file (see below) **MUST** be provided (please note that this file is not necessary for resting state scans). For multiband acquisitions, one **MAY** also save the single-band reference image as type **sbref** (e.g. `sub-control01_task-nback_sbref.nii.gz`).

Each task has a unique label that **MUST** only consist of letters and/or numbers (other characters, including spaces and underscores, are not allowed). Those labels **MUST** be consistent across subjects and sessions.

If more than one run of the same task has been acquired a key/value pair: `_run-1`, `_run-2`, `_run-3` etc. **MUST** be used. If only one run was acquired the `run-<index>` can be omitted. In the context of functional imaging a run is defined as the same task, but in some cases it can mean different set of stimuli (for example randomized order) and participant responses.

The OPTIONAL `acq-<label>` key/value pair corresponds to a custom label one may use to distinguish different set of parameters used for acquiring the same task. For example this should be used when a study includes two resting state images - one single band and one multiband. In such case two files could have the following names: `sub-01_task-rest_acq-singleband_bold.nii.gz` and `sub-01_task-rest_acq-multiband_bold.nii.gz`, however the user is **MAY** choose any other label than **singleband** and **multiband** as long as they are consistent across subjects and sessions and consist only of the legal label characters.

Similarly the OPTIONAL `ce-<label>` key/value can be used to distinguish sequences using different contrast enhanced images. The label is the name of the contrast agent. The key ContrastBolusIngredient **MAY** be also be added in the JSON file, with the same label.

Similarly the OPTIONAL `rec-<label>` key/value can be used to distinguish different reconstruction algorithms (for example ones using motion correction).

Similarly the OPTIONAL `dir-<label>` and `rec-<label>` key/values can be used to distinguish different phase-encoding directions and reconstruction algorithms (for example ones using motion correction). See **fmap** Case 4 for more information on `dir` field specification.

Multi-echo data **MUST** be split into one file per echo. Each file shares the same name with the exception of the `_echo-<index>` key/value. For example:

```

1 sub-01/
2   func/
3     sub-01_task-cuedSGT_run-1_echo-1_bold.nii.gz
4     sub-01_task-cuedSGT_run-1_echo-1_bold.json
5     sub-01_task-cuedSGT_run-1_echo-2_bold.nii.gz
6     sub-01_task-cuedSGT_run-1_echo-2_bold.json
7     sub-01_task-cuedSGT_run-1_echo-3_bold.nii.gz
8     sub-01_task-cuedSGT_run-1_echo-3_bold.json

```

Please note that the `<index>` denotes the number/index (in a form of an integer) of the echo not the echo time value which needs to be stored in the field EchoTime of the separate JSON file.

Some meta information about the acquisition **MUST** be provided in an additional JSON file.

**Required fields**

| Field name     | Definition  |
|----------------|---|
| RepetitionTime | REQUIRED. The time in seconds between the beginning of an acquisition of one volume and the beginning of acquisition of the volume following it (TR). Please note that this definition includes time between scans (when no data has been acquired) in case of sparse acquisition schemes. This value needs to be consistent with the <code>pixdim[4]</code> field (after accounting for units stored in <code>xyzt_units</code> field) in the NIFTI header. This field is mutually exclusive with <code>VolumeTiming</code> and is derived from DICOM Tag 0018, 0080 and converted to seconds. |
| VolumeTiming   | REQUIRED. The time at which each volume was acquired during the acquisition. It is described using a list of times (in JSON format) referring to the onset of each volume in the BOLD series. The list must have the same length as the BOLD series, and the values must be non-negative and monotonically increasing. This field is mutually exclusive with <code>RepetitionTime</code> and <code>DelayTime</code> . If defined, this requires acquisition time (TA) be defined via either <code>SliceTiming</code> or <code>AcquisitionDuration</code> be defined.                            |
| TaskName       | REQUIRED. Name of the task. No two tasks should have the same name. Task label ( <code>task-</code> ) included in the file name is derived from this field by removing all non alphanumeric ( <code>[a-zA-Z0-9]</code> ) characters. For example task name <code>faces n-back</code> will corresponds to task label <code>facesnback</code> . A RECOMMENDED convention is to name resting state task using labels beginning with <code>rest</code> .  |

For the fields described above and in the following section, the term "Volume" refers to a reconstruction of the object being imaged (e.g., brain or part of a brain). In case of multiple channels in a coil, the term "Volume" refers to a combined image rather than an image from each coil.

**Other RECOMMENDED metadata****Timing Parameters**

| Field name                        | Definition  |
|-----------------------------------|---|
| NumberOfVolumesDiscardedByScanner | RECOMMENDED. Number of volumes ("dummy scans") discarded by the scanner (as opposed to those discarded by the user post hoc) before saving the imaging file. For example, a sequence that automatically discards the first 4 volumes before saving would have this field as 4. A sequence that doesn't discard dummy scans would have this set to 0. Please note that the onsets recorded in the <code>_event.tsv</code> file should always refer to the beginning of the acquisition of the first volume in the corresponding imaging file - independent of the value of <code>NumberOfVolumesDiscardedByScanner</code> field. |
| NumberOfVolumesDiscardedByUser    | RECOMMENDED. Number of volumes ("dummy scans") discarded by the user before including the file in the dataset. If possible, including all of the volumes is strongly recommended. Please note that the onsets recorded in the <code>_event.tsv</code> file should always refer to the beginning of the acquisition of the first volume in the corresponding imaging file - independent of the value of <code>NumberOfVolumesDiscardedByUser</code> field.   |
| DelayTime                         | RECOMMENDED. User specified time (in seconds) to delay the acquisition of data for the following volume. If the field is not present it is assumed to be set to zero. Corresponds to Siemens CSA header field <code>lDelayTimeInTR</code> . This field is REQUIRED for sparse sequences using the <code>RepetitionTime</code> field that do not have the <code>SliceTiming</code> field set to allowed for accurate calculation of "acquisition time". This field is mutually exclusive with <code>VolumeTiming</code> .  |

| Field name          | Definition   |
|---------------------|--|
| AcquisitionDuration | RECOMMENDED. Duration (in seconds) of volume acquisition. Corresponds to DICOM Tag 0018,9073 <b>Acquisition Duration</b> . This field is REQUIRED for sequences that are described with the <b>VolumeTiming</b> field and that do not have the <b>SliceTiming</b> field set to allow for accurate calculation of "acquisition time". This field is mutually exclusive with <b>RepetitionTime</b> . |
| DelayAfterTrigger   | RECOMMENDED. Duration (in seconds) from trigger delivery to scan onset. This delay is commonly caused by adjustments and loading times. This specification is entirely independent of <b>NumberOfVolumesDiscardedByScanner</b> or <b>NumberOfVolumesDiscardedByUser</b> , as the delay precedes the acquisition.   |

The following table recapitulates the different ways that specific fields have to be populated for functional sequences. Note that all these options can be used for non sparse sequences but that only options B, D and E are valid for sparse sequences.

|          | RepetitionTime | SliceTiming | AcquisitionDuration | DelayTime | VolumeTiming |
|----------|----------------|-------------|---------------------|-----------|--------------|
| option A | [ X ]          |             | [ ]                 |           | [ ]          |
| option B | [ ]            | [ X ]       |                     | [ ]       | [ X ]        |
| option C | [ ]            |             | [ X ]               | [ ]       | [ X ]        |
| option D | [ X ]          | [ X ]       | [ ]                 |           | [ ]          |
| option E | [ X ]          |             | [ ]                 | [ X ]     | [ ]          |

### Legend

- [ X ] --> has to be filled
- [ ] --> has to be left empty
- empty cell --> can be specified but not required

### fMRI task information

| Field name      | Definition   |
|-----------------|--|
| Instructions    | RECOMMENDED. Text of the instructions given to participants before the scan. This is especially important in context of resting state fMRI and distinguishing between eyes open and eyes closed paradigms. |
| TaskDescription | RECOMMENDED. Longer description of the task.   |
| CogAtlasID      | RECOMMENDED. URL of the corresponding <a href="#">Cognitive Atlas</a> Task term.   |
| CogPOID         | RECOMMENDED. URL of the corresponding <a href="#">CogPO</a> term.  |

See Common metadata fields for a list of additional terms and their definitions.

Example:

```
1 sub-control01/
2   func/
3     sub-control01_task-nback_bold.json
```

```
1 {
2   "TaskName": "N Back",
3   "RepetitionTime": 0.8,
4   "EchoTime": 0.03,
5   "FlipAngle": 78,
6   "SliceTiming": [0.0, 0.2, 0.4, 0.6, 0.0, 0.2, 0.4, 0.6, 0.0, 0.2, 0.4, 0.6, 0.0, 0.2, 0.4,
7     0.6],
8   "MultibandAccelerationFactor": 4,
9   "ParallelReductionFactorInPlane": 2,
```

```

9   "PhaseEncodingDirection": "j",
10  "InstitutionName": "Stanford University",
11  "InstitutionAddress": "450 Serra Mall, Stanford, CA 94305-2004, USA",
12  "DeviceSerialNumber": "11035"
13 }

```

If this information is the same for all participants, sessions and runs it can be provided in `task-<label>_bold.json` (in the root directory of the dataset). However, if the information differs between subjects/runs it can be specified in the `sub-<label>/func/sub-<label>_task-<label>[_acq-<label>][_run-<index>]_bold.json` file. If both files are specified fields from the file corresponding to a particular participant, task and run takes precedence.

## Diffusion imaging data

Template:

```

1 sub-<label>/[ses-<label>/]
2   dwi/
3     sub-<label>[_ses-<label>][_acq-<label>][_dir-<label>][_run-<index>]_dwi.nii.gz]
4     sub-<label>[_ses-<label>][_acq-<label>][_dir-<label>][_run-<index>]_dwi.bval
5     sub-<label>[_ses-<label>][_acq-<label>][_dir-<label>][_run-<index>]_dwi.bvec
6     sub-<label>[_ses-<label>][_acq-<label>][_dir-<label>][_run-<index>]_dwi.json
7     sub-<label>[_ses-<label>][_acq-<label>][_dir-<label>][_run-<index>]_sbref.nii.gz]
8     sub-<label>[_ses-<label>][_acq-<label>][_dir-<label>][_run-<index>]_sbref.json

```

Diffusion-weighted imaging data acquired for that participant. The OPTIONAL `acq-<label>` key/value pair corresponds to a custom label the user may use to distinguish different set of parameters. For example this should be used when a study includes two diffusion images - one single band and one multiband. In such case two files could have the following names: `sub-01_acq-singleband_dwi.nii.gz` and `sub-01_acq-multiband_dwi.nii.gz`, however the user is free to choose any other label than `singleband` and `multiband` as long as they are consistent across subjects and sessions. For multiband acquisitions, one can also save the single-band reference image as type `sbref` (e.g. `dwi/sub-control01_sbref.nii.gz`) The `bvec` and `bval` files are in the [FSL format](#): The `bvec` files contain 3 rows with `n` space-delimited floating-point numbers (corresponding to the `n` volumes in the relevant NIFTI file). The first row contains the `x` elements, the second row contains the `y` elements and third row contains the `z` elements of a unit vector in the direction of the applied diffusion gradient, where the `i`-th elements in each row correspond together to the `i`-th volume with `[0,0,0]` for non-diffusion-weighted volumes. Inherent to the FSL format for `bvec` specification is the fact that the coordinate system of the `bvecs` is with respect to the participant (i.e., defined by the axes of the corresponding `dwi.nii` file) and not the magnet's coordinate system, which means that any rotations applied to `dwi.nii` also need to be applied to the corresponding `bvec` file.

`bvec` example:

```

1 0 0 0.021828 -0.015425 -0.70918 -0.2465
2 0 0 0.80242 0.22098 -0.00063106 0.1043
3 0 0 -0.59636 0.97516 -0.70503 -0.96351

```

The `bval` file contains the `b`-values (in `s/mm2`) corresponding to the volumes in the relevant NIFTI file), with 0 designating non-diffusion-weighted volumes, space-delimited.

`bval` example:

```

1 0 0 2000 2000 1000 1000

```

`.bval` and `.bvec` files can be saved on any level of the directory structure and thus define those values for all sessions and/or subjects in one place (see Inheritance principle).

See Common metadata fields for a list of additional terms that can be included in the corresponding JSON file.

JSON example:

```

1 {
2   "PhaseEncodingDirection": "j-",
3   "TotalReadoutTime": 0.095
4 }

```

## Fieldmap data

Data acquired to correct for B0 inhomogeneities can come in different forms. The current version of this standard considers four different scenarios. Please note that in all cases fieldmap data can be linked to a specific scan(s) it was acquired for by filling the `IntendedFor` field in the corresponding JSON file. For example:

```
1 {
2   "IntendedFor": "func/sub-01_task-motor_bold.nii.gz"
3 }
```

The `IntendedFor` field may contain one or more filenames with paths relative to the subject subfolder. The path needs to use forward slashes instead of backward slashes. Here's an example with multiple target scans:

```
1 {
2   "IntendedFor": ["ses-pre/func/sub-01_ses-pre_task-motor_run-1_bold.nii.gz",
3                  "ses-post/func/sub-01_ses-post_task-motor_run-1_bold.nii.gz"]
4 }
```

The `IntendedFor` field is OPTIONAL and in case the fieldmaps do not correspond to any particular scans it does not have to be filled.

Multiple fieldmaps can be stored. In such case the `_run-1`, `_run-2` should be used. The OPTIONAL `acq-<label>` key/value pair corresponds to a custom label the user may use to distinguish different set of parameters.

### Case 1: Phase difference image and at least one magnitude image

Template:

```
1 sub-<label>/[ses-<label>/]
2   fmap/
3     sub-<label>[_ses-<label>][_acq-<label>][_run-<index>]_phasediff.nii.gz]
4     sub-<label>[_ses-<label>][_acq-<label>][_run-<index>]_phasediff.json
5     sub-<label>[_ses-<label>][_acq-<label>][_run-<index>]_magnitude1.nii.gz]
```

OPTIONAL

```
1 sub-<label>/[ses-<label>/]
2   fmap/
3     sub-<label>[_ses-<label>][_acq-<label>][_run-<index>]_magnitude2.nii.gz]
```

This is a common output for build in fieldmap sequence on Siemens scanners. In this particular case the sidcar JSON file has to define the Echo Times of the two phase images used to create the difference image. `EchoTime1` corresponds to the shorter echo time and `EchoTime2` to the longer echo time. Similarly `_magnitude1` image corresponds to the shorter echo time and the OPTIONAL `_magnitude2` image to the longer echo time. For example:

```
1 {
2   "EchoTime1": 0.00600,
3   "EchoTime2": 0.00746,
4   "IntendedFor": "func/sub-01_task-motor_bold.nii.gz"
5 }
```

### Case 2: Two phase images and two magnitude images

Template:

```
1 sub-<label>/[ses-<label>/]
2   fmap/
3     sub-<label>[_ses-<label>][_acq-<label>][_run-<index>]_phase1.nii.gz]
4     sub-<label>[_ses-<label>][_acq-<label>][_run-<index>]_phase1.json
5     sub-<label>[_ses-<label>][_acq-<label>][_run-<index>]_phase2.nii.gz]
6     sub-<label>[_ses-<label>][_acq-<label>][_run-<index>]_phase2.json
7     sub-<label>[_ses-<label>][_acq-<label>][_run-<index>]_magnitude1.nii.gz]
8     sub-<label>[_ses-<label>][_acq-<label>][_run-<index>]_magnitude2.nii.gz]
```

Similar to the case above, but instead of a precomputed phase difference map two separate phase images are presented. The two sidcar JSON files need to specify corresponding `EchoTime` values. For example:

```
1 {
2   "EchoTime": 0.00746,
3   "IntendedFor": "func/sub-01_task-motor_bold.nii.gz"
4 }
```

### Case 3: A real fieldmap image

Template:

```
1 sub-<label>/[ses-<label>/]
2   fmap/
3     sub-<label>[_ses-<label>][_acq-<label>][_run-<index>]_magnitude.nii.gz]
4     sub-<label>[_ses-<label>][_acq-<label>][_run-<index>]_fieldmap.nii.gz]
5     sub-<label>[_ses-<label>][_acq-<label>][_run-<index>]_fieldmap.json
```

In some cases (for example GE) the scanner software will output a precomputed fieldmap denoting the B0 inhomogeneities along with a magnitude image used for coregistration. In this case the sidcar JSON file needs to include the units of the fieldmap. The possible options are: `Hz`, `rad/s`, or `Tesla`. For example:

```
1 {
2   "Units": "rad/s",
3   "IntendedFor": "func/sub-01_task-motor_bold.nii.gz"
4 }
```

### Case 4: Multiple phase encoded directions ("pepolar")

Template:

```
1 sub-<label>/[ses-<label>/]
2   fmap/
3     sub-<label>[_ses-<label>][_acq-<label>][_ce-<label>]_dir-<label>[_run-<index>]_epi.nii.gz]
4     sub-<label>[_ses-<label>][_acq-<label>][_ce-<label>]_dir-<label>[_run-<index>]_epi.json
```

The phase-encoding polarity (PEpolar) technique combines two or more Spin Echo EPI scans with different phase encoding directions to estimate the underlying inhomogeneity/deformation map. Examples of tools using this kind of images are FSL TOPUP, AFNI 3dqwarp and SPM. In such a case, the phase encoding direction is specified in the corresponding JSON file as one of: `i`, `j`, `k`, `i-`, `j-`, `k-`. For these differentially phase encoded sequences, one also needs to specify the Total Readout Time defined as the time (in seconds) from the center of the first echo to the center of the last echo (aka "FSL definition" - see [here](#) and [here](#) how to calculate it). For example

```
1 {
2   "PhaseEncodingDirection": "j-",
3   "TotalReadoutTime": 0.095,
4   "IntendedFor": "func/sub-01_task-motor_bold.nii.gz"
5 }
```

`label` value of `_dir-` can be set to arbitrary alphanumeric label (`[a-zA-Z0-9]+` for example `LR` or `AP`) that can help users to distinguish between different files, but should not be used to infer any scanning parameters (such as phase encoding directions) of the corresponding sequence. Please rely only on the JSON file to obtain scanning parameters. `_epi` files can be a 3D or 4D - in the latter case all timepoints share the same scanning parameters. To indicate which run is intended to be used with which functional or diffusion scan the `IntendedFor` field in the JSON file should be used.

## Chapter 6

# Magnetoencephalography (MEG)

Support for MEG was developed as a BIDS Extension Proposal. Please cite the following paper when referring to this part of the standard in context of the academic literature:

Niso Galan, J.G., Gorgolewski, K.J., Bock, E., Brooks, T.L., Flandin, G., Gramfort, A., Henson, R.N., Jas, M., Litvak, V., Moreau, J., Oostenveld, R., Schoffelen, J.-M., Tadel, F., Wexler, J., Baillet, S. (2018). **MEG-BIDS, the brain imaging data structure extended to magnetoencephalography**. Scientific data, 5. doi: [10.1038/sdata.2018.110](https://doi.org/10.1038/sdata.2018.110)

## MEG recording data

Template:

```
1 sub-<label>/
2   [ses-<label>]/
3     meg/
4       sub-<label>[_ses-<label>]_task-<label>[_run-<index>] [_proc-<label>]_meg.<
5         manufacturer_specific_extension>
6       [sub-<label>[_ses-<label>]_task-<label>[_run-<index>] [_proc-<label>]_meg.json]
```

Unprocessed MEG data MUST be stored in the native file format of the MEG instrument with which the data was collected. With the MEG specification of BIDS, we wish to promote the adoption of good practices in the management of scientific data. Hence, the emphasis is not to impose a new, generic data format for the modality, but rather to standardize the way data is stored in repositories. Further, there is currently no widely accepted standard file format for MEG, but major software applications, including free and open-source solutions for MEG data analysis, provide readers of such raw files.

Some software readers may skip important metadata that is specific to MEG system manufacturers. It is therefore RECOMMENDED that users provide additional meta information extracted from the manufacturer raw data files in a sidecar JSON file. This allows for easy searching and indexing of key metadata elements without the need to parse files in proprietary data format. Other relevant files MAY be included alongside the MEG data; examples are provided below.

This template is for MEG data of any kind, including but not limited to task-based, resting-state, and noise recordings. If multiple Tasks were performed within a single Run, the task description can be set to **task-multitask**. The `_meg.json` SHOULD contain details on the Tasks. Some manufacturers data storage conventions use folders which contain data files of various nature: e.g., CTF's .ds format, or BTi/4D. Please refer to Appendix VI for examples from a selection of MEG manufacturers.

The `proc` label is analogous to `rec` for MR and denotes a variant of a file that was a result of particular processing performed on the device. This is useful for files produced in particular by Elekta's MaxFilter (e.g. sss, tsss, trans, quat, mc, etc.), which some installations impose to be run on raw data because of active shielding software corrections before the MEG data can actually be exploited.

## Sidecar JSON (`*_meg.json`)

Generic fields MUST be present:

| Field name | Definition   |
|------------|--|
| TaskName   | REQUIRED. Name of the task (for resting state use the <b>rest</b> prefix). Different Tasks SHOULD NOT have the same name. The Task label is derived from this field by removing all non alphanumeric ([a-zA-Z0-9]) characters. |

SHOULD be present: For consistency between studies and institutions, we encourage users to extract the values of these fields from the actual raw data. Whenever possible, please avoid using ad-hoc wording.

| Field name             | Definition   |
|------------------------|--|
| InstitutionName        | RECOMMENDED. The name of the institution in charge of the equipment that produced the composite instances.   |
| InstitutionAddress     | RECOMMENDED. The address of the institution in charge of the equipment that produced the composite instances.  |
| Manufacturer           | RECOMMENDED. Manufacturer of the MEG system (CTF, Elekta/Neuromag, BTi/4D, KIT/Yokogawa, ITAB, KRIS, Other). See Appendix VII with preferred names   |
| ManufacturersModelName | RECOMMENDED. Manufacturer's designation of the MEG scanner model (e.g. CTF-275). See Appendix VII with preferred names   |
| SoftwareVersions       | RECOMMENDED. Manufacturer's designation of the acquisition software.   |
| TaskDescription        | RECOMMENDED. Description of the task.  |
| Instructions           | RECOMMENDED. Text of the instructions given to participants before the scan. This is not only important for behavioral or cognitive tasks but also in resting state paradigms (e.g. to distinguish between eyes open and eyes closed).                     |
| CogAtlasID             | RECOMMENDED. URL of the corresponding Cognitive Atlas term that describes the task (e.g. Resting State with eyes closed " <a href="http://www.cognitiveatlas.org/task/id/trm_54e69c642d89b">http://www.cognitiveatlas.org/task/id/trm_54e69c642d89b</a> ") |
| CogPOID                | RECOMMENDED. URL of the corresponding CogPO term that describes the task (e.g. Rest " <a href="http://wiki.cogpo.org/index.php?title=Rest">http://wiki.cogpo.org/index.php?title=Rest</a> ")   |
| DeviceSerialNumber     | RECOMMENDED. The serial number of the equipment that produced the composite instances. A pseudonym can also be used to prevent the equipment from being identifiable, as long as each pseudonym is unique within the dataset.                              |

Specific MEG fields MUST be present:

| Field name          | Definition   |
|---------------------|--|
| SamplingFrequency   | REQUIRED. Sampling frequency (in Hz) of all the data in the recording, regardless of their type (e.g., 2400)   |
| PowerLineFrequency  | REQUIRED. Frequency (in Hz) of the power grid at the geographical location of the MEG instrument (i.e. 50 or 60)   |
| DewarPosition       | REQUIRED. Position of the dewar during the MEG scan: <b>upright</b> , <b>supine</b> or <b>degrees</b> of angle from vertical: for example on CTF systems, upright=15°, supine = 90°.   |
| SoftwareFilters     | REQUIRED. A JSON object of temporal software filters applied, or "n/a" if the data is not available. Each key:value pair in the JSON object is a name of the filter and an object in which its parameters are defined as key:value pairs. E.g.,<br>{"SSS": {"frame": "head", "badlimit": 7}, "SpatialCompensation": {"GradientOrder": "Order of the gradient compensation"}} |
| DigitizedLandmarks  | REQUIRED. Boolean ("true" or "false") value indicating whether anatomical landmark points (i.e. fiducials) are contained within this recording.  |
| DigitizedHeadPoints | REQUIRED. Boolean (true or false) value indicating whether head points outlining the scalp/face surface are contained within this recording.   |



SHOULD be present:

| Field name                 | Definition   |
|----------------------------|--|
| MEGChannelCount            | RECOMMENDED. Number of MEG channels (e.g. 275)   |
| MEGREFChannelCount         | RECOMMENDED. Number of MEG reference channels (e.g. 23). For systems without such channels (e.g. Neuromag Vectorview), <code>MEGREFChannelCount=0</code>   |
| EEGChannelCount            | RECOMMENDED. Number of EEG channels recorded simultaneously (e.g. 21)  |
| ECOGChannelCount           | RECOMMENDED. Number of ECoG channels   |
| SEEGChannelCount           | RECOMMENDED. Number of SEEG channels   |
| EOGChannelCount            | RECOMMENDED. Number of EOG channels  |
| ECGChannelCount            | RECOMMENDED. Number of ECG channels  |
| EMGChannelCount            | RECOMMENDED. Number of EMG channels  |
| MiscChannelCount           | RECOMMENDED. Number of miscellaneous analog channels for auxiliary signals   |
| TriggerChannelCount        | RECOMMENDED. Number of channels for digital (TTL bit level) triggers   |
| RecordingDuration          | RECOMMENDED. Length of the recording in seconds (e.g. 3600)  |
| RecordingType              | RECOMMENDED. Defines whether the recording is <b>continuous</b> or <b>epoched</b> ; this latter limited to time windows about events of interest (e.g., stimulus presentations, subject responses etc.)  |
| EpochLength                | RECOMMENDED. Duration of individual epochs in seconds (e.g. 1) in case of epoched data   |
| ContinuousHeadLocalization | RECOMMENDED. Boolean ( <b>true</b> or <b>false</b> ) value indicating whether continuous head localisation was performed.  |
| HeadCoilFrequency          | RECOMMENDED. List of frequencies (in Hz) used by the head localisation coils ('HLC' in CTF systems, 'HPI' in Elekta, 'COH' in BTi/4D) that track the subject's head position in the MEG helmet (e.g. [293, 307, 314, 321])   |
| MaxMovement                | RECOMMENDED. Maximum head movement (in mm) detected during the recording, as measured by the head localisation coils (e.g., 4.8)   |
| SubjectArtefactDescription | RECOMMENDED. Freeform description of the observed subject artefact and its possible cause (e.g. "Vagus Nerve Stimulator", "non-removable implant"). If this field is set to <b>n/a</b> , it will be interpreted as absence of major source of artifacts except cardiac and blinks.   |
| AssociatedEmptyRoom        | RECOMMENDED. Relative path in BIDS folder structure to empty-room file associated with the subject's MEG recording. The path needs to use forward slashes instead of backward slashes (e.g. <code>sub-emptyroom/ses-/meg/sub-emptyroom_ses-_task-noise_run-_meg.ds</code> ).   |
| HardwareFilters            | RECOMMENDED. A <b>JSON object</b> of temporal hardware filters applied, or <b>"n/a"</b> if the data is not available. Each key:value pair in the JSON object is a name of the filter and an object in which its parameters are defined as key:value pairs. E.g., <code>{"Highpass RC filter": {"Half amplitude cutoff (Hz)": 0.0159, "Roll-off": "6dB/Octave"}}</code> |

Specific EEG fields (if recorded with MEG) SHOULD be present:

| Field name                | Definition  |
|---------------------------|---|
| EEGPlacementScheme        | OPTIONAL. Placement scheme of EEG electrodes. Either the name of a standardised placement system (e.g., "10-20") or a list of standardised electrode names (e.g. ["Cz", "Pz"]). |
| CapManufacturer           | OPTIONAL. Manufacturer of the EEG cap (e.g. EasyCap)  |
| CapManufacturersModelName | OPTIONAL. Manufacturer's designation of the EEG cap model (e.g., M10)   |
| EEGReference              | OPTIONAL. Description of the type of EEG reference used (e.g., M1 for left mastoid, <b>average</b> , or <b>longitudinal bipolar</b> ).  |

By construct, EEG when recorded simultaneously with the same MEG system , should have the same

SamplingFrequency as MEG. Note that if EEG is recorded with a separate amplifier, it should be stored separately under a new /eeg data type (see the EEG specification).

Example:

```

1 {
2   "InstitutionName": "Stanford University",
3   "InstitutionAddress": "450 Serra Mall, Stanford, CA 94305-2004, USA",
4   "Manufacturer": "CTF",
5   "ManufacturersModelName": "CTF-275",
6   "DeviceSerialNumber": "11035",
7   "SoftwareVersions": "Acq 5.4.2-linux-20070507",
8   "PowerLineFrequency": 60,
9   "SamplingFrequency": 2400,
10  "MEGChannelCount": 270,
11  "MEGREFChannelCount": 26,
12  "EEGChannelCount": 0,
13  "EOGChannelCount": 2,
14  "ECGChannelCount": 1,
15  "EMGChannelCount": 0,
16  "DewarPosition": "upright",
17  "SoftwareFilters": {
18    "SpatialCompensation": {"GradientOrder": "3rd"}
19  },
20  "RecordingDuration": 600,
21  "RecordingType": "continuous",
22  "EpochLength": 0,
23  "TaskName": "rest",
24  "ContinuousHeadLocalization": true,
25  "HeadCoilFrequency": [1470,1530,1590],
26  "DigitizedLandmarks": true,
27  "DigitizedHeadPoints": true
28 }
```

Note that the date and time information SHOULD be stored in the Study key file (`scans.tsv`), see Scans file. As it is indicated there, date time information MUST be expressed in the following format YYYY-MM-DDThh:mm:ss (ISO8601 date-time format). For example: 2009-06-15T13:45:30. It does not need to be fully detailed, depending on local REB/IRB ethics board policy.

## Channels description (\*\_channels.tsv)

Template:

```

1 sub-<label>/
2   [ses-<label>]/
3   meg/
4   [sub-<label>[_ses-<label>]_task-<label>[_run-<index>] [_proc-<label>]_channels.tsv]
```

This file is RECOMMENDED as it provides easily searchable information across BIDS datasets for e.g., general curation, response to queries or batch analysis. To avoid confusion, the channels SHOULD be listed in the order they appear in the MEG data file. Missing values MUST be indicated with n/a.

The columns of the Channels description table stored in \*\_channels.tsv are:

MUST be present:

| Column name | Definition   |
|-------------|--|
| name        | REQUIRED. Channel name (e.g., MRT012, MEG023)  |
| type        | REQUIRED. Type of channel; MUST use the channel types listed below.  |
| units       | REQUIRED. Physical unit of the data values recorded by this channel in SI (see Appendix V: Units for allowed symbols). |

SHOULD be present:

| Column name        | Definition  |
|--------------------|---|
| description        | OPTIONAL. Brief free-text description of the channel, or other information of interest. See examples below.   |
| sampling_frequency | OPTIONAL. Sampling rate of the channel in Hz.   |
| low_cutoff         | OPTIONAL. Frequencies used for the high-pass filter applied to the channel in Hz. If no high-pass filter applied, use n/a.  |
| high_cutoff        | OPTIONAL. Frequencies used for the low-pass filter applied to the channel in Hz. If no low-pass filter applied, use n/a. Note that hardware anti-aliasing in A/D conversion of all MEG/EEG electronics applies a low-pass filter; specify its frequency here if applicable. |
| notch              | OPTIONAL. Frequencies used for the notch filter applied to the channel, in Hz. If no notch filter applied, use n/a.   |
| software_filters   | OPTIONAL. List of temporal and/or spatial software filters applied (e.g. "SSS", "SpatialCompensation"). Note that parameters should be defined in the general MEG sidcar .json file. Indicate n/a in the absence of software filters applied.                               |
| status             | OPTIONAL. Data quality observed on the channel (good/bad). A channel is considered <b>bad</b> if its data quality is compromised by excessive noise. Description of noise type SHOULD be provided in [status_description].  |
| status_description | OPTIONAL. Freeform text description of noise or artifact affecting data quality on the channel. It is meant to explain why the channel was declared bad in [status].  |

Example:

```

1 name type units description sampling_frequency low_cutoff high_cutoff notch software_filters
  status
2 UDI0001 TRIG V analogue trigger 1200 0.1 300 0 n/a good
3 MLC11 MEGGRADAXIAL T sensor 1st-order grad 1200 0 n/a 50 SSS bad

```

Restricted keyword list for field type

| Keyword          | Definition                       |
|------------------|----------------------------------|
| MEGMAG           | MEG magnetometer                 |
| MEGGRADAXIAL     | MEG axial gradiometer            |
| MEGGRADPLANAR    | MEG planargrاديometer            |
| MEGREFMAG        | MEG reference magnetometer       |
| MEGREFGRADAXIAL  | MEG reference axial gradiometer  |
| MEGREFGRADPLANAR | MEG reference planar gradiometer |
| MEGOTHER         | Any other type of MEG sensor     |
| EEG              | Electrode channel                |
| ECOG             | Electrode channel                |
| SEEG             | Electrode channel                |
| DBS              | Electrode channel                |
| VEOG             | Vertical EOG (electrooculogram)  |
| HEOG             | Horizontal EOG                   |
| EOG              | Generic EOG channel              |
| ECG              | ElectroCardioGram (heart)        |
| EMG              | ElectroMyoGram (muscle)          |
| TRIG             | System Triggers                  |
| AUDIO            | Audio signal                     |
| PD               | Photodiode                       |
| EYEGAZE          | Eye Tracker gaze                 |
| PUPIL            | Eye Tracker pupil diameter       |
| MISC             | Miscellaneous                    |

| Keyword  | Definition   |
|----------|--|
| SYSCLOCK | System time showing elapsed time since trial started |
| ADC      | Analog to Digital input                              |
| DAC      | Digital to Analog output                             |
| HLU      | Measured position of head and head coils             |
| FITERR   | Fit error signal from each head localization coil    |
| OTHER    | Any other type of channel                            |

Example of free text for field **description**

- stimulus, response, vertical EOG, horizontal EOG, skin conductance, sats, intracranial, eyetracker

Example:

```

1 name type units description
2 VEOG VEOG V vertical EOG
3 FDI EMG V left first dorsal interosseous
4 UDI0001 TRIG V analog trigger signal
5 UADC001 AUDIO V envelope of audio signal presented to participant

```

## Coordinate System JSON (\*\_coordsystem.json)

Template:

```

1 sub-<label>/
2   [ses-<label>]/
3   meg/
4   [sub-<label>[_ses-<label>][_acq-<label>]_coordsystem.json]

```

OPTIONAL. A JSON document specifying the coordinate system(s) used for the MEG, EEG, head localization coils, and anatomical landmarks.

MEG and EEG sensors:

| Field name                     | Description   |
|--------------------------------|---|
| MEGCoordinateSystem            | REQUIRED. Defines the coordinate system for the MEG sensors. See Appendix VIII: preferred names of Coordinate systems. If <b>Other</b> , provide definition of the coordinate system in [MEGCoordinateSystemDescription]. |
| MEGCoordinateUnits             | REQUIRED. Units of the coordinates of MEGCoordinateSystem. MUST be m, cm, or mm.  |
| MEGCoordinateSystemDescription | OPTIONAL. Freeform text description or link to document describing the MEG coordinate system system in detail.  |
| EEGCoordinateSystem            | OPTIONAL. Describes how the coordinates of the EEG sensors are to be interpreted.   |
| EEGCoordinateUnits             | OPTIONAL. Units of the coordinates of EEGCoordinateSystem. MUST be m, cm, or mm.  |
| EEGCoordinateSystemDescription | OPTIONAL. Freeform text description or link to document describing the EEG coordinate system system in detail.  |

Head localization coils:

| Field name                          | Description   |
|-------------------------------------|---|
| HeadCoilCoordinates                 | OPTIONAL. Key:value pairs describing head localization coil labels and their coordinates, interpreted following the <b>HeadCoilCoordinateSystem</b> , e.g., {NAS: [12.7,21.3,13.9], LPA: [5.2,11.3,9.6], RPA: [20.2,11.3,9.1]}. Note that coils are not always placed at locations that have a known anatomical name (e.g. for Elekta, Yokogawa systems); in that case generic labels can be used (e.g. {coil1: [12.2,21.3,12.3], coil2: [6.7,12.3,8.6], coil3: [21.9,11.0,8.1]} ). |
| HeadCoilCoordinateSystem            | OPTIONAL. Defines the coordinate system for the coils. See Appendix VIII: preferred names of Coordinate systems. If "Other", provide definition of the coordinate system in <b>HeadCoilCoordinateSystemDescription</b> .  |
| HeadCoilCoordinateUnits             | OPTIONAL. Units of the coordinates of <b>HeadCoilCoordinateSystem</b> . MUST be m, cm, or mm.   |
| HeadCoilCoordinateSystemDescription | OPTIONAL. Freeform text description or link to document describing the Head Coil coordinate system system in detail.  |

## Digitized head points:

| Field name                                     | Description  |
|--|--|
| DigitizedHeadPoints                            | OPTIONAL. Relative path to the file containing the locations of digitized head points collected during the session (e.g., <b>sub-01_headshape.pos</b> ). RECOMMENDED for all MEG systems, especially for CTF and BTi/4D. For Elekta/Neuromag the head points will be stored in the fif file. |
| DigitizedHeadPointsCoordinateSystem            | OPTIONAL. Defines the coordinate system for the digitized head points. See Appendix VIII: preferred names of Coordinate systems. If Other, provide definition of the coordinate system in <b>DigitizedHeadPointsCoordinateSystemDescription</b> .  |
| DigitizedHeadPointsCoordinateUnits             | OPTIONAL. Units of the coordinates of <b>DigitizedHeadPointsCoordinateSystem</b> . MUST be m, cm, or mm.   |
| DigitizedHeadPointsCoordinateSystemDescription | OPTIONAL. Freeform text description or link to document describing the Digitized head Points coordinate system system in detail.   |

## Anatomical MRI:

| Field name  | Description  |
|-------------|--|
| IntendedFor | OPTIONAL. Path or list of path relative to the subject subfolder pointing to the structural MRI, possibly of different types if a list is specified, to be used with the MEG recording. The path(s) need(s) to use forward slashes instead of backward slashes (e.g. <b>ses-/anat/sub-01_T1w.nii.gz</b> ). |

## Anatomical landmarks:

| Field name                                    | Description  |
|---|--|
| AnatomicalLandmarkCoordinates                 | OPTIONAL. Key:value pairs of the labels and 3-D digitized locations of anatomical landmarks, interpreted following the <code>AnatomicalLandmarkCoordinateSystem</code> , e.g., <code>{"NAS": [12.7, 21.3, 13.9], "LPA": [5.2, 11.3, 9.6], "RPA": [20.2, 11.3, 9.1]}</code> . |
| AnatomicalLandmarkCoordinateSystem            | OPTIONAL. Defines the coordinate system for the anatomical landmarks. See Appendix VIII: preferred names of Coordinate systems. If <code>Other</code> , provide definition of the coordinate system in <code>AnatomicalLandmarkCoordinateSystemDescription</code> .          |
| AnatomicalLandmarkCoordinateUnits             | OPTIONAL. Units of the coordinates of <code>AnatomicalLandmarkCoordinateSystem</code> . MUST be m, cm, or mm.  |
| AnatomicalLandmarkCoordinateSystemDescription | OPTIONAL. Freeform text description or link to document describing the Head Coil coordinate system in detail.  |

It is also RECOMMENDED that the MRI voxel coordinates of the actual anatomical landmarks for co-registration of MEG with structural MRI are stored in the `AnatomicalLandmarkCoordinates` field in the JSON sidecar of the corresponding T1w MRI anatomical data of the subject seen in the MEG session (see here) - for example: `sub-01/ses-mri/anat/sub-01_ses-mri_acq-mprage_T1w.json`

In principle, these locations are those of absolute anatomical markers. However, the marking of NAS, LPA and RPA is more ambiguous than that of e.g., AC and PC. This may result in some variability in their 3-D digitization from session to session, even for the same participant. The solution would be to use only one T1w file and populate the `AnatomicalLandmarkCoordinates` field with session-specific labels e.g., `"NAS-session1": [127, 213, 139], "NAS-session2": [123, 220, 142]`, etc.

Fiducials information:

| Field name           | Description  |
|----------------------|--|
| FiducialsDescription | OPTIONAL. A freeform text field documenting the anatomical landmarks that were used and how the head localization coils were placed relative to these. This field can describe, for instance, whether the true anatomical locations of the left and right pre-auricular points were used and digitized, or rather whether they were defined as the intersection between the tragus and the helix (the entry of the ear canal), or any other anatomical description of selected points in the vicinity of the ears. |

For more information on the definition of anatomical landmarks, please visit: [http://www.fieldtriptoolbox.org/faq/how\\_are\\_the\\_lpa\\_and\\_rpa\\_points\\_defined](http://www.fieldtriptoolbox.org/faq/how_are_the_lpa_and_rpa_points_defined)

For more information on typical coordinate systems for MEG-MRI coregistration: [http://www.fieldtriptoolbox.org/faq/how\\_are\\_the\\_different\\_head\\_and\\_mri\\_coordinate\\_systems\\_defined](http://www.fieldtriptoolbox.org/faq/how_are_the_different_head_and_mri_coordinate_systems_defined), or: <http://neuroimage.usc.edu/brainstorm/CoordinateSystems>

## Landmark photos (\*\_photo.jpg)

Photos of the anatomical landmarks and/or head localization coils (\*\_photo.jpg)

Template:


```

1 sub-<label>/
2   [ses-<label>]/
3   meg/
4   [sub-<label>[_ses-<label>][_acq-<label>]_photo.jpg]
```

Photos of the anatomical landmarks and/or head localization coils on the subject's head are RECOMMENDED. If the coils are not placed at the location of actual anatomical landmarks, these latter may be marked with a piece of felt-tip taped to the skin. Please note that the photos may need to be cropped or blurred to conceal identifying features prior to sharing, depending on the terms of the consent given by the participant.

The `acq` parameter can be used to indicate acquisition of different photos of the same face (or other body part in different angles to show, for example, the location of the nasion (NAS) as opposed to the right periauricular point (RPA)).

Example of the NAS fiducial placed between the eyebrows, rather than at the actual anatomical nasion: `sub-0001_ses-001_acq-NAS_photo.jpg`

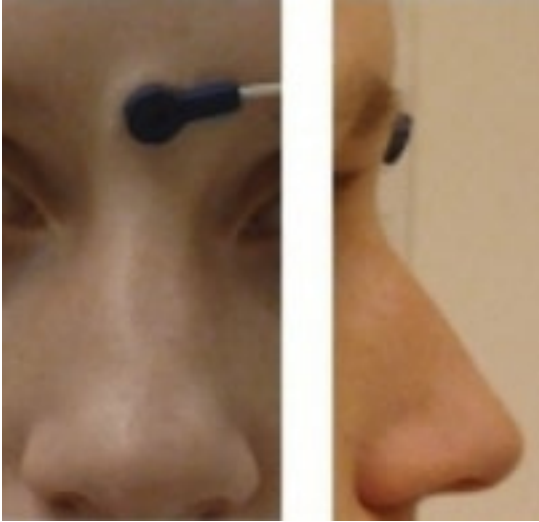

## Head shape and electrode description (\*\_headshape.<ext>)

Template:

```
1 sub-<label>/
2   [ses-<label>]/
3   meg/
4   [sub-<label>[_ses-<label>][_acq-<label>]_headshape.<manufacturer_specific_extension>]
```

This file is RECOMMENDED.

The 3-D locations of points that describe the head shape and/or EEG electrode locations can be digitized and stored in separate files. The `*_acq-<label>` can be used when more than one type of digitization is done for a session, for example when the head points are in a separate file from the EEG locations. These files are stored in the specific format of the 3-D digitizer's manufacturer (see Appendix VI).

Example:

```
1 sub-control01
2   ses-01
3       sub-control01_ses-01_acq-HEAD_headshape.pos
4       sub-control01_ses-01_acq-ECG_headshape.pos
```

Note that the `*_headshape` file(s) is shared by all the runs and tasks in a session. If the subject needs to be taken out of the scanner and the head-shape has to be updated, then for MEG it could be considered to be a new session.

## Empty-room MEG recordings

Empty-room MEG recordings capture the environment and system noise. Their collection is RECOMMENDED, before/during/after each session. This data is stored inside a subject folder named `sub-emptyroom`. The `session label` SHOULD be that of the date of the empty-room recording (e.g. `ses-YYYYMMDD`). The `scans.tsv` file containing the date/time of the acquisition SHOULD also be included. Hence, users will be able to

retrieve the empty-room recording that best matches a particular session with a participant, based on date/time of recording.

Example:

```
1 sub-control01/  
2 sub-control02/  
3 sub-emptyroom/  
4   ses-20170801/  
5     sub-emptyroom_ses-20170801_scans.tsv  
6     meg/  
7       sub-emptyroom_ses-20170801_task-noise_meg.ds  
8       sub-emptyroom_ses-20170801_task-noise_meg.json
```

TaskName in the \*\_meg.json file should be set to "noise".



## Chapter 7

# Electroencephalography (EEG)

Support for EEG was developed as a BIDS Extension Proposal. Please cite the following paper when referring to this part of the standard in context of the academic literature:

Pernet, C. R., Appelhoff, S., Gorgolewski, K.J., Flandin, G., Phillips, C., Delorme, A., Oostenveld, R. (2019). **EEG-BIDS, an extension to the brain imaging data structure for electroencephalography**. Scientific data, 6. doi: [10.1038/s41597-019-0104-8](https://doi.org/10.1038/s41597-019-0104-8)

The following example EEG datasets have been formatted using this specification and can be used for practical guidance when curating a new dataset.

- Single session per subject: [eeg\\_matchingpennies](#)
- Multiple sessions per subject: [eeg\\_rishikesh](#)
- Combined with fMRI: [eeg\\_rest\\_fmri](#)

Further datasets are available from the [BIDS examples repository](#).

## EEG recording data

Template:

```
1 sub-<label>/
2   [ses-<label>]/
3   eeg/
4     sub-<label>[_ses-<label>]_task-<label>[_run-<index>]_eeg.<
       manufacturer_specific_extension>
5     sub-<label>[_ses-<label>]_task-<label>[_run-<index>]_eeg.json
```

The EEG community uses a variety of formats for storing raw data, and there is no single standard that all researchers agree on. For BIDS, EEG data **MUST** be stored in one of the following formats:

- [European data format](#) ([.edf](#))
- [BrainVision Core Data Format](#) ([.vhdr](#), [.vmrk](#), [.eeg](#)) by Brain Products GmbH
- The format used by the MATLAB toolbox [EEGLAB](#) ([.set](#) and [.fdt](#) files)
- [Biosemi](#) data format ([.bdf](#))

It is **RECOMMENDED** to use the European data format, or the BrainVision data format. It is furthermore discouraged to use the other accepted formats over these **RECOMMENDED** formats, particularly because there are conversion scripts available in most commonly used programming languages to convert data into the **RECOMMENDED** formats. The data in their original format, if different from the supported formats, can be stored in the `/sourcedata` directory.

The original data format is especially valuable in case conversion elicits the loss of crucial metadata specific to manufacturers and specific EEG systems. We also encourage users to provide additional meta information extracted from the manufacturer specific data files in the sidecar JSON file. Other relevant files **MAY** be included alongside the original EEG data in `/sourcedata`.

Note the `RecordingType`, which depends on whether the data stream on disk is interrupted or not. Continuous data is by definition 1 segment without interruption. Epoched data consists of multiple segments that all have the same length (e.g., corresponding to trials) and that have gaps in between. Discontinuous data consists of multiple segments of different length, for example due to a pause in the acquisition.

Note that for proper documentation of EEG recording metadata it is important to understand the difference between electrode and channel: An EEG electrode is attached to the skin, whereas a channel is the combination of the analog differential amplifier and analog-to-digital converter that result in a potential (voltage) difference that is stored in the EEG dataset. We employ the following short definitions:

- **Electrode** = A single point of contact between the acquisition system and the recording site (e.g., scalp, neural tissue, ...). Multiple electrodes can be organized as caps (for EEG), arrays, grids, leads, strips, probes, shafts, etc.
- **Channel** = A single analog-to-digital converter in the recording system that regularly samples the value of a transducer, which results in the signal being represented as a time series in the digitized data. This can be connected to two electrodes (to measure the potential difference between them), a magnetic field or magnetic gradient sensor, temperature sensor, accelerometer, etc.

Although the *reference* and *ground* electrodes are often referred to as channels, they are in most common EEG systems not recorded by themselves. Therefore they are not represented as channels in the data. The type of referencing for all channels and optionally the location of the reference electrode and the location of the ground electrode MAY be specified.

## Sidecar JSON (\*\_eeg.json)

Generic fields MUST be present:

| Field name | Definition   |
|------------|--|
| TaskName   | REQUIRED. Name of the task (for resting state use the <b>rest</b> prefix). Different Tasks SHOULD NOT have the same name. The Task label is derived from this field by removing all non alphanumeric ([a-zA-Z0-9]) characters. |

SHOULD be present: For consistency between studies and institutions, we encourage users to extract the values of these fields from the actual raw data. Whenever possible, please avoid using ad hoc wording.

| Field name             | Definition   |
|------------------------|--|
| InstitutionName        | RECOMMENDED. The name of the institution in charge of the equipment that produced the composite instances.   |
| InstitutionAddress     | RECOMMENDED. The address of the institution in charge of the equipment that produced the composite instances.  |
| Manufacturer           | RECOMMENDED. Manufacturer of the EEG system (e.g., Biosemi, Brain Products, Neuroscan).  |
| ManufacturersModelName | RECOMMENDED. Manufacturer's designation of the EEG system model (e.g., BrainAmp DC).   |
| SoftwareVersions       | RECOMMENDED. Manufacturer's designation of the acquisition software.   |
| TaskDescription        | RECOMMENDED. Description of the task.  |
| Instructions           | RECOMMENDED. Text of the instructions given to participants before the scan. This is not only important for behavioral or cognitive tasks but also in resting state paradigms (e.g., to distinguish between eyes open and eyes closed).                                      |
| CogAtlasID             | RECOMMENDED. URL of the corresponding <a href="#">Cognitive Atlas</a> term that describes the task (e.g., Resting State with eyes closed " <a href="http://www.cognitiveatlas.org/task/id/trm_54e69c642d89b">http://www.cognitiveatlas.org/task/id/trm_54e69c642d89b</a> "). |
| CogPOID                | RECOMMENDED. URL of the corresponding <a href="#">CogPO</a> term that describes the task (e.g., Rest " <a href="http://wiki.cogpo.org/index.php?title=Rest">http://wiki.cogpo.org/index.php?title=Rest</a> ") .  |

| Field name         | Definition  |
|--------------------|---|
| DeviceSerialNumber | RECOMMENDED. The serial number of the equipment that produced the composite instances. A pseudonym can also be used to prevent the equipment from being identifiable, as long as each pseudonym is unique within the dataset. |

Specific EEG fields **MUST** be present:

| Field name         | Definition   |
|--------------------|--|
| EEGReference       | REQUIRED. General description of the reference scheme used and (when applicable) of location of the reference electrode in the raw recordings (e.g., "left mastoid", "Cz", "CMS"). If different channels have a different reference, this field should have a general description and the channel specific reference should be defined in the <code>_channels.tsv</code> file. |
| SamplingFrequency  | REQUIRED. Sampling frequency (in Hz) of all the data in the recording, regardless of their type (e.g., 2400).  |
| PowerLineFrequency | REQUIRED. Frequency (in Hz) of the power grid at the geographical location of the EEG instrument (i.e., 50 or 60).   |
| SoftwareFilters    | REQUIRED. A <b>JSON object</b> of temporal software filters applied, or <b>"n/a"</b> if the data is not available. Each key:value pair in the JSON object is a name of the filter and an object in which its parameters are defined as key:value pairs. E.g., <code>{"Anti-aliasing filter": {"half-amplitude cutoff (Hz)": 500, "Roll-off": "6dB/Octave"}}</code> .           |

**SHOULD** be present:

| Field name                | Definition  |
|---------------------------|---|
| CapManufacturer           | RECOMMENDED. Name of the cap manufacturer (e.g., "EasyCap").  |
| CapManufacturersModelName | RECOMMENDED. Manufacturer's designation of the EEG cap model (e.g., "actiCAP 64 Ch Standard-2").  |
| EEGChannelCount           | RECOMMENDED. Number of EEG channels included in the recording (e.g., 128).  |
| ECGChannelCount           | RECOMMENDED. Number of ECG channels.  |
| EMGChannelCount           | RECOMMENDED. Number of EMG channels.  |
| EOGChannelCount           | RECOMMENDED. Number of EOG channels.  |
| MiscChannelCount          | RECOMMENDED. Number of miscellaneous analog channels for auxiliary signals.   |
| TriggerChannelCount       | RECOMMENDED. Number of channels for digital (TTL bit level) trigger.  |
| RecordingDuration         | RECOMMENDED. Length of the recording in seconds (e.g., 3600).   |
| RecordingType             | RECOMMENDED. Defines whether the recording is <b>continuous</b> , <b>discontinuous</b> or <b>epoched</b> .  |
| EpochLength               | RECOMMENDED. Duration of individual epochs in seconds (e.g., 1) in case of epoched data.  |
| EEGGround                 | RECOMMENDED. Description of the location of the ground electrode (e.g., "placed on right mastoid (M2)").  |
| HeadCircumference         | RECOMMENDED. Circumference of the participants head, expressed in cm (e.g., 58).  |
| EEGPlacementScheme        | RECOMMENDED. Placement scheme of EEG electrodes. Either the name of a standardized placement system (e.g., "10-20") or a list of standardized electrode names (e.g., ["Cz", "Pz"]). |

| Field name                 | Definition  |
|----------------------------|---|
| HardwareFilters            | RECOMMENDED. A <a href="#">JSON object</a> of temporal hardware filters applied, or "n/a" if the data is not available. Each key:value pair in the JSON object is a name of the filter and an object in which its parameters are defined as key:value pairs. E.g., {"Highpass RC filter": {"Half amplitude cutoff (Hz)": 0.0159, "Roll-off": "6dB/Octave"}} |
| SubjectArtefactDescription | RECOMMENDED. Free-form description of the observed subject artifact and its possible cause (e.g., "Vagus Nerve Stimulator", "non-removable implant"). If this field is set to n/a, it will be interpreted as absence of major source of artifacts except cardiac and blinks.  |

Example:

```

1 {
2   "TaskName": "Seeing stuff",
3   "TaskDescription": "Subjects see various images for which phase, amplitude spectrum, and color
   vary continuously",
4   "Instructions": "Your task is to detect images when they appear for the 2nd time, only then
   press the response button with your right/left hand (counterbalanced across subjects)",
5   "InstitutionName": "The world best university, 10 Beachfront Avenue, Papeete",
6   "SamplingFrequency": 2400,
7   "Manufacturer": "Brain Products",
8   "ManufacturersModelName": "BrainAmp DC",
9   "CapManufacturer": "EasyCap",
10  "CapManufacturersModelName": "M1-ext",
11  "EEGChannelCount": 87,
12  "EOGChannelCount": 2,
13  "ECGChannelCount": 1,
14  "EMGChannelCount": 0,
15  "MiscChannelCount": 0,
16  "TriggerChannelCount": 1,
17  "PowerLineFrequency": 50,
18  "EEGPlacementScheme": "10 percent system",
19  "EEGReference": "single electrode placed on FCz",
20  "EEGGround": "placed on AFz",
21  "SoftwareFilters": {
22    "Anti-aliasing filter": {
23      "half-amplitude cutoff (Hz)": 500,
24      "Roll-off": "6dB/Octave"
25    }
26  },
27  "HardwareFilters": {
28    "ADC's decimation filter (hardware bandwidth limit)": {
29      "-3dB cutoff point (Hz)": 480,
30      "Filter order sinc response": 5
31    }
32  },
33  "RecordingDuration": 600,
34  "RecordingType": "continuous"
35 }

```

Note that the date and time information SHOULD be stored in the Study key file (`scans.tsv`). As it is indicated there, date time information MUST be expressed in the following format `YYYY-MM-DDThh:mm:ss` ([ISO8601](#) date-time format). For example: 2009-06-15T13:45:30. It does not need to be fully detailed, depending on local REB/IRB ethics board policy.

## Channels description (\*\_channels.tsv)

Template:

```
1 sub-<label>/
2   [ses-<label>]/
3   eeg/
4   [sub-<label>[_ses-<label>]_task-<label>[_run-<index>]]_channels.tsv]
```

This file is RECOMMENDED as it provides easily searchable information across BIDS datasets for e.g., general curation, response to queries or batch analysis. The required columns are channel **name**, **type** and **units** in this specific order. To avoid confusion, the channels SHOULD be listed in the order they appear in the EEG data file. Any number of additional columns may be added to provide additional information about the channels. Note that electrode positions SHOULD NOT be added to this file, but to \*\_electrodes.tsv.

The columns of the Channels description table stored in \*\_channels.tsv are:

MUST be present:

| Column name | Definition   |
|-------------|--|
| name        | REQUIRED. Channel name (e.g., FC1, Cz)   |
| type        | REQUIRED. Type of channel; MUST use the channel types listed below.  |
| units       | REQUIRED. Physical unit of the data values recorded by this channel in SI units (see Appendix V: Units for allowed symbols). |

SHOULD be present:

| Column name        | Definition  |
|--------------------|---|
| description        | OPTIONAL. Free-form text description of the channel, or other information of interest. See examples below.  |
| sampling_frequency | OPTIONAL. Sampling rate of the channel in Hz.   |
| reference          | OPTIONAL. Name of the reference electrode(s) (not needed when it is common to all channels, in that case it can be specified in *_eeg.json as EEGReference).  |
| low_cutoff         | OPTIONAL. Frequencies used for the high-pass filter applied to the channel in Hz. If no high-pass filter applied, use n/a.  |
| high_cutoff        | OPTIONAL. Frequencies used for the low-pass filter applied to the channel in Hz. If no low-pass filter applied, use n/a. Note that hardware anti-aliasing in A/D conversion of all EEG electronics applies a low-pass filter; specify its frequency here if applicable. |
| notch              | OPTIONAL. Frequencies used for the notch filter applied to the channel, in Hz. If no notch filter applied, use n/a.   |
| status             | OPTIONAL. Data quality observed on the channel (good/bad). A channel is considered <b>bad</b> if its data quality is compromised by excessive noise. Description of noise type SHOULD be provided in [status_description].  |
| status_description | OPTIONAL. Free-form text description of noise or artifact affecting data quality on the channel. It is meant to explain why the channel was declared bad in [status].   |

Restricted keyword list for field **type** in alphabetic order (shared with the MEG and iEEG modality; however, only the types that are common in EEG data are listed here):

| Keyword | Description  |
|---------|--|
| AUDIO   | Audio signal   |
| EEG     | Electroencephalogram channel                                 |
| EOG     | Generic electrooculogram (eye), different from HEOG and VEOG |

| Keyword  | Description  |
|----------|--|
| ECG      | Electrocardiogram (heart)                            |
| EMG      | Electromyogram (muscle)                              |
| EYEGAZE  | Eye tracker gaze                                     |
| GSR      | Galvanic skin response                               |
| HEOG     | Horizontal EOG (eye)                                 |
| MISC     | Miscellaneous  |
| PUPIL    | Eye tracker pupil diameter                           |
| REF      | Reference channel                                    |
| RESP     | Respiration  |
| SYSCLOCK | System time showing elapsed time since trial started |
| TEMP     | Temperature  |
| TRIG     | System triggers                                      |
| VEOG     | Vertical EOG (eye)                                   |

Example of free-form text for field **description**

- n/a, stimulus, response, skin conductance, battery status

Example:

| 1 | name    | type | units  | description                    | status | status_description   |
|---|---------|------|--------|--------------------------------|--------|----------------------|
| 2 | VEOG    | VEOG | microV | n/a                            | good   | n/a                  |
| 3 | FDI     | EMG  | microV | left first dorsal interosseous | good   | n/a                  |
| 4 | Cz      | EEG  | microV | n/a                            | bad    | high frequency noise |
| 5 | UADC001 | MISC | n/a    | envelope of audio signal       | good   | n/a                  |

## Electrodes description (\*\_electrodes.tsv)

Template:

```

1 sub-<label>/
2   [ses-<label>]/
3   eeg/
4   [sub-<label>[_ses-<label>][_acq-<label>][_run-<index>]_electrodes.tsv]
```

File that gives the location of EEG electrodes. Note that coordinates are expected in cartesian coordinates according to the EEGCoordinateSystem and EEGCoordinateSystemUnits fields in \*\_coordsystem.json. **If an \*\_electrodes.tsv file is specified, a \*\_coordsystem.json file MUST be specified as well.** The order of the required columns in the \*\_electrodes.tsv file MUST be as listed below.

MUST be present:

| Column name | Definition                                   |
|-------------|--|
| name        | REQUIRED. Name of the electrode              |
| x           | REQUIRED. Recorded position along the x-axis |
| y           | REQUIRED. Recorded position along the y-axis |
| z           | REQUIRED. Recorded position along the z-axis |

SHOULD be present:

| Column name | Definition  |
|-------------|---|
| type        | RECOMMENDED. Type of the electrode (e.g., cup, ring, clip-on, wire, needle) |
| material    | RECOMMENDED. Material of the electrode, e.g., Tin, Ag/AgCl, Gold            |
| impedance   | RECOMMENDED. Impedance of the electrode in kOhm                             |

Example:

| 1 | name | x       | y       | z       | type    | material |
|---|------|---------|---------|---------|---------|----------|
| 2 | A1   | -0.0707 | 0.0000  | -0.0707 | clip-on | Ag/AgCl  |
| 3 | F3   | -0.0567 | 0.0677  | 0.0469  | cup     | Ag/AgCl  |
| 4 | Fz   | 0.0000  | 0.0714  | 0.0699  | cup     | Ag/AgCl  |
| 5 | REF  | -0.0742 | -0.0200 | -0.0100 | cup     | Ag/AgCl  |
| 6 | GND  | 0.0742  | -0.0200 | -0.0100 | cup     | Ag/AgCl  |

The `acq` parameter can be used to indicate acquisition of the same data. For example, this could be the recording of electrode positions with a different electrode position recording device, or repeated digitization before and after the recording.

## Coordinate System JSON (\*\_coordsystem.json)

Template:

```
1 sub-<label>/
2   [ses-<label>]/
3   eeg/
4   [sub-<label>[_ses-<label>][_acq-<label>]_coordsystem.json]
```

A `*_coordsystem.json` file is used to specify the fiducials, the location of anatomical landmarks, and the coordinate system and units in which the position of electrodes and landmarks is expressed. **The `*_coordsystem.json` is REQUIRED if the optional `*_electrodes.tsv` is specified.** If a corresponding anatomical MRI is available, the locations of landmarks and fiducials according to that scan should also be stored in the `*_T1w.json` file which goes alongside the MRI data.

For disambiguation, we employ the following definitions for fiducials and anatomical landmarks respectively:

- Fiducials = objects with a well defined location used to facilitate the localization of electrodes and coregistration with other geometric data such as the participant's own T1 weighted magnetic resonance head image, a T1 weighted template head image, or a spherical head model. Commonly used fiducials are vitamin-E pills, which show clearly in an MRI, or reflective spheres that are localized with an infrared optical tracking system.
- Anatomical landmarks = locations on a research subject such as the nasion, which is the intersection of the frontal bone and two nasal bones of the human skull.

Fiducials are typically used in conjunction with anatomical landmarks. An example would be the placement of vitamin-E pills on top of anatomical landmarks, or the placement of LEDs on the nasion and preauricular points to triangulate the position of other LED-lit electrodes on a research subject's head.

- For more information on the definition of anatomical landmarks, please visit: [http://www.fieldtriptoolbox.org/faq/how\\_are\\_the\\_lpa\\_and\\_rpa\\_points\\_defined](http://www.fieldtriptoolbox.org/faq/how_are_the_lpa_and_rpa_points_defined)
- For more information on coordinate systems for coregistration, please visit: [http://www.fieldtriptoolbox.org/faq/how\\_are\\_the\\_different\\_head\\_and\\_mri\\_coordinate\\_systems\\_defined](http://www.fieldtriptoolbox.org/faq/how_are_the_different_head_and_mri_coordinate_systems_defined)

General fields:

| Keyword     | Description  |
|-------------|--|
| IntendedFor | OPTIONAL. Relative path to associate the electrodes, landmarks and fiducials to an MRI/CT. |

Fields relating to the EEG electrode positions:

| Keyword             | Description  |
|---------------------|--|
| EEGCoordinateSystem | REQUIRED. Refers to the coordinate system in which the EEG electrode positions are to be interpreted (see Appendix VIII).                  |
| EEGCoordinateUnits  | REQUIRED. Units in which the coordinates that are listed in the field <code>EEGCoordinateSystem</code> are represented (e.g., "mm", "cm"). |

| Keyword                        | Description   |
|--------------------------------|---|
| EEGCoordinateSystemDescription | RECOMMENDED. Free-form text description of the coordinate system. May also include a link to a documentation page or paper describing the system in greater detail. |

Fields relating to the position of fiducials measured during an EEG session/run:

| Keyword                              | Description   |
|--------------------------------------|---|
| FiducialsDescription                 | OPTIONAL. Free-form text description of how the fiducials such as vitamin-E capsules were placed relative to anatomical landmarks, and how the position of the fiducials were measured (e.g., both with Polhemus and with T1w MRI).                             |
| FiducialsCoordinates                 | RECOMMENDED. Key:value pairs of the labels and 3-D digitized position of anatomical landmarks, interpreted following the <code>FiducialsCoordinateSystem</code> (e.g., <code>{"NAS": [12.7,21.3,13.9], "LPA": [5.2,11.3,9.6], "RPA": [20.2,11.3,9.1]}</code> ). |
| FiducialsCoordinateSystem            | RECOMMENDED. Refers to the coordinate space to which the landmarks positions are to be interpreted - preferably the same as the <code>EEGCoordinateSystem</code> .  |
| FiducialsCoordinateUnits             | RECOMMENDED. Units in which the coordinates that are listed in the field <code>AnatomicalLandmarkCoordinateSystem</code> are represented (e.g., "mm", "cm").  |
| FiducialsCoordinateSystemDescription | RECOMMENDED. Free-form text description of the coordinate system. May also include a link to a documentation page or paper describing the system in greater detail.   |

Fields relating to the position of anatomical landmark measured during an EEG session/run:

| Keyword                                       | Description  |
|---|--|
| AnatomicalLandmarkCoordinates                 | RECOMMENDED. Key:value pairs of the labels and 3-D digitized position of anatomical landmarks, interpreted following the <code>AnatomicalLandmarkCoordinateSystem</code> (e.g., <code>{"NAS": [12.7,21.3,13.9], "LPA": [5.2,11.3,9.6], "RPA": [20.2,11.3,9.1]}</code> ). |
| AnatomicalLandmarkCoordinateSystem            | RECOMMENDED. Refers to the coordinate space to which the landmarks positions are to be interpreted - preferably the same as the <code>EEGCoordinateSystem</code> .   |
| AnatomicalLandmarkCoordinateUnits             | RECOMMENDED. Units in which the coordinates that are listed in the field <code>AnatomicalLandmarkCoordinateSystem</code> are represented (e.g., "mm", "cm").   |
| AnatomicalLandmarkCoordinateSystemDescription | RECOMMENDED. Free-form text description of the coordinate system. May also include a link to a documentation page or paper describing the system in greater detail.  |

If the position of anatomical landmarks is measured using the same system or device used to measure electrode positions, and if thereby the anatomical landmarks are expressed in the same coordinates, the coordinates of the anatomical landmarks can be specified in `electrodes.tsv`. The same applies to the coordinates of the fiducials.

Anatomical landmarks or fiducials measured on an anatomical MRI that match the landmarks or fiducials during an EEG session/run, must be stored separately in the corresponding `*_T1w.json` or `*_T2w.json` file and



should be expressed in voxels (starting from [0, 0, 0]).

Example:

```

1 {
2   "IntendedFor":"/sub-01/ses-01/anat/sub-01_T1w.nii",
3   "EEGCoordinateSystem":"Other",
4   "EEGCoordinateUnits":"mm",
5   "EEGCoordinateSystemDescription":"RAS orientation: Origin halfway between LPA and RPA,
      positive x-axis towards RPA, positive y-axis orthogonal to x-axis through Nasion, z-axis
      orthogonal to xy-plane, pointing in superior direction.",
6   "FiducialsDescription":"Electrodes and fiducials were digitized with Polhemus, fiducials were
      recorded as the centre of vitamin E capsules stuck on the left/right pre-auricular and
      on the nasion, these are also visible on the T1w MRI"
7 }

```

## Landmark photos (\*\_photo.jpg)

Photos of the anatomical landmarks and/or fiducials.

Template:

```

1 sub-<label>/
2   [ses-<label>]/
3   eeg/
4   [sub-<label>[_ses-<label>][_acq-<label>]_photo.jpg]

```

Photos of the anatomical landmarks and/or fiducials are OPTIONAL. Please note that the photos may need to be cropped or blurred to conceal identifying features prior to sharing, depending on the terms of the consent given by the participant.

The `acq` parameter can be used to indicate acquisition of different photos of the same face (or other body part) in different angles to show, for example, the location of the nasion (NAS) as opposed to the right periauricular point (RPA).

Example:

Picture of a NAS fiducial placed between the eyebrows, rather than at the actual anatomical nasion: `sub-0001_ses-001_acq-NAS_photo.jpg`



## Chapter 8

# intracranial Electroencephalography (iEEG)

Support for iEEG was developed as a BIDS Extension Proposal. Please cite the following paper when referring to this part of the standard in context of the academic literature:

Holdgraf, C., Appelhoff, S., Bickel, S., Bouchard, K., D'Ambrosio, S., David, O., Devinsky, O., Dichter, B., Flinker, A., Foster, B. L., Gorgolewski, K. J., Groen, I., Groppe, D., Gunduz, A., Hamilton, L., Honey, C. J., Jas, M., Knight, R., Lauchaux, J.-P., Lau, J. C., Lee-Messer, C., Lundstrom, B. N., Miller, K. J., Ojemann, J. G., Oostenveld, R., Petridou, N., Piantoni, G., Pigorini, A., Pouratian, N., Ramsey, N. F., Stolk, A., Swann, N. C., Tadel, F., Voytek, B., Wandell, B. A., Winawer, J., Whitaker, K., Zehl, L., Hermes, D. (2019). **iEEG-BIDS, extending the Brain Imaging Data Structure specification to human intracranial electrophysiology**. Scientific data, 6. doi: [10.1038/s41597-019-0105-7](https://doi.org/10.1038/s41597-019-0105-7)

## iEEG recording data

Template:

```
1 sub-<label>/
2   [ses-<label>]/
3     ieeg/
4       sub-<label>[_ses-<label>]_task-<task_label>[_run-<index>]_ieeg.<
5         manufacturer_specific_extension>
6       sub-<label>[_ses-<label>]_task-<task_label>[_run-<index>]_ieeg.json
```

The iEEG community uses a variety of formats for storing raw data, and there is no single standard that all researchers agree on. For BIDS, iEEG data **MUST** be stored in one of the following formats:

- [European Data Format \(.edf\)](#)
- [BrainVision Core Data Format \(.vhdr, .eeg, .vmrk\)](#) by Brain Products GmbH
- The format used by the MATLAB toolbox [EEGLAB](#) (.set and .fdt files)
- [Neurodata Without Borders \(.nwb\)](#)
- [MEF3 \(.mef\)](#)

It is **RECOMMENDED** to use the European data format, or the BrainVision data format. It is furthermore discouraged to use the other accepted formats over these **RECOMMENDED** formats, particularly because there are conversion scripts available in most commonly used programming languages to convert data into the **RECOMMENDED** formats.

Future versions of BIDS may extend this list of supported file formats. File formats for future consideration **MUST** have open access documentation, **MUST** have open source implementation for both reading and writing in at least two programming languages and **SHOULD** be widely supported in multiple software packages. Other formats that may be considered in the future should have a clear added advantage over the existing formats and should have wide adoption in the BIDS community.

The data format in which the data was originally stored is especially valuable in case conversion elicits the loss of crucial metadata specific to manufacturers and specific iEEG systems. We also encourage users to provide additional meta information extracted from the manufacturer-specific data files in the sidecar JSON file. Other relevant files MAY be included alongside the original iEEG data in the `/sourcedata` directory.

Note the `RecordingType`, which depends on whether the data stream on disk is interrupted or not. Continuous data is by definition 1 segment without interruption. Epoched data consists of multiple segments that all have the same length (e.g., corresponding to trials) and that have gaps in between. Discontinuous data consists of multiple segments of different length, for example due to a pause in the acquisition.

## Terminology: Electrodes vs. Channels

For proper documentation of iEEG recording metadata it is important to understand the difference between electrode and channel: An iEEG electrode is placed on or in the brain, whereas a channel is the combination of the analog differential amplifier and analog-to-digital converter that result in a potential (voltage) difference that is stored in the iEEG dataset. We employ the following short definitions:

- **Electrode** = A single point of contact between the acquisition system and the recording site (e.g., scalp, neural tissue, ...). Multiple electrodes can be organized as arrays, grids, leads, strips, probes, shafts, caps (for EEG), etc.
- **Channel** = A single analog-to-digital converter in the recording system that regularly samples the value of a transducer, which results in the signal being represented as a time series in the digitized data. This can be connected to two electrodes (to measure the potential difference between them), a magnetic field or magnetic gradient sensor, temperature sensor, accelerometer, etc.

Although the *reference* and *ground* electrodes are often referred to as channels, they are in most common iEEG systems not recorded by themselves. Therefore they are not represented as channels in the data. The type of referencing for all channels and optionally the location of the reference electrode and the location of the ground electrode MAY be specified.

## Sidecar JSON (`*_ieeg.json`)

For consistency between studies and institutions, we encourage users to extract the values of metadata fields from the actual raw data. Whenever possible, please avoid using ad hoc wording.

Generic fields **MUST** be present:

| Field name | Definition  |
|------------|---|
| TaskName   | REQUIRED. Name of the task (for resting state use the "rest" prefix). No two tasks should have the same name. Task label is derived from this field by removing all non alphanumeric ([a-zA-Z0-9]) characters. Note this does not have to be a "behavioral task" that subjects perform, but can reflect some information about the conditions present when the data was acquired (e.g., "rest", "sleep", or "seizure"). |

**SHOULD** be present: For consistency between studies and institutions, we encourage users to extract the values of these fields from the actual raw data. Whenever possible, please avoid using ad hoc wording.

| Field name             | Definition  |
|------------------------|---|
| InstitutionName        | RECOMMENDED. The name of the institution in charge of the equipment that produced the composite instances.    |
| InstitutionAddress     | RECOMMENDED. The address of the institution in charge of the equipment that produced the composite instances. |
| Manufacturer           | RECOMMENDED. Manufacturer of the amplifier system (e.g., "TDT, Blackrock").                                   |
| ManufacturersModelName | RECOMMENDED. Manufacturer's designation of the iEEG amplifier model.  |
| SoftwareVersions       | RECOMMENDED. Manufacturer's designation of the acquisition software.  |
| TaskDescription        | RECOMMENDED. Longer description of the task.  |

| Field name         | Definition  |
|--------------------|---|
| Instructions       | RECOMMENDED. Text of the instructions given to participants before the recording. This is especially important in context of resting state and distinguishing between eyes open and eyes closed paradigms.                    |
| CogAtlasID         | RECOMMENDED. URL of the corresponding <a href="#">Cognitive Atlas Task</a> term.  |
| CogPOID            | RECOMMENDED. URL of the corresponding <a href="#">CogPO</a> term.   |
| DeviceSerialNumber | RECOMMENDED. The serial number of the equipment that produced the composite instances. A pseudonym can also be used to prevent the equipment from being identifiable, as long as each pseudonym is unique within the dataset. |

Specific iEEG fields MUST be present:

| Field name         | Definition  |
|--------------------|---|
| iEEGReference      | REQUIRED. General description of the reference scheme used and (when applicable) of location of the reference electrode in the raw recordings (e.g., "left mastoid", "bipolar", "T01" for electrode with name T01, "intracranial electrode on top of a grid, not included with data", "upside down electrode"). If different channels have a different reference, this field should have a general description and the channel specific reference should be defined in the <code>_channels.tsv</code> file. |
| SamplingFrequency  | REQUIRED. Sampling frequency (in Hz) of all the iEEG channels in the recording (e.g., 2400). All other channels should have frequency specified as well in the <code>channels.tsv</code> file.  |
| PowerLineFrequency | REQUIRED. Frequency (in Hz) of the power grid where the iEEG recording was done (i.e., 50 or 60).   |
| SoftwareFilters    | REQUIRED. A <a href="#">JSON object</a> of temporal software filters applied, or "n/a" if the data is not available. Each key:value pair in the JSON object is a name of the filter and an object in which its parameters are defined as key:value pairs. E.g., <code>{"HighPass": {"HalfAmplitudeCutOffHz": 1, "RollOff": "6dB/Octave"}}</code>  |

Specific iEEG fields SHOULD be present:

| Field name                      | Definition   |
|---------------------------------|--|
| DCOffsetCorrection              | RECOMMENDED. A description of the method (if any) used to correct for a DC offset. If the method used was subtracting the mean value for each channel, use "mean".   |
| HardwareFilters                 | RECOMMENDED. A <a href="#">JSON object</a> of temporal hardware filters applied, or "n/a" if the data is not available. Each key:value pair in the JSON object is a name of the filter and an object in which its parameters are defined as key:value pairs. E.g., <code>{"Highpass RC filter": {"Half amplitude cutoff (Hz)": 0.0159, "Roll-off": "6dB/Octave"}}</code> |
| ElectrodeManufacturer           | RECOMMENDED. can be used if all electrodes are of the same manufacturer (e.g., AD-TECH, DIXI). If electrodes of different manufacturers are used, please use the corresponding table in the <code>_electrodes.tsv</code> file.   |
| ElectrodeManufacturersModelName | RECOMMENDED. If different electrode types are used, please use the corresponding table in the <code>_electrodes.tsv</code> file.   |
| ECOGChannelCount                | RECOMMENDED. Number of iEEG surface channels included in the recording (e.g., 120).  |
| SEEGChannelCount                | RECOMMENDED. Number of iEEG depth channels included in the recording (e.g., 8).  |
| EEGChannelCount                 | RECOMMENDED. Number of scalp EEG channels recorded simultaneously (e.g., 21).  |
| EOGChannelCount                 | RECOMMENDED. Number of EOG channels.   |

| Field name                 | Definition   |
|----------------------------|--|
| ECGChannelCount            | RECOMMENDED. Number of ECG channels.   |
| EMGChannelCount            | RECOMMENDED. Number of EMG channels.   |
| MiscChannelCount           | RECOMMENDED. Number of miscellaneous analog channels for auxiliary signals.  |
| TriggerChannelCount        | RECOMMENDED. Number of channels for digital (TTL bit level) triggers.  |
| RecordingDuration          | RECOMMENDED. Length of the recording in seconds (e.g., 3600).  |
| RecordingType              | RECOMMENDED. Defines whether the recording is "continuous", "discontinuous" or "epoched"; this latter limited to time windows about events of interest (e.g., stimulus presentations, subject responses etc.)  |
| EpochLength                | RECOMMENDED. Duration of individual epochs in seconds (e.g., 1) in case of epoched data. If recording was continuous or discontinuous, leave out the field.  |
| iEEGGround                 | RECOMMENDED. Description of the location of the ground electrode ("placed on right mastoid (M2)").   |
| iEEGPlacementScheme        | RECOMMENDED. Freeform description of the placement of the iEEG electrodes. Left/right/bilateral/depth/surface (e.g., "left frontal grid and bilateral hippocampal depth" or "surface strip and STN depth" or "clinical indication bitemporal, bilateral temporal strips and left grid"). |
| iEEGElectrodeGroups        | RECOMMENDED. Field to describe the way electrodes are grouped into strips, grids or depth probes e.g., {'grid1': "10x8 grid on left temporal pole", 'strip2': "1x8 electrode strip on xxx"}.   |
| SubjectArtefactDescription | RECOMMENDED. Freeform description of the observed subject artefact and its possible cause (e.g., "door open", "nurse walked into room at 2 min", "seizure at 10 min"). If this field is left empty, it will be interpreted as absence of artifacts.                                      |

Specific iEEG fields MAY be present:

| Field name                      | Definition  |
|---------------------------------|---|
| ElectricalStimulation           | OPTIONAL. Boolean field to specify if electrical stimulation was done during the recording (options are "true" or "false"). Parameters for event-like stimulation should be specified in the <code>_events.tsv</code> file (see example below). |
| ElectricalStimulationParameters | OPTIONAL. Free form description of stimulation parameters, such as frequency, shape etc. Specific onsets can be specified in the <code>_events.tsv</code> file. Specific shapes can be described here in freeform text.                         |

Example:

```

1 {
2   "TaskName": "visual",
3   "InstitutionName": "Stanford Hospital and Clinics",
4   "InstitutionAddress": "300 Pasteur Dr, Stanford, CA 94305",
5   "Manufacturer": "Tucker Davis Technologies",
6   "ManufacturersModelName": "n/a",
7   "TaskDescription": "visual gratings and noise patterns",
8   "Instructions": "look at the dot in the center of the screen and press the button when it
   changes color",
9   "iEEGReference": "left mastoid",
10  "SamplingFrequency": 1000,
11  "PowerLineFrequency": 60,
12  "SoftwareFilters": "n/a",
13  "DCOffsetCorrection": 0,

```

```

14  "HardwareFilters":{"Highpass RC filter": {"Half amplitude cutoff (Hz)": 0.0159, "Roll-off": "6
    dB/octave"}},
15  "ElectrodeManufacturer":"AdTech",
16  "ECOGChannelCount":120,
17  "SEEGChannelCount":0,
18  "EEGChannelCount":0,
19  "EOGChannelCount":0,
20  "ECGChannelCount":0,
21  "EMGChannelCount":0,
22  "MiscChannelCount":0,
23  "TriggerChannelCount":0,
24  "RecordingDuration":233.639,
25  "RecordingType":"continuous",
26  "iEEGGround":"placed on the right mastoid",
27  "iEEGPlacementScheme":"right occipital temporal surface",
28  "ElectricalStimulation":false
29 }

```

Note that the date and time information **SHOULD** be stored in the Study key file (`scans.tsv`). As it is indicated there, date time information **MUST** be expressed in the following format `YYYY-MM-DDThh:mm:ss` (ISO8601 date-time format). For example: 2009-06-15T13:45:30. It does not need to be fully detailed, depending on local REB/IRB ethics board policy.

## Channels description (\*\_channels.tsv)

Template:

```

1 sub-<label>/
2   [ses-<label>]/
3   ieeg/
4   [sub-<label>[_ses-<label>]_task-<label>[_run-<index>]]_channels.tsv]

```

A channel represents one time series recorded with the recording system (for example, there can be a bipolar channel, recorded from two electrodes or contact points on the tissue). Although this information can often be extracted from the iEEG recording, listing it in a simple `.tsv` document makes it easy to browse or search (e.g., searching for recordings with a sampling frequency of  $\geq 1000$  Hz). Hence, the `channels.tsv` is RECOMMENDED. The two required columns are channel **name** and **type**. Channels **SHOULD** appear in the table in the same order they do in the iEEG data file. Any number of additional columns may be provided to provide additional information about the channels. Note that electrode positions **SHOULD NOT** be added to this file but to `*_electrodes.tsv`.

The columns of the Channels description table stored in `*_channels.tsv` are:

**MUST** be present:

| Column name | Definition  |
|-------------|---|
| name        | REQUIRED. Label of the channel. The label must correspond to <code>_electrodes.tsv</code> name and all ieeg type channels are required to have a position. The reference channel name <b>MAY</b> be provided in the reference column.         |
| type        | REQUIRED. Type of channel, see below for adequate keywords in this field.   |
| units       | REQUIRED. Physical unit of the value represented in this channel, e.g., V for Volt, specified according to the <a href="#">SI unit symbol</a> and possibly prefix symbol (e.g., mV, V). For guidelines for Units and Prefixes see Appendix V. |
| low_cutoff  | REQUIRED. Frequencies used for the low pass filter applied to the channel in Hz. If no low pass filter was applied, use <b>n/a</b> . Note that anti-alias is a low pass filter, specify its frequencies here if applicable.                   |
| high_cutoff | REQUIRED. Frequencies used for the high pass filter applied to the channel in Hz. If no high pass filter applied, use <b>n/a</b> .  |

**SHOULD** be present:

| Column name        | Definition  |
|--------------------|---|
| reference          | OPTIONAL. Specification of the reference (e.g., 'mastoid', 'ElectrodeName01', 'intracranial', 'CAR', 'other', 'n/a'). If the channel is not an electrode channel (e.g., a microphone channel) use <b>n/a</b> .  |
| group              | OPTIONAL. Which group of channels (grid/strip/seeg/depth) this channel belongs to. This is relevant because one group has one cable-bundle and noise can be shared. This can be a name or number. Note that any groups specified in <code>_electrodes.tsv</code> must match those present here. |
| sampling_frequency | OPTIONAL. Sampling rate of the channel in Hz.   |
| description        | OPTIONAL. Brief free-text description of the channel, or other information of interest (e.g., position (e.g., "left lateral temporal surface", etc.).   |
| notch              | OPTIONAL. Frequencies used for the notch filter applied to the channel, in Hz. If no notch filter applied, use <b>n/a</b> .   |
| status             | OPTIONAL. Data quality observed on the channel (good/bad). A channel is considered bad if its data quality is compromised by excessive noise. Description of noise type <b>SHOULD</b> be provided in <code>[status_description]</code> .  |
| status_description | OPTIONAL. Freeform text description of noise or artifact affecting data quality on the channel. It is meant to explain why the channel was declared bad in <code>[status]</code> .  |

**Example** sub-01\_channels.tsv:

```

1 name type units low_cutoff high_cutoff status status_description
2 LT01 ECG V 300 0.11 good n/a
3 LT02 ECG V 300 0.11 bad broken
4 H01 SEEG V 300 0.11 bad line_noise
5 ECG1 ECG V n/a 0.11 good n/a
6 TR1 TRIG n/a n/a n/a good n/a

```

Restricted keyword list for field type in alphabetic order (shared with the MEG and EEG modality; however, only types that are common in iEEG data are listed here):

| Keyword  | Description  |
|----------|--|
| EEG      | Electrode channel from electroencephalogram                            |
| ECOG     | Electrode channel from electrocorticogram (intracranial)               |
| SEEG     | Electrode channel from stereo-electroencephalogram (intracranial)      |
| DBS      | Electrode channel from deep brain stimulation electrode (intracranial) |
| VEOG     | Vertical EOG (electrooculogram)  |
| HEOG     | Horizontal EOG   |
| EOG      | Generic EOG channel if HEOG or VEOG information not available          |
| ECG      | ElectroCardioGram (heart)  |
| EMG      | ElectroMyoGram (muscle)  |
| TRIG     | System Triggers  |
| AUDIO    | Audio signal   |
| PD       | Photodiode   |
| EYEGAZE  | Eye Tracker gaze   |
| PUPIL    | Eye Tracker pupil diameter   |
| MISC     | Miscellaneous  |
| SYSCLOCK | System time showing elapsed time since trial started                   |
| ADC      | Analog to Digital input  |
| DAC      | Digital to Analog output   |
| REF      | Reference channel  |
| OTHER    | Any other type of channel  |

The free text field for the channel description can for example be specified as intracranial, stimulus, response,

vertical EOG, horizontal EOG, skin conductance, eyetracker, etc.

## Electrode description (\*\_electrodes.tsv)

Template:

```
1 sub-<label>/
2   [ses-<label>]/
3   i EEG/
4   sub-<label>[_ses-<label>][_space-<label>]_electrodes.tsv
```

File that gives the location, size and other properties of iEEG electrodes. Note that coordinates are expected in cartesian coordinates according to the `iEEGCoordinateSystem` and `iEEGCoordinateSystemUnits` fields in `*_coordsystem.json`. If an `*_electrodes.tsv` file is specified, a `*_coordsystem.json` file **MUST** be specified as well.

The optional space label (`*[_space-<label>]_electrodes.tsv`) can be used to indicate the way in which electrode positions are interpreted. The space label needs to be taken from the list in Appendix VIII

For examples:

- `_space-MNI152Lin` (electrodes are coregistered and scaled to a specific MNI template)
- `_space-Talairach` (electrodes are coregistered and scaled to Talairach space)

When referring to the `*_electrodes.tsv` file in a certain *space* as defined above, the `space-<label>` of the accompanying `*_coordsystem.json` **MUST** correspond.

For example:

- `sub-01_space-Talairach_electrodes.tsv`
- `sub-01_space-Talairach_coordsystem.json`

The order of the required columns in the `*_electrodes.tsv` file **MUST** be as listed below.

**MUST** be present:

| Column name | Definition   |
|-------------|--|
| name        | REQUIRED. Name of the electrode contact point.   |
| x           | REQUIRED. X position. The positions of the center of each electrode in xyz space. Units are in millimeters or pixels and are specified in <code>__space-_electrode.json</code> . |
| y           | REQUIRED. Y position.  |
| z           | REQUIRED. Z position. If electrodes are in 2D space this should be a column of n/a values.   |
| size        | REQUIRED. Surface area of the electrode, in $\text{mm}^2$ .  |

**SHOULD** be present:

| Column name  | Definition   |
|--------------|--|
| material     | OPTIONAL. Material of the electrodes.  |
| manufacturer | OPTIONAL. Recommended field to specify the manufacturer for each electrode. Can be used if electrodes were manufactured by more than one company.                                |
| group        | OPTIONAL. Optional field to specify the group that the electrode is a part of. Note that any group specified here should match a group specified in <code>_channels.tsv</code> . |
| hemisphere   | OPTIONAL. Optional field to specify the hemisphere in which the electrode is placed, one of ['L' or 'R'] (use capital).  |

**MAY** be present:

| Column name | Definition  |
|-------------|---|
| type        | OPTIONAL. Optional type of the electrode, e.g., cup, ring, clip-on, wire, needle, ... |



| Column name | Definition  |
|-------------|---|
| impedance   | OPTIONAL. Impedance of the electrode in kOhm.   |
| dimension   | OPTIONAL. Size of the group (grid/strip/probe) that this electrode belongs to. Must be of form [AxB] with the smallest dimension first (e.g., [1x8]). |

Example:

```

1 name  x   y   z   size  manufacturer
2 LT01  19  -39 -16  2.3   Integra
3 LT02  23  -40 -19  2.3   Integra
4 H01   27  -42 -21  5     AdTech

```

## Coordinate System JSON (\*\_coordsystem.json)

Template:

```

1 sub-<label>/
2   [ses-<label>]/
3     ieeg/
4       sub-<label>[_ses-<label>][_space-<label>]_coordsystem.json

```

This `_coordsystem.json` file contains the coordinate system in which electrode positions are expressed. The associated MRI, CT, X-Ray, or operative photo can also be specified.

General fields:

| Field name  | Definition   |
|-------------|--|
| IntendedFor | RECOMMENDED. This can be an MRI/CT or a file containing the operative photo, x-ray or drawing with path relative to the project folder. If only a surface reconstruction is available, this should point to the surface reconstruction file. Note that this file should have the same coordinate system specified in <code>iEEGCoordinateSystem</code> . For example, <b>T1</b> :<br>sub-<label>/ses-<label>/anat/sub-01_T1w.nii.gz <b>Surface</b> : /derivatives/surfaces/sub-<label>/ses-<label>/anat/sub-01_T1w_pial.R.surf.gii<br><b>Operative photo</b> :<br>/sub-<label>/ses-<label>/ieeg/sub-0001_ses-01_acq-photo1_photo.jpg<br><b>Talairach</b> : /derivatives/surfaces/sub-Talairach/ses-01/anat/sub-Talairach_T1w_pial.R.surf.gii |

Fields relating to the iEEG electrode positions:

| Field name                      | Definition  |
|---------------------------------|---|
| iEEGCoordinateSystem            | REQUIRED. Defines the coordinate system for the iEEG electrodes. See Appendix VIII for a list of restricted keywords. If positions correspond to pixel indices in a 2D image (of either a volume-rendering, surface-rendering, operative photo, or operative drawing), this must be "Pixels". For more information, see the section on 2D coordinate systems      |
| iEEGCoordinateUnits             | REQUIRED. Units of the <code>_electrodes.tsv</code> , MUST be "m", "mm", "cm" or "pixels".  |
| iEEGCoordinateSystemDescription | RECOMMENDED. Freeform text description or link to document describing the iEEG coordinate system system in detail (e.g., "Coordinate system with the origin at anterior commissure (AC), negative y-axis going through the posterior commissure (PC), z-axis going to a mid-hemispheric point which lies superior to the AC-PC line, x-axis going to the right"). |

| Field name                          | Definition  |
|-------------------------------------|---|
| iEEGCoordinateProcessingDescription | RECOMMENDED. Has any post-processing (such as projection) been done on the electrode positions (e.g., "surface_projection", "none").                      |
| iEEGCoordinateProcessingReference   | RECOMMENDED. A reference to a paper that defines in more detail the method used to localize the electrodes and to post-process the electrode positions. . |

## Recommended 3D coordinate systems

It is preferred that electrodes are localized in a 3D coordinate system (with respect to a pre- and/or post-operative anatomical MRI or CT scans or in a standard space as specified in BIDS Appendix VIII about preferred names of coordinate systems, such as ACPC).

## Allowed 2D coordinate systems

If electrodes are localized in 2D space (only x and y are specified and z is n/a), then the positions in this file must correspond to the locations expressed in pixels on the photo/drawing/rendering of the electrodes on the brain. In this case, coordinates must be (row,column) pairs, with (0,0) corresponding to the upper left pixel and (N,0) corresponding to the lower left pixel.

## Multiple coordinate systems

If electrode positions are known in multiple coordinate systems (e.g., MRI, CT and MNI), these spaces can be distinguished by the optional [\_space-<label>] field, see the [\\*\\_electrodes.tsv-section](#) for more information. Note that the [\_space-<label>] fields must correspond between \*\_electrodes.tsv and \*\_coordsystem.json if they refer to the same data.

Example:

```

1 {
2   "IntendedFor": "/sub-01/ses-01/anat/sub-01_T1w.nii.gz",
3   "iEEGCoordinateSystem": "ACPC",
4   "iEEGCoordinateUnits": "mm",
5   "iEEGCoordinateSystemDescription": "Coordinate system with the origin at anterior commissure
      (AC), negative y-axis going through the posterior commissure (PC), z-axis going to a
      mid-hemispheric point which lies superior to the AC-PC line, x-axis going to the right",
6   "iEEGCoordinateProcessingDescription": "surface_projection",
7   "iEEGCoordinateProcessingReference": "Hermes et al., 2010 JNeuroMeth"
8 }
```

## Photos of the electrode positions (\*\_photo.jpg)

Template:

```

1 sub-<label>/
2   [ses-<label>]/
3     ieeg/
4       sub-<label>[_ses-<label>][_acq-<label>]_photo.json
```

These can include photos of the electrodes on the brain surface, photos of anatomical features or landmarks (such as sulcal structure), and fiducials. Photos can also include an X-ray picture, a flatbed scan of a schematic drawing made during surgery, or screenshots of a brain rendering with electrode positions. The photos may need to be cropped and/or blurred to conceal identifying features or entirely omitted prior to sharing, depending on obtained consent.

If there are photos of the electrodes, the acquisition field should be specified with:


- \*\_photo.jpg in case of an operative photo
- \*\_acq-xray#\_photo.jpg in case of an x-ray picture

- \*\_acq-drawing#\_photo.jpg in case of a drawing or sketch of electrode placements
- \*\_acq-render#\_photo.jpg in case of a rendering

The session label may be used to specify when the photo was taken.


Example of the operative photo of ECoG electrodes (here is an annotated example in which electrodes and vasculature are marked, taken from Hermes et al., JNeuroMeth 2010).

```
1 sub-0001_ses-01_acq-photo1_photo.jpg
2 sub-0001_ses-01_acq-photo2_photo.jpg
```

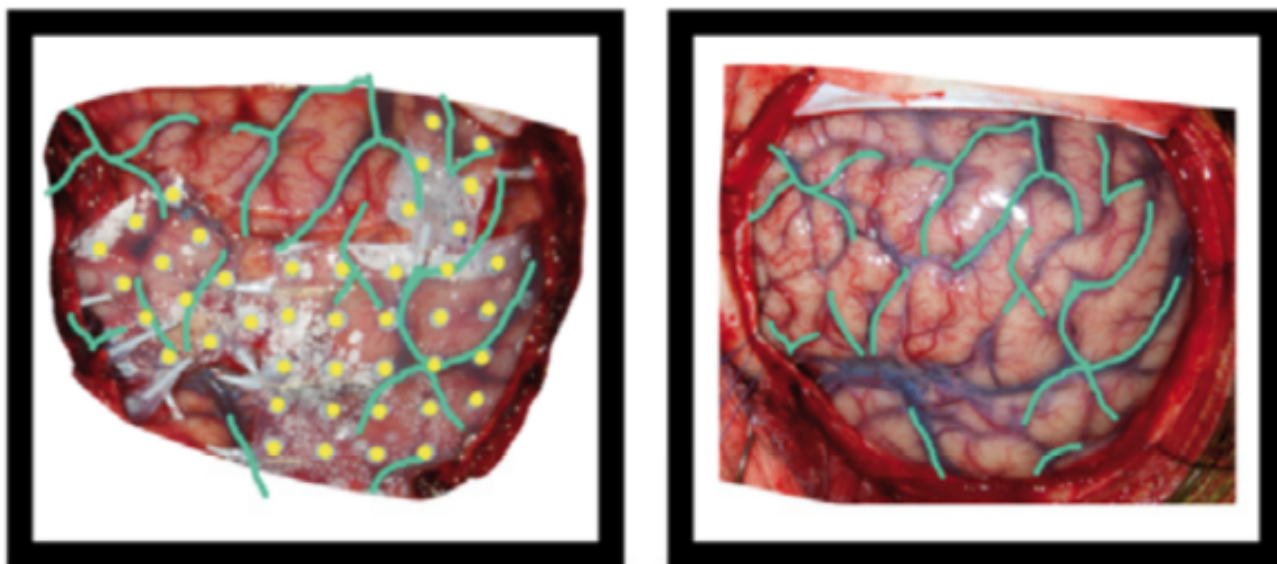

Below is an example of a volume rendering of the cortical surface with a superimposed subdural electrode implantation. This map is often provided by the EEG technician and provided to the epileptologists (e.g., see Burneo JG et al. 2014 <https://doi.org/10.1016/j.clineuro.2014.03.020>).

```
1 sub-0002_ses-01_acq-render_photo.jpg
```

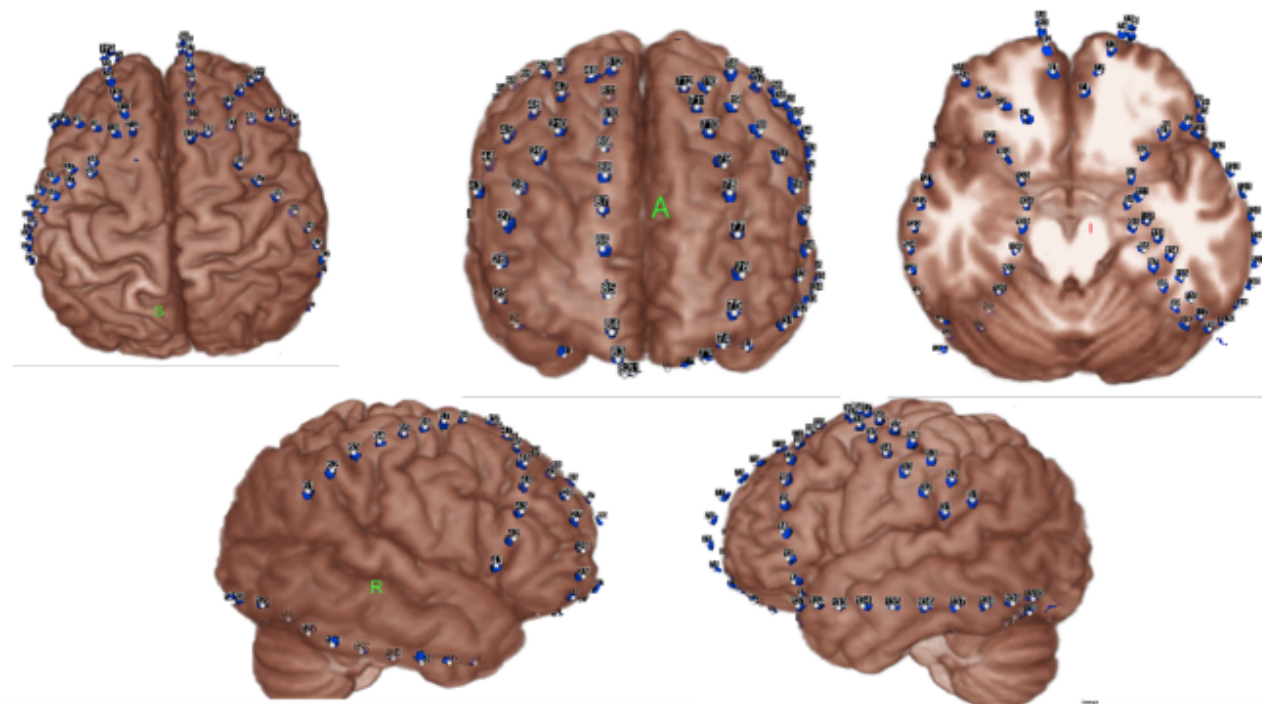

## Electrical stimulation

In case of electrical stimulation of brain tissue by passing current through the iEEG electrodes, and the electrical stimulation has an event structure (on-off, onset, duration), the `_events.tsv` file can contain the electrical stimulation parameters in addition to other events. Note that these can be intermixed with other task events. Electrical stimulation parameters can be described in columns called `electrical_stimulation_<label>`, with labels chosen by the researcher and optionally defined in more detail in an accompanying `_events.json` file (as per the main BIDS spec). Functions for complex stimulation patterns can, similar as when a video is presented, be stored in a folder in the `/stimuli/` folder. For example: `/stimuli/electrical_stimulation_functions/biphasic.tsv`

Example:

|   | onset | duration | trial_type                     | electrical_stimulation_type | electrical_stimulation_site |
|---|-------|----------|--------------------------------|-----------------------------|-----------------------------|
|   |       |          | electrical_stimulation_current |                             |                             |
| 1 | 1.2   | 0.001    | electrical_stimulation         | biphasic                    | LT01-LT02                   |
| 2 |       | 0.005    |                                |                             |                             |
| 3 | 1.3   | 0.001    | electrical_stimulation         | biphasic                    | LT01-LT02                   |
|   |       | 0.005    |                                |                             |                             |
| 4 | 2.2   | 0.001    | electrical_stimulation         | biphasic                    | LT02-LT03                   |
|   |       | 0.005    |                                |                             |                             |
| 5 | 4.2   | 1        | electrical_stimulation         | complex                     | LT02-LT03                   |
|   | a     |          |                                |                             | n/                          |
| 6 | 15.2  | 3        | auditory_stimulus              | n/a                         | n/a                         |
|   | a     |          |                                |                             | n/                          |

## Chapter 9

# Task events

Template:

```
1 sub-<label>/[ses-<label>]
2   func/
3     <matches>_events.tsv
4     <matches>_events.json
```

Where `<matches>` corresponds to task file name. For example: `sub-control01_task-nback`. It is also possible to have a single `_events.tsv` file describing events for all participants and runs (see Inheritance Principle). As with all other tabular data, `_events` files may be accompanied by a JSON file describing the columns in detail (see here).

The purpose of this file is to describe timing and other properties of events recorded during the scan. Events MAY be either stimuli presented to the participant or participant responses. A single event file MAY include any combination of stimuli and response events. Events MAY overlap in time. Please mind that this does not imply that only so called "event related" study designs are supported (in contrast to "block" designs) - each "block of events" can be represented by an individual row in the `_events.tsv` file (with a long duration). Each task events file **REQUIRES** a corresponding task imaging data file (but a single events file MAY be shared by multiple imaging data files - see Inheritance principle). The tabular files consists of one row per event and a set of **REQUIRED** and **OPTIONAL** columns:

| Column name   | Description   |
|---------------|---|
| onset         | REQUIRED. Onset (in seconds) of the event measured from the beginning of the acquisition of the first volume in the corresponding task imaging data file. If any acquired scans have been discarded before forming the imaging data file, ensure that a time of 0 corresponds to the first image stored. In other words negative numbers in "onset" are allowed <sup>5</sup> .  |
| duration      | REQUIRED. Duration of the event (measured from onset) in seconds. Must always be either zero or positive. A "duration" value of zero implies that the delta function or event is so short as to be effectively modeled as an impulse.   |
| sample        | OPTIONAL. Onset of the event according to the sampling scheme of the recorded modality (i.e., referring to the raw data file that the <code>events.tsv</code> file accompanies).  |
| trial_type    | OPTIONAL. Primary categorisation of each trial to identify them as instances of the experimental conditions. For example: for a response inhibition task, it could take on values "go" and "no-go" to refer to response initiation and response inhibition experimental conditions.   |
| response_time | OPTIONAL. Response time measured in seconds. A negative response time can be used to represent preemptive responses and "n/a" denotes a missed response.  |
| stim_file     | OPTIONAL. Represents the location of the stimulus file (image, video, sound etc.) presented at the given onset time. There are no restrictions on the file formats of the stimuli files, but they should be stored in the <code>/stimuli</code> folder (under the root folder of the dataset; with optional subfolders). The values under the <code>stim_file</code> column correspond to a path relative to <code>/stimuli</code> . For example <code>"images/cat03.jpg"</code> will be translated to <code>/stimuli/images/cat03.jpg</code> . |
| value         | OPTIONAL. Marker value associated with the event (e.g., the value of a TTL trigger that was recorded at the onset of the event).  |

| Column name | Description  |
|-------------|--|
| HED         | OPTIONAL. Hierarchical Event Descriptor (HED) Tag. See Appendix III for details. |

5 For example in case there is an in scanner training phase that begins before the scanning sequence has started events from this sequence should have negative onset time counting down to the beginning of the acquisition of the first volume.

An arbitrary number of additional columns can be added. Those allow describing other properties of events that could be later referred in modelling and hypothesis extensions of BIDS.

In case of multi-echo task run, a single `_events.tsv` file will suffice for all echoes.

Example:

```
1 sub-control01/
2   func/
3     sub-control01_task-stopsignal_events.tsv

1 onset duration trial_type response_time stim_file
2 1.2 0.6 go 1.435 images/red_square.jpg
3 5.6 0.6 stop 1.739 images/blue_square.jpg
```

References to existing databases can also be encoded using additional columns. Example 2 includes references to the Karolinska Directed Emotional Faces (KDEF) database<sup>6</sup>:

<sup>6</sup><http://www.emotionlab.se/resources/kdef>

Example:

```
1 sub-control01/
2   func/
3     sub-control01_task-emoface_events.tsv

1 onset duration trial_type identifier database response_time
2 1.2 0.6 afraid AF01AFAP kdef 1.435
3 5.6 0.6 angry AM01AFAN kdef 1.739
4 5.6 0.6 sad AF01ANSA kdef 1.739
```

For multi-echo files `events.tsv` file is applicable to all echos of particular run:

```
1 sub-01_task-cuedSGT_run-1_events.tsv
2 sub-01_task-cuedSGT_run-1_echo-1_bold.nii.gz
3 sub-01_task-cuedSGT_run-1_echo-2_bold.nii.gz
4 sub-01_task-cuedSGT_run-1_echo-3_bold.nii.gz
```

## Chapter 10

# Physiological and other continuous recordings

Template:

```
1 sub-<label>/[ses-<label>/]
2   func/
3     <matches>[_recording-<label>]_physio.tsv.gz
4     <matches>[_recording-<label>]_physio.json
5     <matches>[_recording-<label>]_stim.tsv.gz
6     <matches>[_recording-<label>]_stim.json
```

Optional: Yes

Where `<matches>` corresponds to task file name without the `_bold.nii[.gz]` suffix. For example: `sub-control01_task-nback_run-1`. If the same continuous recording has been used for all subjects (for example in the case where they all watched the same movie) one file can be used and placed in the root directory. For example: `task-movie_stim.tsv.gz`

Physiological recordings such as cardiac and respiratory signals and other continuous measures (such as parameters of a film or audio stimuli) can be specified using two files: a gzip compressed TSV file with data (without header line) and a JSON for storing the following metadata fields:

| Field name        | Definition  |
|-------------------|---|
| SamplingFrequency | REQUIRED. Sampling frequency in Hz of all columns in the file.  |
| StartTime         | REQUIRED. Start time in seconds in relation to the start of acquisition of the first data sample in the corresponding neural dataset (negative values are allowed). |
| Columns           | REQUIRED. Names of columns in file.   |

Additional metadata may be included as in any TSV file to specify, for example, the units of the recorded time series. Please note that in contrast to other TSV files this one does not include a header line. Instead the name of columns are specified in the JSON file. This is to improve compatibility with existing software (FSL PNM) as well as make support for other file formats possible in the future. Recordings with different sampling frequencies and/or starting times should be stored in separate files. The following naming conventions should be used for column names:

| Column name | Definition   |
|-------------|--|
| cardiac     | continuous pulse measurement                         |
| respiratory | continuous breathing measurement                     |
| trigger     | continuous measurement of the scanner trigger signal |

Any combination of those three can be included as well as any other stimuli related continuous variables (such as low level image properties in a video watching paradigm).

Physiological recordings (including eye tracking) should use the `_physio` suffix, and signals related to the stimulus should use `_stim` suffix. For motion parameters acquired from scanner side motion correction please use `_physio` suffix.

More than one continuous recording file can be included (with different sampling frequencies). In such case use different labels. For example: `_recording-contrast`, `_recording-saturation`. The full file name could then look like this: `sub-control01_task-nback_run-2_recording-movie_stim.tsv.gz`

For multi-echo data, `physio.tsv` file is applicable to all echos of particular run. For eg:

```
1 sub-01_task-cuedSGT_run-1_physio.tsv.gz
2 sub-01_task-cuedSGT_run-1_echo-1_bold.nii.gz
3 sub-01_task-cuedSGT_run-1_echo-2_bold.nii.gz
4 sub-01_task-cuedSGT_run-1_echo-3_bold.nii.gz
```

Example:

```
1 sub-control01/
2   func/
3     sub-control01_task-nback_physio.tsv.gz
```

(after decompression)

```
1 34 110 0
2 44 112 0
3 23 100 1
```

```
1 sub-control01/
2   func/
3     sub-control01_task-nback_physio.json
```

```
1 {
2   "SamplingFrequency": 100.0,
3   "StartTime": -22.345,
4   "Columns": ["cardiac", "respiratory", "trigger"],
5   "cardiac": {
6     "Units": "mV"
7   }
8 }
```



## Chapter 11

# Behavioral experiments (with no MRI)

Template:

```
1 sub-<label>/[ses-<label>/]
2   beh/
3     sub-<label>[_ses-<label>]_task-<task_name>_events.tsv
4     sub-<label>[_ses-<label>]_task-<task_name>_events.json
5     sub-<label>[_ses-<label>]_task-<task_name>_beh.tsv
6     sub-<label>[_ses-<label>]_task-<task_name>_beh.json
7     sub-<label>[_ses-<label>]_task-<task_name>_physio.tsv.gz
8     sub-<label>[_ses-<label>]_task-<task_name>_physio.json
9     sub-<label>[_ses-<label>]_task-<task_name>_stim.tsv.gz
10    sub-<label>[_ses-<label>]_task-<task_name>_stim.json
```

In addition to logs from behavioral experiments performed along imaging data acquisitions one can also include data from experiments performed outside of the scanner. The results of those experiments can be stored in the `beh` folder using the same formats for event timing (`_events.tsv`), metadata (`_events.json`), physiological (`_physio.tsv.gz`, `_physio.json`) and other continuous recordings (`_stim.tsv.gz`, `_stim.json`) as for tasks performed during MRI acquisitions. Additionally, events files that do not include the mandatory `onset` and `duration` columns can still be included, but should be labelled `_beh.tsv` rather than `_events.tsv`.

## Chapter 12

# Longitudinal studies with multiple sessions (visits)

Multiple sessions (visits) are encoded by adding an extra layer of directories and file names in the form of **ses-*<label>***. Session label can consist only of alphanumeric characters [a-zA-Z0-9] and should be consistent across subjects. If numbers are used in session labels we recommend using zero padding (for example **ses-01**, **ses-11** instead of **ses-1**, **ses-11**). This makes results of alphabetical sorting more intuitive. Acquisition time of session can be defined in the sessions file (see below for details).

The extra session layer (at least one **/ses-*<label>*** subfolder) should be added for all subjects if at least one subject in the dataset has more than one session. Skipping the session layer for only some subjects in the dataset is not allowed. If a **/ses-*<label>*** subfolder is included as part of the directory hierarchy, then the same **ses-*<label>*** tag must also be included as part of the file names themselves.

```
1 sub-control01/
2   ses-predrug/
3     anat/
4       sub-control01_ses-predrug_T1w.nii.gz
5       sub-control01_ses-predrug_T1w.json
6       sub-control01_ses-predrug_T2w.nii.gz
7       sub-control01_ses-predrug_T2w.json
8     func/
9       sub-control01_ses-predrug_task-nback_bold.nii.gz
10      sub-control01_ses-predrug_task-nback_bold.json
11      sub-control01_ses-predrug_task-nback_events.tsv
12      sub-control01_ses-predrug_task-nback_cont-physio.tsv.gz
13      sub-control01_ses-predrug_task-nback_cont-physio.json
14      sub-control01_ses-predrug_task-nback_sbref.nii.gz
15     dwi/
16       sub-control01_ses-predrug_dwi.nii.gz
17       sub-control01_ses-predrug_dwi.bval
18       sub-control01_ses-predrug_dwi.bvec
19     fmap/
20       sub-control01_ses-predrug_phasediff.nii.gz
21       sub-control01_ses-predrug_phasediff.json
22       sub-control01_ses-predrug_magnitude1.nii.gz
23     sub-control01_ses-predrug_scans.tsv
24   ses-postdrug/
25     func/
26       sub-control01_ses-postdrug_task-nback_bold.nii.gz
27       sub-control01_ses-postdrug_task-nback_bold.json
28       sub-control01_ses-postdrug_task-nback_events.tsv
29       sub-control01_ses-postdrug_task-nback_cont-physio.tsv.gz
30       sub-control01_ses-postdrug_task-nback_cont-physio.json
31       sub-control01_ses-postdrug_task-nback_sbref.nii.gz
32     fmap/
```

```

33         sub-control01_ses-postdrug_phasediff.nii.gz
34         sub-control01_ses-postdrug_phasediff.json
35         sub-control01_ses-postdrug_magnitude1.nii.gz
36     sub-control01_ses-postdrug_scans.tsv
37     sub-control01_sessions.tsv
38 participants.tsv
39 dataset_description.json
40 README
41 CHANGES

```

## Sessions file

Template:

```

1 sub-<label>/
2     sub-<label>_sessions.tsv

```

Optional: Yes

In case of multiple sessions there is an option of adding additional participant key files describing variables changing between sessions. In such case one file per participant should be added. These files need to include compulsory `session_id` column and describe each session by one and only one row. Column names in per participant key files have to be different from group level participant key column names.

`_sessions.tsv` example:

```

1 session_id  acq_time  systolic_blood_pressure
2 ses-predrug 2009-06-15T13:45:30 120
3 ses-postdrug 2009-06-16T13:45:30 100
4 ses-followup 2009-06-17T13:45:30 110

```

## Multi-site or multi-center studies

This version of the BIDS specification does not explicitly cover studies with data coming from multiple sites or multiple centers (such extension is planned in BIDS 2.0.0). There are however ways to model your data without any loss in terms of metadata.

### Treat each site/center as a separate dataset

The simplest way of dealing with multiple sites is to treat data from each site as a separate and independent BIDS dataset with a separate `participants.tsv` and other metadata files. This way you can feed each dataset individually to BIDS Apps and everything should just work.

### Option 2: Combining sites/centers into one dataset

Alternatively you can combine data from all sites into one dataset. To identify which site each subjects comes from you can add a `site` column in the `participants.tsv` file indicating the source site. This solution allows you to analyze all of the subjects together in one dataset. One caveat is that subjects from all sites will have to have unique labels. To enforce that and improve readability you can use a subject label prefix identifying the site. For example `sub-NUY001`, `sub-MIT002`, `sub-MPG002` etc. Remember that hyphens and underscores are not allowed in subject labels.

## Chapter 13

# BIDS Extension Proposals

The BIDS specification can be extended in a backwards compatible way and will evolve over time. This is accomplished with BIDS Extension Proposals (BEPs), which are community-driven processes.

On the [BIDS homepage](#) you can find a list of extension protocols that are currently being worked on.

## Chapter 14

# Appendix I: Contributors

Legend (source: <https://github.com/kentcdodds/all-contributors>)

| Emoji | Represents  |
|-------|---|
|       | Answering Questions (on the mailing list, NeuroStars, GitHub, or in person) |
|       | Bug reports   |
|       | Blogposts   |
|       | Code  |
|       | Documentation and specification   |
|       | Design  |
|       | Examples  |
|       | Event Organizers  |
|       | Financial Support   |
|       | Funding/Grant Finders   |
|       | Ideas & Planning  |
|       | Infrastructure (Hosting, Build-Tools, etc)                                  |
|       | Plugin/utility libraries  |
|       | Reviewed Pull Requests  |
|       | Tools   |
|       | Translation   |
|       | Tests   |
|       | Tutorials   |
|       | Talks   |
|       | Videos  |

The following individuals have contributed to the Brain Imaging Data Structure ecosystem (in alphabetical order). If you contributed to the BIDS ecosystem and your name is not listed, please add it.

- Fidel Alfaro Almagro
- Stefan Appelhoff
- Tibor Auer
- Sylvain Baillet
- Stephan Bickel
- Elizabeth Bock
- Kristofer Bouchard
- Eric Bridgeford
- Teon L. Brooks
- Suyash Bhogawar
- Vince D. Calhoun
- Alexander L. Cohen
- R. Cameron Craddock
- Sasha D'Ambrosio
- Samir Das
- Olivier David

- Orrin Devinsky
- Alejandro de la Vega
- Arnaud Delorme
- Benjamin Dichter
- Eugene P. Duff
- Elizabeth DuPre
- Eric A. Earl
- Anders Eklund
- Oscar Esteban
- Franklin W. Feingold
- Guillaume Flandin
- Adeen Flinker
- Brett L. Foster
- Remi Gau
- Satrajit S. Ghosh
- Tristan Glatard
- Mathias Goncalves
- Krzysztof J. Gorgolewski
- Alexandre Gramfort
- Jeffrey S. Grethe
- Iris Groen
- David Groppe
- Aysegul Gunduz
- Yaroslav O. Halchenko
- Liberty Hamilton
- Daniel A. Handwerker
- Michael Hanke
- Michael P. Harms
- Richard N. Henson
- Peer Herholz
- Dora Hermes
- Katja Heuer
- Chris Holdgraf
- Christopher J. Honey
- Jean-Christophe Houde
- International Neuroinformatics Coordinating Facility
- Andrew Janke
- Mainak Jas
- David Keator
- James Kent
- Gregory Kiar
- Robert Knight
- Jean-Philippe Lachaux
- Pamela LaMontagne
- Kevin Larcher
- Jonathan C. Lau
- Laura and John Arnold Foundation
- Christopher Lee-Messer
- Jon Haitz Legarreta
- Xiangrui Li
- Vladimir Litvak
- Brian N. Lundstrom
- Dan Lurie
- Camille Maumet
- Christopher J. Markiewicz
- Kai J. Miller
- Jeremy Moreau
- Zachary Michael
- Ezequiel Mikulan
- Michael P. Milham
- Henk Mutsaerts

- National Institute of Mental Health
- Mikael Naveau
- B. Nolan Nichols
- Thomas E. Nichols
- Dylan Nielson
- Guiomar Niso
- Michael P. Notter
- Jeffrey G. Ojemann
- Robert Oostenveld
- Dimitri Papadopoulos Orfanos
- Patrick Park
- Dianne Patterson
- John Pellman
- Cyril Pernet
- Franco Pestilli
- Natalia Petridou
- Dmitry Petrov
- Christophe Phillips
- Gio Piantoni
- Andrea Pigorini
- Russell A. Poldrack
- Jean-Baptiste Poline
- Wouter V. Potters
- Nader Pouratian
- Pradeep Reddy Raamana
- Vasudev Raguram
- Nick F. Ramsey
- Alex Rockhill
- Ariel Rokem
- Matt Sanderson
- Gunnar Schaefer
- Jan-Mathijs Schoffelen
- Vanessa Sochat
- Arjen Stolk
- Nicole C. Swann
- François Tadel
- Roberto Toro
- Sébastien Tourbier
- William Triplett
- Jessica A. Turner
- Bradley Voytek
- Brian A. Wandell
- Joseph Wexler
- Kirstie Whitaker
- Jonathan Winawer
- Gaël Varoquaux
- Tal Yarkoni
- Lyuba Zehl
- Sébastien Tourbier

## Chapter 15

# Appendix II: Licenses

This section lists a number of common licenses for datasets and defines suggested abbreviations for use in the dataset metadata specifications.

Please note that this list only serves to provide some examples for possible licenses. The terms of any license should be consistent with the informed consent obtained from participants and any institutional limitations on distribution.

| Identifier | License name   | Description   |
|------------|--|---|
| PD         | Public Domain  | No license required for any purpose; the work is not subject to copyright in any jurisdiction.  |
| PDDL       | Open Data Commons<br>Public Domain<br>Dedication and License | License to assign public domain like permissions without giving up the copyright: <a href="http://opendatacommons.org/licenses/pddl/">http://opendatacommons.org/licenses/pddl/</a>   |
| CC0        | Creative Commons Zero<br>1.0 Universal.                      | Use this if you are a holder of copyright or database rights, and you wish to waive all your interests in your work worldwide:<br><a href="http://opendatacommons.org/licenses/cc0/">http://opendatacommons.org/licenses/cc0/</a> |



## Chapter 16

# Appendix III: Hierarchical Event Descriptor (HED) Tags

HED is a controlled vocabulary of terms describing events in a behavioral paradigm. HED was originally developed with EEG in mind, but is applicable to all behavioral experiments. Each level of a hierarchical tag is delimited with a forward slash (/). An HED string contains one or more HED tags separated by commas (,). Parentheses (brackets, ()) group tags and enable specification of multiple items and their attributes in a single HED string (see section 2.4 in [HED Tagging Strategy Guide](#)). For more information about HED and tools available to validate and match HED strings, please visit [www.hedtags.org](http://www.hedtags.org). Since dedicated fields already exist for the overall task classification in the sidecar JSON files (**CogAtlasID** and **CogPOID**), HED tags from the **Paradigm** HED subcategory should not be used to annotate events.

There are several ways to associate HED annotations with events within the BIDS framework. The most direct way is to use the HED column of the `__events.tsv` file to annotate events:

Example:

```
1 onset    duration    HED
2 1.1 n/a Event/Category/Experimental stimulus, Event/Label/CrossFix, Sensory presentation/Visual
   , Item/Object/2D Shape/Cross
3 1.3 n/a Event/Category/Participant response, Event/Label/ButtonPress, Action/Button press
4 ...
```

The direct approach requires that each line in the events file must be annotated. Since there are typically thousands of events in each experiment, this method of annotation is usually not convenient unless the annotations are automatically generated. In many experiments, the event instances fall into a much smaller number of categories, and often these categories are labeled with numerical codes or short names. It is therefore more convenient to associate the HED annotations with these categories and allow the analysis tools to make the association with individual event instances during analysis. To use this approach, your `__events.tsv` file should associate a category (often called an event code) with each event instance. Since BIDS allows an arbitrary number of columns to be included in an `__events.tsv` file, you can make this association by including columns representing various types of event categories in your `__events.tsv` file.

Example:

```
1 onset    duration    mycodes
2 1.1 n/a Fixation
3 1.3 n/a Button
4 1.8 n/a Target
5 ...
```

If you provide an `__events.json` file somewhere in your data hierarchy that has an HED mapping for **mycodes**, the HED tags associated with a given **mycodes** value can then be associated with the event instances in that category. You may provide a HED column and multiple category columns. The union of the relevant HED tags will then be associated with the event instance.

Example:

```

1 {
2   "mycodes": {
3     "HED": {
4       "Fixation": "Event/Category/Experimental stimulus, Event/Label/CrossFix, Event/
        Description/A cross appears at screen center to serve as a fixation point,
        Sensory presentation/Visual, Item/Object/2D Shape/Cross, Attribute/Visual/
        Fixation point, Attribute/Visual/Rendering type/Screen, Attribute/Location/
        Screen/Center",
5       "Target": "Event/Label/Target image, Event/Description/A white airplane as the RSVP
        target superimposed on a satellite image is displayed., Event/Category/
        Experimental stimulus, (Item/Object/Vehicle/Aircraft/Airplane, Participant/
        Effect/Cognitive/Target, Sensory presentation/Visual/Rendering type/Screen/2D),
        (Item/Natural scene/Arial/Satellite, Sensory presentation/Visual/Rendering type/
        Screen/2D)",
6       "Button": "...",
7     }
8   }
9 }

```

The tags in the HED column are often specific to the event instances, while the common properties associated with categories such as `mycodes` are encapsulated in the `__events.json` dictionary. Downstream tools should not distinguish between tags specified using the different mechanisms. Further, the normal BIDS inheritance principle applies so these data dictionaries can appear higher in the BIDS hierarchy.

# Chapter 17

## Appendix IV: Entity table

This section compiles the entities (key-value pairs) described throughout this specification, and establishes a common order within a filename. For example, if a file has an acquisition and reconstruction label, the acquisition entity must precede the reconstruction entity. REQUIRED and OPTIONAL entities for a given file type are denoted. Entity formats indicate whether the value is alphanumeric (<label>) or numeric (<index>).

A general introduction to entities is given in the section on file name structure

| Entity   | Subject  | Session  | Task     | Acquisition | Contrast Enhancing Agent | Reconstruction | Phase-Encoding Direction | Run      | Corresponding modality | Echo     | Recording | Processed (on device) | Space    |
|--|----------|----------|----------|-------------|--------------------------|----------------|--------------------------|----------|------------------------|----------|-----------|-----------------------|----------|
| Format   | sub      | ses      | task     | acq         | ce-<label>               | rec            | dir                      | run      | mod                    | echo     | record    | proc                  | space    |
|  | -<label> | -<label> | -<label> | -<label>    | -<label>                 | -<label>       | -<label>                 | -<label> | -<label>               | -<label> | -<label>  | -<label>              | -<label> |
|  | >        | >        | >        | >           | >                        | >              | >                        | >        | >                      | >        | >         | >                     | >        |
| anat(T1w T2w T1rho T1map T2map T2star FLAIR FLASH PD PDmap PDT2 inplaneT1 inplaneT2 angio) | REQUIRED | REQUIRED | OPTIONAL | OPTIONAL    | OPTIONAL                 | OPTIONAL       |                          |          |                        |          |           |                       |          |
| anat(defacemask)   | REQUIRED | OPTIONAL | OPTIONAL | OPTIONAL    | OPTIONAL                 |                |                          |          | OPTIONAL               |          |           |                       |          |
| func(bold cbv phase sbref events)  | REQUIRED | REQUIRED | REQUIRED | OPTIONAL    | OPTIONAL                 | OPTIONAL       | OPTIONAL                 | OPTIONAL |                        | OPTIONAL |           |                       |          |
| func(physio stim)  | REQUIRED | REQUIRED | REQUIRED | OPTIONAL    | OPTIONAL                 | OPTIONAL       | OPTIONAL                 | OPTIONAL |                        | OPTIONAL |           |                       |          |
| dwi(dwi bvec bval)   | REQUIRED | REQUIRED | OPTIONAL | OPTIONAL    | OPTIONAL                 | OPTIONAL       | OPTIONAL                 | OPTIONAL |                        |          |           |                       |          |
| fmap(phasediff phase1 phase2 magnitude1 magnitude2 magnitude fieldmap)                     | REQUIRED | REQUIRED | OPTIONAL | OPTIONAL    | OPTIONAL                 |                |                          | OPTIONAL |                        |          |           |                       |          |
| fmap(epi)  | REQUIRED | REQUIRED | OPTIONAL | OPTIONAL    | OPTIONAL                 | REQUIRED       | OPTIONAL                 |          |                        |          |           |                       |          |
| beh(events stim physio)  | REQUIRED | REQUIRED | REQUIRED |             |                          |                |                          |          |                        |          |           |                       |          |
| meg  | REQUIRED | REQUIRED | REQUIRED | OPTIONAL    |                          |                |                          | OPTIONAL |                        |          |           | OPTIONAL              |          |
| eeg  | REQUIRED | REQUIRED | REQUIRED | OPTIONAL    |                          |                |                          | OPTIONAL |                        |          |           |                       |          |
| ieeg   | REQUIRED | REQUIRED | REQUIRED | OPTIONAL    |                          |                |                          | OPTIONAL |                        |          |           |                       |          |
| channels(meg/eeg/ieeg)   | REQUIRED | REQUIRED | REQUIRED |             |                          |                |                          | OPTIONAL |                        |          |           |                       |          |
| headshape(meg)   | REQUIRED | REQUIRED | OPTIONAL | OPTIONAL    |                          |                |                          |          |                        |          |           |                       | OPTIONAL |
| markers(meg)   | REQUIRED | REQUIRED | OPTIONAL | OPTIONAL    |                          |                |                          |          |                        |          |           |                       | OPTIONAL |
| photo(meg/eeg/ieeg)  | REQUIRED | REQUIRED | OPTIONAL | OPTIONAL    |                          |                |                          |          |                        |          |           |                       |          |
| electrodes(eeg/ieeg)   | REQUIRED | REQUIRED | OPTIONAL | OPTIONAL    |                          |                |                          |          |                        |          |           |                       | OPTIONAL |
| events(meg/eeg/ieeg)   | REQUIRED | REQUIRED | REQUIRED |             |                          |                |                          | OPTIONAL |                        |          |           |                       |          |

## Chapter 18

# Appendix V: Units

Following the International System of Units (SI, abbreviated from the French *Système international (d'unités)*)

| Unit name      | Unit symbol        | Quantity name                              |
|----------------|--------------------|--|
| metre          | m                  | length                                     |
| kilogram       | kg                 | mass                                       |
| second         | s                  | time                                       |
| ampere         | A                  | electric current                           |
| kelvin         | K                  | thermodynamic temperature                  |
| mole           | mol                | amount of substance                        |
| candela        | cd                 | luminous intensity                         |
| radian         | rad                | angle                                      |
| steradian      | sr                 | solid angle                                |
| hertz          | Hz                 | frequency                                  |
| newton         | N                  | force, weight                              |
| pascal         | Pa                 | pressure, stress                           |
| joule          | J                  | energy, work, heat                         |
| watt           | W                  | power, radiant flux                        |
| coulomb        | C                  | electric charge or quantity of electricity |
| volt           | V                  | voltage (electrical potential), emf        |
| farad          | F                  | capacitance                                |
| ohm            | $\Omega$           | resistance, impedance, reactance           |
| siemens        | S                  | electrical conductance                     |
| weber          | Wb                 | magnetic flux                              |
| tesla          | T                  | magnetic flux density                      |
| henry          | H                  | inductance                                 |
| degree Celsius | $^{\circ}\text{C}$ | temperature relative to 273.15 K           |
| lumen          | lm                 | luminous flux                              |
| lux            | lx                 | illuminance                                |
| becquerel      | Bq                 | radioactivity (decays per unit time)       |
| gray           | Gy                 | absorbed dose (of ionizing radiation)      |
| sievert        | Sv                 | equivalent dose (of ionizing radiation)    |
| katal          | kat                | catalytic activity                         |

## Prefixes

### Multiples

| Prefix name | Prefix symbol | Factor          |
|-------------|---------------|-----------------|
| deca        | da            | 10 <sup>1</sup> |
| hecto       | h             | 10 <sup>2</sup> |
| kilo        | k             | 10 <sup>3</sup> |
| mega        | M             | 10 <sup>6</sup> |

| Prefix name           | Prefix symbol | Factor           |
|-----------------------|---------------|------------------|
| <a href="#">giga</a>  | G             | 10 <sup>9</sup>  |
| <a href="#">tera</a>  | T             | 10 <sup>12</sup> |
| <a href="#">peta</a>  | P             | 10 <sup>15</sup> |
| <a href="#">exa</a>   | E             | 10 <sup>18</sup> |
| <a href="#">zetta</a> | Z             | 10 <sup>21</sup> |
| <a href="#">yotta</a> | Y             | 10 <sup>24</sup> |

## Submultiples

| Prefix name           | Prefix symbol | Factor            |
|-----------------------|---------------|-------------------|
| <a href="#">deci</a>  | d             | 10 <sup>-1</sup>  |
| <a href="#">centi</a> | c             | 10 <sup>-2</sup>  |
| <a href="#">milli</a> | m             | 10 <sup>-3</sup>  |
| <a href="#">micro</a> | μ             | 10 <sup>-6</sup>  |
| <a href="#">nano</a>  | n             | 10 <sup>-9</sup>  |
| <a href="#">pico</a>  | p             | 10 <sup>-12</sup> |
| <a href="#">femto</a> | f             | 10 <sup>-15</sup> |
| <a href="#">atto</a>  | a             | 10 <sup>-18</sup> |
| <a href="#">zepto</a> | z             | 10 <sup>-21</sup> |
| <a href="#">yocto</a> | y             | 10 <sup>-24</sup> |

## Chapter 19

# Appendix VI: MEG file formats

Each MEG system brand has specific file organization and data formats. RECOMMENDED values for `manufacturer_specific_extensions`:

| Value             | Definition   |
|-------------------|--|
| <code>ctf</code>  | CTF (folder with <code>.ds</code> extension)   |
| <code>fif</code>  | Neuromag / Elekta / MEGIN and BabyMEG (file with extension <code>.fif</code> )   |
| <code>4d</code>   | BTi / 4D Neuroimaging (folder containing multiple files without extensions)  |
| <code>kit</code>  | KIT / Yokogawa / Ricoh (file with extension <code>.sqd</code> , <code>.con</code> , <code>.raw</code> , <code>.ave</code> or <code>.mrk</code> ) |
| <code>kdf</code>  | KRISS (file with extension <code>.kdf</code> )   |
| <code>itab</code> | Chieti system (file with extension <code>.raw</code> and <code>.mhd</code> )   |

Below are specifications for each system brand.

### CTF

Each experimental run with a CTF system yields a folder with a `.ds` extension, containing several files. The OPTIONAL digitized positions of the head points are usually stored in a separate `.pos` file, not necessarily within the `.ds` folder.

```
1 [sub-<label>[_ses-<label>]_headshape.pos]
2 sub-<label>[_ses-<label>]_task-<label>[_run-<index>]_meg.ds>
```

CTF's data storage is therefore via directories containing multiple files. The files contained within a `.ds` directory are named such that they match the parent directory, but preserve the original file extension (e.g., `.meg4`, `.res4`, etc.). The renaming of CTF datasets SHOULD be done with a specialized software such as the CTF newDs command-line application or [MNE-BIDS](#).

Example:

```
1 sub-control01/
2   ses-001/
3     sub-control01_ses-001_scans.tsv
4     meg/
5       sub-control01_ses-001_coordsystem.json
6       sub-control01_ses-001_headshape.pos
7       sub-control01_ses-001_task-rest_run-01_meg.ds
8       sub-control01_ses-001_task-rest_run-01_meg.json
9       sub-control01_ses-001_task-rest_run-01_channels.tsv
```

To learn more about CTF's data organization: [http://www.fieldtriptoolbox.org/getting\\_started/ctf](http://www.fieldtriptoolbox.org/getting_started/ctf)

## Neuromag/Elekta/MEGIN

Neuromag/Elekta/MEGIN data and Tristan Technologies BabyMEG data is stored with file extension `.fif`. The digitized positions of the head points are saved inside the fif file along with the MEG data, with typically no `_headshape` file.

```
1 sub-<label>[_ses-<label>]_task-<label>[_run-<index>]_meg.fif
```

Note that we do not provide specifications for cross-talk and fine-calibration matrix files in the current version of BIDS.

Example:

```
1 sub-control01/
2   ses-001/
3     sub-control01_ses-001_scans.tsv
4     meg/
5       sub-control01_ses-001_coordsystem.json
6       sub-control01_ses-001_task-rest_run-01_meg.fif
7       sub-control01_ses-001_task-rest_run-01_meg.json
8       sub-control01_ses-001_task-rest_run-01_channels.tsv
```

After applying the MaxFilter pre-processing tool, files should be renamed with the corresponding label (e.g., `proc-sss`) and placed into a `derivatives` subfolder.

Example:

```
1 sub-control01_ses-001_task-rest_run-01_proc-sss_meg.fif
2 sub-control01_ses-001_task-rest_run-01_proc-sss_meg.json
```

In the case of data runs exceeding 2Gb, the data is stored in two separate files:

```
1 sub-control01_ses-001_task-rest_run-01_meg.fif
2 sub-control01_ses-001_task-rest_run-01_meg-1.fif
```

Each of these two files has a pointer to the next file. In some software applications, like [MNE](https://mne.tools/), one can simply specify the name of the first file, and data will be read in both files via this pointer. For this reason, it is RECOMMENDED to rename and write back the file once read, to avoid the persistence of a pointer associated with the old file name.

Naming convention:

```
1 sub-control01_ses-001_task-rest_run-01_part-01_meg.fif
2 sub-control01_ses-001_task-rest_run-01_part-02_meg.fif
```

More about the Neuromag/Elekta/MEGIN data organization at: [http://www.fieldtriptoolbox.org/getting\\_started/neuromag](http://www.fieldtriptoolbox.org/getting_started/neuromag) And BabyMEG : [http://www.fieldtriptoolbox.org/getting\\_started/babysquid](http://www.fieldtriptoolbox.org/getting_started/babysquid)

## BTi/4D neuroimaging

Each experimental run on a 4D neuroimaging/BTi system results in a folder containing multiple files without extensions.

```
1 [sub-<label>[_ses-<label>]_headshape.pos]
2 sub-<label>[_ses-<label>]_task-<label>[_run-<index>]_meg>
```

One SHOULD rename/create a father run specific directory and keep the original files for each run inside (e.g., `c,rfhp0.1Hz`, `config` and `hs_file`).

Example:

```
1 sub-control01/
2   ses-001/
3     sub-control01_ses-001_scans.tsv
4     meg/
5       sub-control01_ses-001_coordsystem.json
```

```

6      sub-control01_ses-001_headshape.pos
7      sub-control01_ses-001_task-rest_run-01_meg
8      sub-control01_ses-001_task-rest_run-01_meg.json
9      sub-control01_ses-001_task-rest_run-01_channels.tsv

```

Where:

```

1 sub-control01_ses-001_task-rest_run-01_meg/
2     config
3     hs_file
4     e,rfhp1.0Hz.COH
5     c,rfDC

```

More about the 4D neuroimaging/BTi data organization at: [http://www.fieldtriptoolbox.org/getting\\_started/bti](http://www.fieldtriptoolbox.org/getting_started/bti)

## KIT/Yokogawa/Ricoh

Each experimental run on a KIT/Yokogawa/Ricoh system yields a raw (.sqd, .con) file with its associated marker coil file(s) (.sqd, .mrk), which contains coil positions in the acquisition system's native space. Head points and marker points in head space are acquired using third-party hardware.

Example:

```

1 sub-control01/
2     ses-001/
3         sub-control01_ses-001_scans.tsv
4         meg/
5             sub-control01_ses-001_coordsystem.json
6             sub-control01_ses-001_headshape.txt
7             sub-control01_ses-001_task-rest_run-01_meg
8             sub-control01_ses-001_task-rest_run-01_meg.json
9             sub-control01_ses-001_task-rest_run-01_channels.tsv
10            sub-control01_ses-001_task-rest[_acq-<label>]_run-01_markers.<mrk,sqd>
11            sub-control01_ses-001_task-rest_run-01_meg.<con,sqd>

```

If there are files with multiple marker coils, the marker files must have the `acq-<label>` parameter and no more than two marker files may be associated with one raw data file. While the acquisition parameter can take any value, it is RECOMMENDED that if the two marker measurements occur before and after the raw data acquisition, `pre` and `post` are used to differentiate the two situations.

More about the KIT/Yokogawa/Ricoh data organization at: [http://www.fieldtriptoolbox.org/getting\\_started/yokogawa](http://www.fieldtriptoolbox.org/getting_started/yokogawa)

## KRISS

Each experimental run on the KRISS system produces a file with extension .kdf. Additional files can be available in the same folder: the digitized positions of the head points (\\_digitizer.txt), the position of the center of the MEG coils (.chn) and the event markers (.trg).

```

1 [sub-<label>[_ses-<label>]_headshape.txt]
2 sub-<label>[_ses-<label>]_task-<label>[_run-<index>]_meg.kdf
3 sub-<label>[_ses-<label>]_task-<label>[_run-<index>]_meg.chn
4 sub-<label>[_ses-<label>]_task-<label>[_run-<index>]_meg.trg
5 sub-<label>[_ses-<label>]_task-<label>[_acq-<label>]_digitizer.txt

```

Example:

```

1 sub-control01/
2     ses-001/
3         sub-control01_ses-001_scans.tsv
4         meg/
5             sub-control01_ses-001_coordsystem.json

```



```
6      sub-control01_ses-001_headshape.txt
7      sub-control01_ses-001_task-rest_run-01_meg
8      sub-control01_ses-001_task-rest_run-01_meg.json
9      sub-control01_ses-001_task-rest_run-01_channels.tsv
10     sub-control01_ses-001_task-rest_run-01_meg.chn
11     sub-control01_ses-001_task-rest_run-01_meg.kdf
12     sub-control01_ses-001_task-rest_run-01_meg.trg
13     sub-control01_ses-001_task-rest_digitizer.txt
```

## ITAB

Each experimental run on a ITAB-ARGOS153 system yields a raw (.raw) data file plus an associated binary header file (.mhd). The raw data file has an ASCII header that contains detailed information about the data acquisition system, followed by binary data. The associated binary header file contains part of the information from the ASCII header, specifically the one needed to process data, plus other information on offline preprocessing performed after data acquisition (e.g., sensor position relative to subject's head, head markers, stimulus information).

Example:

```
1 sub-control01/
2   ses-001/
3     sub-control01_ses-001_coordsystem.json
4     sub-control01_ses-001_headshape.txt
5     sub-control01_ses-001_task-rest_run-01_meg
6     sub-control01_ses-001_task-rest_run-01_meg.json
7     sub-control01_ses-001_task-rest_run-01_channels.tsv
8     sub-control01_ses-001_task-rest_run-01_meg.raw
9     sub-control01_ses-001_task-rest_run-01_meg.raw.mhd
```

## Aalto MEG–MRI

For stand-alone MEG data, the Aalto hybrid device uses the standard .fif data format and follows the conventions of Elekta/Neuromag as described above. The .fif files may contain unreconstructed MRI data. The inclusion of MRI data and information for accurate reconstruction will be fully standardized at a later stage.

## Chapter 20

# Appendix VII: preferred names of MEG systems

Restricted keywords for Manufacturer field in the \*\_meg.json file:

- CTF
- Neuromag/Elekta/Megin
- BTi/4D
- KIT/Yokogawa/Ricoh
- KRISS
- ITAB
- [Aalto/MEG-MRI](#)
- Other

Restricted keywords for ManufacturersModelName field in the \*\_meg.json file:

| System Model Name         | Manufacturer          | Details   |
|---------------------------|-----------------------|---|
| CTF-64                    | CTF                   |   |
| CTF-151                   | CTF                   | <a href="https://www.ctf.com/products">https://www.ctf.com/products</a>                                 |
| CTF-275                   | CTF                   | CTF-275: OMEGA 2000   |
| Neuromag-122              | Neuromag/Elekta/Megin |   |
| ElektaVectorview          | Neuromag/Elekta/Megin | 102 magnetometers + 204 planar gradiometers   |
| ElektaTRIUX               | Neuromag/Elekta/Megin | <a href="https://www.elekta.com/diagnostic-solutions/">https://www.elekta.com/diagnostic-solutions/</a> |
| 4D-Magnes-WH2500          | BTi/4D                |   |
| 4D-Magnes-WH3600          | BTi/4D                |   |
| KIT-157                   | KIT/Yokogawa          |   |
| KIT-160                   | KIT/Yokogawa          |   |
| KIT-208                   | KIT/Yokogawa          |   |
| ITAB-ARGOS153             | ITAB                  |   |
| Aalto-MEG-MRI-<br>YYYY/MM | Aalto/MEG-MRI         | YYYY-MM (year, month; or major version)   |

## Chapter 21

# Appendix VIII: Preferred Names of Coordinate Systems

To interpret a coordinate (x, y, z), it is required that you know relative to which origin the coordinates are expressed, you have to know the interpretation of the three axes, and you have to know the units in which the numbers are expressed. This information is sometimes called the coordinate system.

These letters help describe the coordinate system definition:

- A/P means anterior/posterior
- L/R means left/right
- S/I means superior/inferior

For example: **RAS** means that the first dimension (X) points towards the right hand side of the head, the second dimension (Y) points towards the Anterior aspect of the head, and the third dimension (Z) points towards the top of the head.

Besides coordinate systems, defined by their origin and direction of the axes, BIDS defines "spaces" as an artificial frame of reference, created to describe different anatomies in a unifying manner (see e.g., <https://doi.org/10.1016/j.neuroimage.2012.01.024>). The "space" and all coordinates expressed in this space are by design a transformation of the real world geometry, and nearly always different from the individual subject space that it stems from. An example is the Talairach-Tournoux space, which is constructed by piecewise linear scaling of an individual's brain to that of the Talairach-Tournoux 1988 atlas. In the Talairach-Tournoux space, the origin of the coordinate system is at the AC and units are expressed in mm.

The coordinate systems below all relate to neuroscience and therefore to the head or brain coordinates. Please be aware that all data acquisition starts with "device coordinates" (scanner), which does not have to be identical to the initial "file format coordinates" (DICOM), which are again different from the "head" coordinates (e.g., NIFTI). Not only do device coordinate vary between hardware manufacturers, but also the head coordinates differ, mostly due to different conventions used in specific software packages developed by different (commercial or academic) groups.

## Coordinate Systems applicable to MEG, EEG, and iEEG

Generally, across the MEG, EEG, and iEEG modalities, the first two pieces of information (origin, orientation) are specified in `XXXCoordinateSystem`, and the units are specified in `XXXCoordinateSystemUnits`.

Allowed values for the `XXXCoordinateSystem` field come from a list of restricted keywords, as listed in the sections below. If no value from the list of restricted keywords fits, there is always the option to specify the value as follows:

- **Other:** Use this for other coordinate systems and specify further details in the `XXXCoordinateSystemDescription` field

## MEG Specific Coordinate Systems

Restricted keywords for the `XXXCoordinateSystem` field in the `coordinatesystem.json` file for MEG datasets:

- **CTF**: ALS orientation and the origin between the ears
- **ElektaNeuromag**: RAS orientation and the origin between the ears
- **4DBti**: ALS orientation and the origin between the ears
- **KitYokogawa**: ALS orientation and the origin between the ears
- **ChietiItab**: RAS orientation and the origin between the ears

Note that the short descriptions above do not capture all details, There are detailed extensive descriptions of these EEG coordinate systems on the [FieldTrip toolbox web page](#)

## EEG Specific Coordinate Systems

Restricted keywords for the `XXXCoordinateSystem` field in the `coordsystem.json` file for EEG datasets:

- **BESA**: Although natively this is a spherical coordinate system, the electrode positions should be expressed in Cartesian coordinates, with a RAS orientation. The X axis is the T8-T7 line, positive at T8. The Y axis is the Oz-Fpz line, positive at Fpz. The origin of the sphere fitted to the electrodes is approximately 4 cm above the point between the ears.
- **Captrak**: RAS orientation and the origin between the ears

Note that the short descriptions above do not capture all details, There are detailed extensive descriptions of these EEG coordinate systems on the [FieldTrip toolbox web page](#) and on the [BESA wiki](#).

## iEEG Specific Coordinate Systems

Restricted keywords for the `XXXCoordinateSystem` field in the `coordsystem.json` file for iEEG datasets:

- **Pixels**: If electrodes are localized in 2D space (only x and y are specified and z is n/a), then the positions in this file must correspond to the locations expressed in pixels on the photo/drawing/rendering of the electrodes on the brain. In this case, coordinates must be (row,column) pairs, with (0,0) corresponding to the upper left pixel and (N,0) corresponding to the lower left pixel.
- **ACPC**: The origin of the coordinate system is at the Anterior Commissure and the negative y-axis is passing through the Posterior Commissure. The positive z-axis is passing through a mid-hemispheric point in the superior direction. The anatomical landmarks are determined in the individual's anatomical scan and no scaling or deformations have been applied to the individual's anatomical scan. For more information, see the [ACPC site](#) on the FieldTrip toolbox wiki.

## Template Based Coordinate Systems

The transformation of the real world geometry to an artificial frame of reference is described in `XXXCoordinateSystem`. Unless otherwise specified below, the origin is at the AC and the orientation of the axes is RAS. Unless specified explicitly in the sidecar file in the `XXCoordinateSystemUnits` field, the units are assumed to be mm.

### Standard template identifiers

| Coordinate System | Description   |
|-------------------|---|
| ICBM452AirSpace   | Reference space defined by the "average of 452 T1-weighted MRIs of normal young adult brains" with "linear transforms of the subjects into the atlas space using a 12-parameter affine transformation"<br><a href="https://www.loni.usc.edu/research/atlas">https://www.loni.usc.edu/research/atlas</a> |
| ICBM452Warp5Space | Reference space defined by the "average of 452 T1-weighted MRIs of normal young adult brains" "based on a 5th order polynomial transformation into the atlas space" <a href="https://www.loni.usc.edu/research/atlas">https://www.loni.usc.edu/research/atlas</a>                                       |
| IXI549Space       | Reference space defined by the average of the "549 (...) subjects from the IXI dataset" linearly transformed to ICBM MNI 452. Used by SPM12.<br><a href="http://www.brain-development.org/">http://www.brain-development.org/</a>   |

| Coordinate System             | Description  |
|-------------------------------|--|
| fsaverage[ Sym]               | The <b>fsaverage</b> is a <b>dual template</b> providing both volumetric and surface coordinates references. The volumetric template corresponds to a FreeSurfer variant of MNI305 space. The <b>fsaverage</b> atlas also defines a surface reference system (formerly described as fsaverage[3 4 5 6 sym]). The <b>fsaverageSym</b> is the symmetric counterpart of <b>fsaverage</b> .  |
| fsLR                          | The <b>fsLR</b> is a <b>dual template</b> providing both volumetric and surface coordinates references. The volumetric template corresponds to MNI152NLin6Asym. Surface templates are given at several sampling densities: 164k (used by HCP pipelines for 3T and 7T anatomical analysis), 59k (used by HCP pipelines for 7T MRI bold and DWI analysis), 32k (used by HCP pipelines for 3T MRI bold and DWI analysis), or 4k (used by HCP pipelines for MEG analysis) fsaverage_LR surface reconstructed from the T1w image. |
| MNIColin27                    | Average of 27 T1 scans of a single subject<br><a href="http://www.bic.mni.mcgill.ca/ServicesAtlases/Colin27Highres">http://www.bic.mni.mcgill.ca/ServicesAtlases/Colin27Highres</a>  |
| MNI152Lin                     | Also known as ICBM (version with linear coregistration)<br><a href="http://www.bic.mni.mcgill.ca/ServicesAtlases/ICBM152Lin">http://www.bic.mni.mcgill.ca/ServicesAtlases/ICBM152Lin</a>   |
| MNI152NLin2009[a-c][Sym Asym] | Also known as ICBM (non-linear coregistration with 40 iterations, released in 2009). It comes in either three different flavours each in symmetric or asymmetric version.<br><a href="http://www.bic.mni.mcgill.ca/ServicesAtlases/ICBM152NLin2009">http://www.bic.mni.mcgill.ca/ServicesAtlases/ICBM152NLin2009</a>   |
| MNI152NLin6Sym                | Also known as symmetric ICBM 6th generation (non-linear coregistration). Used by SPM99 - SPM8.<br><a href="http://www.bic.mni.mcgill.ca/ServicesAtlases/ICBM152NLin6">http://www.bic.mni.mcgill.ca/ServicesAtlases/ICBM152NLin6</a>  |
| MNI152NLin6Asym               | A variation of MNI152NLin6Sym built by A. Janke that is released as the <i>MNI template</i> of FSL. Volumetric templates included with <a href="#">HCP-Pipelines</a> correspond to this template too. See <a href="https://doi.org/10.1016/j.neuroimage.2012.01.024">10.1016/j.neuroimage.2012.01.024</a> .  |
| MNI305                        | Also known as avg305.  |
| NIHPD[Sym Asym]               | Pediatric templates generated from the NIHPD sample. Available for different age groups (4.5–18.5 y.o., 4.5–8.5 y.o., 7–11 y.o., 7.5–13.5 y.o., 10–14 y.o., 13–18.5 y.o. This template also comes in either -symmetric or -asymmetric flavor. <a href="http://www.bic.mni.mcgill.ca/ServicesAtlases/NIHPD-obj1">http://www.bic.mni.mcgill.ca/ServicesAtlases/NIHPD-obj1</a>  |
| OASIS30AntsOASISAnts          | <a href="https://figshare.com/articles/ANTs_ANTsR_Brain_Templates/915436">https://figshare.com/articles/ANTs_ANTsR_Brain_Templates/915436</a>  |
| OASIS30Atropos                | <a href="https://mindboggle.info/data.html">https://mindboggle.info/data.html</a>  |
| Talairach                     | Piecewise linear scaling of the brain is implemented as described in TT88.<br><a href="http://www.talairach.org/">http://www.talairach.org/</a>  |
| UNCInfant                     | Infant Brain Atlases from Neonates to 1- and 2-year-olds.<br><a href="https://www.nitrc.org/projects/pediatricatlas">https://www.nitrc.org/projects/pediatricatlas</a>   |

The following template identifiers are retained for backwards compatibility of BIDS implementations. Their use is NOT RECOMMENDED for new BIDS datasets and tooling, but their presence MUST NOT produce a validation error.

| Coordinate System           | Recommended identifier |
|-----------------------------|------------------------|
| fsaverage[3 4 5 6 sym]      | fsaverage[ Sym]        |
| UNCInfant[0 1 2 V 21 22 23] | UNCInfant              |

# Chapter 22

## Changelog

### v1.2.2 (2020-02-12)

- BEP 18 suggestions [#398](#) (effigies)
- [MISC] Add Sébastien Tourbier to contributors [#394](#) (sebastientourbier)
- [MISC] moved list of extension proposals to the main BIDS website [#389](#) (robertoostenveld)
- [FIX] Typos and clarifications [#386](#) (apjanke)
- ENH: Add watermark to drafts [#383](#) (effigies)
- [MISC] Teon Brooks retiring moderator duties for BEP021 [#381](#) (teonbrooks)
- [FIX] clarify that string is expected for HowToAcknowledge field in dataset\_description.json [#380](#) (sappelhoff)
- [MISC] Typo and style [#378](#) (TheChymera)
- [FIX] divide readme into 3 parts [#374](#) (sappelhoff)
- FIX: Refer to BIDS consistently, instead of "<Modality>-BIDS" [#366](#) (sappelhoff)
- [FIX] Change recommended anonymization date from 1900 to 1925 [#363](#) (alexrockhill)
- [FIX] Minor fixups of inconsistencies while going through a PDF version [#362](#) (yarikoptic)
- specify further the pipeline following [#345](#) [#358](#) (CPernet)
- FIX: clarify that filters should be specified as object of objects [#348](#) (sappelhoff)
- FIX: Clarify channels.tsv is RECOMMENDED consistently across ephys [#347](#) (sappelhoff)
- [FIX] Typo fix (contract -> contrast) in events documentation [#346](#) (snastase)
- RM: rm TOC.md - seems no longer pertinent/used [#341](#) (yarikoptic)
- [MISC] Move the PR template to a separate folder and improve contents [#338](#) (jhlegarreta)
- [INFRA] Find npm requirements file in Circle [#336](#) (franklin-feingold)
- ENH: Clarify the position toward non-compliant derivative datasets and files [#334](#) (effigies)
- [ENH] Clarify phenotypic and assessment data in new section [#331](#) (sappelhoff)
- [MISC] add information about continuous integration checks to PR template [#330](#) (sappelhoff)
- [FIX] Fix Common principles Key/value files section level [#328](#) (jhlegarreta)
- [INFRA] Set the maximum heading length lint check to false [#325](#) (jhlegarreta)
- [FIX] Number explicitly all cases in MRI field map section headers [#323](#) (jhlegarreta)
- [FIX] Add SoftwareFilters to EEG sidecar example [#322](#) (Remi-Gau)
- [MISC] Fixing Travis errors with Remark [#320](#) (franklin-feingold)
- [INFRA] Link to doc builds in CI checks [#315](#) (jasmainak)
- [MISC] Add BEP027 - BIDS Execution to BEP list [#314](#) (effigies)
- [FIX] Add CBV and phase to Entity table [#312](#) (tsalo)
- [FIX] Separate out imaging-specific "common derivatives" [#310](#) (effigies)
- [FIX] Revise template-generated coordinate systems [#307](#) (oesteban)
- [FIX] Normalization of template-generated standard spaces [#306](#) (oesteban)
- [ENH] Release protocol notes [#304](#) (franklin-feingold)
- [INFRA] Adding contributor appendix sentence to PR template [#299](#) (franklin-feingold)
- [ENH] Added discontinuous datatype for EEG and iEEG [#286](#) (wouterpotters)
- [FIX] Clarify paragraph about custom data types [#264](#) (effigies)

### v1.2.1 (2019-08-14)

- FIX: repair link in anatomical MRI table [#297](#) (sappelhoff)

- [ENH] Clarify requirements in Release Protocol [#294](#) ([franklin-feingold](#))
- [INFRA+FIX] Use linkchecker (from a dedicated docker image) to check all URLs [#293](#) ([yarikoptic](#))
- [ENH] Adding Contributors and updating contributions [#284](#) ([franklin-feingold](#))
- [MISC] update Code of Conduct contact [#281](#) ([franklin-feingold](#))
- [ENH] Update contributing guide and README to make discussion forums easy to find [#279](#) ([emdupre](#))
- [ENH] Starter Kit dropdown menu [#278](#) ([franklin-feingold](#))
- [ENH] BEP Update [#277](#) ([franklin-feingold](#))
- [INFRA] Update pipenv [#274](#) ([sappelhoff](#))
- [INFRA] Transpose the entity table and link to text anchors describing each entity [#272](#) ([sappelhoff](#))
- [ENH] Add Twitter badge to README and link to website to landing page [#268](#) ([franklin-feingold](#))
- [ENH] adding release guidelines [#267](#) ([franklin-feingold](#))
- [FIX] Common principles: Fix filename in inheritance principle [#261](#) ([Lestropie](#))
- [MISC] update modality references [#258](#) ([sappelhoff](#))
- [INFRA] adding logo to RTD [#256](#) ([franklin-feingold](#))
- [INFRA] add footer, replacing mkdocs/material advert with Github link [#250](#) ([sappelhoff](#))
- [MISC] rename logo files, add a README of where they come from, fix favicon [#249](#) ([sappelhoff](#))
- [MISC] updating MEG doc links, manufacturer names, and adding a missing MEG example [#248](#) ([sappelhoff](#))
- [ENH] Add favicon to RTD [#246](#) ([franklin-feingold](#))
- [MISC] Update Authors in BEP025 [#241](#) ([josator2](#))
- [MISC] Document BEPs that are not active anymore, but have not been merged [#240](#) ([sappelhoff](#))
- [FIX] remove ManufacturersAmplifierModelName (again) [#236](#) ([robertoostenveld](#))
- [INFRA] Update release protocol [#235](#) ([effigies](#))
- [INFRA] Enable version panel for quickly finding previous versions [#232](#) ([effigies](#))
- [FIX] use \<label> for \_desc- - not some \<value> + clarify \<value> [#224](#) ([yarikoptic](#))
- [FIX] Clarify Appendix II: The list of licenses only lists examples [#222](#) ([sappelhoff](#))
- [FIX] Trivial column header fix [#220](#) ([nicholst](#))
- [INFRA] Add clarification on merge methods to DECISION\_MAKING [#217](#) ([sappelhoff](#))
- [INFRA] Enable permalink urls to appear at (sub)section headings [#214](#) ([yarikoptic](#))
- [INFRA] bump up mkdocs-materials version [#211](#) ([sappelhoff](#))
- [ENH] Various proposed changes to diffusion derivatives [#205](#) ([Lestropie](#))
- [MISC] Fix github username for @chrisgorgo [#204](#) ([chrisgorgo](#))
- [FIX] clarify example 3 in common principles (inheritance) [#202](#) ([sappelhoff](#))
- [MISC] Expand entity table for MEG/EEG/iEEG specific files [#198](#) ([sappelhoff](#))
- [FIX] make iEEG ToC more consistent with MEG and EEG [#191](#) ([robertoostenveld](#))
- [FIX] Clarify use of acq and task parameters in EEG, MEG, and iEEG [#188](#) ([sappelhoff](#))
- [FIX] clarify use of tools for CTF data renaming [#187](#) ([sappelhoff](#))
- [MISC] Add bep006 and bep010 to completed beps and fix links [#186](#) ([sappelhoff](#))
- [FIX] change file for definition of electrical stimulation labels from \_electrodes.json to \_events.json [#185](#) ([ezemikulan](#))
- [ENH] relax ieeg channel name requirements of letters and numbers only [#182](#) ([sappelhoff](#))
- [FIX] make MEG section headings and ToC consistent to the EEG and iEEG specs [#181](#) ([robertoostenveld](#))
- [FIX] make section headings and ToC consistent between meg and eeg specs [#180](#) ([robertoostenveld](#))
- [MISC] Spelling fixes [#179](#) ([DimitriPapadopoulos](#))
- [ENH] Alternative folder organization for raw, derived, and source data [#178](#) ([chrisgorgo](#))
- [INFRA] Adding instructions for naming PRs [#177](#) ([chrisgorgo](#))
- [MISC] Introducing Stefan Appelhoff as the first Maintainer [#176](#) ([chrisgorgo](#))
- [FIX] Clarify name of "BrainVision" format [#175](#) ([JegouA](#))
- [FIX] Fixes spelling of continuous [#171](#) ([emdupre](#))
- [FIX] Clarify continuous recording metadata fields [#167](#) ([effigies](#))
- [FIX] changed reference of dcm2nii to dcm2niix [#166](#) ([DimitriPapadopoulos](#))
- [FIX] Removing a leftover file [#162](#) ([chrisgorgo](#))
- [ENH] Derived (processed) MR data [#109](#) ([chrisgorgo](#))
- [FIX] Specify marker file names for KIT data (MEG) [#62](#) ([monkeyman192](#))
- [FIX] Remove father-level for meg filetypes other than BTi/4D data [#19](#) ([teonbrooks](#))

## v1.2.0 (2019-03-04)

- [MISC] Adding Dimitri Papadopoulos Orfanos to the list of contributors [#157](#) ([DimitriPapadopoulos](#))



- RF: use "specification" not "protocol" to refer to BIDS [#156](#) (yarikoptic)
- Fix example misalignment [#155](#) (DimitriPapadopoulos)
- Update Pipfile.lock [#144](#) (franklin-feingold)
- [MRG] clarify decimal sep and numerical notation convention [#143](#) (sappelhoff)
- [MRG] clarify encoding of README, CHANGES, TSV, and JSON files [#140](#) (sappelhoff)
- Update site\_name and release protocol [#137](#) (franklin-feingold)
- BF: Example for IntendedFor was missing session indicator in the filename [#129](#) (yarikoptic)
- [ENH] Add "\_phase" suffix to func datatype for functional phase data [#128](#) (tsalo)
- Update to Release\_Protocol.md [#126](#) (franklin-feingold)
- Update tag naming convention [#123](#) (chrisgorgo)
- [MRG] Merge bep006 and bep010 [#108](#) (sappelhoff)
- Adding formal decision-making rules [#104](#) (chrisgorgo)
- [MRG] number of small corrections to the specification [#98](#) (robertoostenveld)

## v1.1.2 (2019-01-10)

- Update 01-contributors.md [#120](#) (oesteban)
- Global fields in data dictionaries [#117](#) (chrisgorgo)
- Propose BEP026 MER [#116](#) (greydongilmore)
- Remove duplicate entries in MEG table [#113](#) (franklin-feingold)
- Propose BEP025 MIDS [#110](#) (josator2)
- repair links [#106](#) (sappelhoff)
- Autogenerate CHANGES.md [#103](#) (franklin-feingold)
- Added contributor information [#100](#) (jgrethe)
- ENH: First(?) good practice recommendation. No excessive overrides in Inheritance principle [#99](#) (yarikoptic)
- adding extensions page [#97](#) (choldgraf)
- Fix up some urls (as detected to be broken/inconsistent [#95](#) (yarikoptic)
- Change BEP numbers to include MRS [#94](#) (Hboni)
- RF harmonize and thus shorten templates etc [#93](#) (yarikoptic)
- put links and some text into README [#91](#) (sappelhoff)
- Add extension proposal in 01-introduction.md [#88](#) (Hboni)
- additional table to recap 'volume acquisition timing' [#87](#) (Remi-Gau)
- Small typo in "scanning sequence" DICOM tag [#84](#) (Remi-Gau)
- Update 01-contributors.md [#83](#) (teonbrooks)
- Added CBV contrast [#82](#) (TheChymera)
- Add CC-BY 4.0 license [#81](#) (KirstieJane)
- Fix Travis break [#80](#) (franklin-feingold)
- ENH: allow \_dir for other EPI (func, dwi) sequences [#78](#) (yarikoptic)
- Added appendix to mkdocs and added some internal links [#77](#) (franklin-feingold)
- DOC: added JC Houde as contributor. [#76](#) (jchoude)
- Updated my contributions [#75](#) (nicholst)
- update HED appendix [#74](#) (sappelhoff)
- [MRG] unicode: replace greek mu and omega by micro and ohm signs [#73](#) (sappelhoff)
- Update 01-contributors.md [#72](#) (francopestilli)
- add ce-<label> for fmri data [#70](#) (dasturge)
- pin pip version [#68](#) (chrisgorgo)
- Fix link in index [#46](#) (chrisgorgo)
- edit contributing md [#44](#) (Park-Patrick)
- Mkdocs configuration and RTD setup [#42](#) (choldgraf)
- Move definitions, compulsory, and raw/derivatives sections to principles [#40](#) (chrisgorgo)
- Remove duplicate section [#39](#) (chrisgorgo)
- mkdocs rendering [#36](#) (chrisgorgo)
- Style consistency [#35](#) (chrisgorgo)
- Renaming files to conform with style guide [#34](#) (chrisgorgo)
- enable travis cache [#32](#) (chrisgorgo)
- FIX - corrected link that is shown for CC0 [#31](#) (robertoostenveld)
- [WIP] added linter integration via travis [#30](#) (chrisgorgo)
- Cleanup [#29](#) (chrisgorgo)
- [MRG] split intro, commons, mr, and meg into folder from specification.md [#28](#) (teonbrooks)



- Add some bids starter kit contributors [#27](#) ([KirstieJane](#))
- Embedded footnotes into text [#25](#) ([franklin-feingold](#))
- Making HED Strategy Guide link prettier [#24](#) ([fake-filo](#))
- Fix/more cleanup [#21](#) ([chrisgorgo](#))
- formatted MEG (8.4) [#17](#) ([franklin-feingold](#))
- small fixes [#16](#) ([chrisgorgo](#))
- Add meg img [#14](#) ([sappelhoff](#))
- [WIP] Cleaning up the specification [#13](#) ([chrisgorgo](#))
- Adding code of conduct [#6](#) ([chrisgorgo](#))
- Renaming the main document [#1](#) ([chrisgorgo](#))

### 1.1.1

- Improved the MEG landmark coordinates description.
- Replaced ManufacturersCapModelName in meg.json with CapManufacturer and CapManufacturersModelName.
- Remove EEGSamplingFrequency and ManufacturersAmplifierModelName from the meg.json.
- Improved the behavioral data description.

### 1.1.0

- Added support for MEG data (merged BEP008)
- Added SequenceName field.
- Added support for describing events with Hierarchical Event Descriptors [\[4.3 Task events\]](#).
- Added VolumeTiming and AcquisitionDuration fields [\[4.1 Task \(including resting state\) imaging data\]](#).
- Added DwellTime field.

### 1.0.2

- Added support for high resolution (anatomical) T2star images [\[4.1 Anatomy imaging data\]](#).
- Added support for multiple defacing masks [\[4.1 Anatomy imaging data\]](#).
- Added optional key and metadata field for contrast enhanced structural scans [\[4.1 Anatomy imaging data\]](#)
- Added DelayTime field [\[4.1 Task \(including resting state\) imaging data\]](#).
- Added support for multi echo BOLD data [\[4.1 Task \(including resting state\) imaging data\]](#).

### 1.0.1

- Added InstitutionName field [\[4.1 Task \(including resting state\) imaging data\]](#).
- Added InstitutionAddress field [\[4.1 Task \(including resting state\) imaging data\]](#).
- Added DeviceSerialNumber field [\[4.1 Task \(including resting state\) imaging data\]](#).
- Added NumberOfVolumesDiscardedByUser and NumberOfVolumesDiscardedByScanner field [\[4.1 Task \(including resting state\) imaging data\]](#).
- Added TotalReadoutTime to functional images metadata list [\[4.1 Task \(including resting state\) imaging data\]](#).

### 1.0.1-rc1

- Added T1 Rho maps [\[4.1 Anatomy imaging data\]](#).
- Added support for phenotypic information split into multiple files [\[3.2 Participant key file\]](#).
- Added recommendations for multi site datasets
- Added SoftwareVersions

- Added run-<run\_index> to the phase encoding maps. Improved the description.
- Added InversionTime metadata key.
- Clarification on the source vs raw language.
- Added trial\_type column to the event files.
- Added missing sub-<participant\_label> in behavioral data file names
- Added ability to store stimuli files.
- Clarified the language describing allowed subject labels.
- Added quantitative proton density maps.

## 1.0.0

- Added ability to specify fieldmaps acquired with multiple parameter sets.
- Added ability to have multiple runs of the same fieldmap.
- Added FLASH anatomical images.

## 1.0.0-rc4

- Replaced links to neurolex with explicit DICOM Tags.
- Added sourcedata.
- Added data dictionaries.
- Be more explicit about contents of JSON files for structural (anatomical) scans.

## 1.0.0-rc3

- Renamed PhaseEncodingDirection values from “x”, “y”, “z” to “i”, “j”, “k” to avoid confusion with FSL parameters
- Renamed SliceEncodingDirection values from “x”, “y”, “z” to “i”, “j”, “k”

## 1.0.0-rc2

- Removed the requirement that TSV files cannot include more than two consecutive spaces.
- Refactor of the definitions sections (copied from the manuscript)
- Make support for uncompressed .nii files more explicit.
- Added BIDSVersion to dataset.json
- Remove the statement that SliceEncodingDirection is necessary for slice time correction
- Change dicom converter recommendation from dcmstack to dcm2nii and dicm2nii following interactions with the community (see <https://github.com/moloney/dcmstack/issues/39> and <https://github.com/neurolabusc/dcm2niix/issues/4>).
- Added section on behavioral experiments with no accompanying MRI acquisition
- Add \_magnitude.nii.gz image for GE type fieldmaps.
- Replaced EchoTimeDifference with EchoTime1 and EchoTime2 (SPM toolbox requires this input).
- Added support for single band reference image for DWI.
- Added DatasetDOI field in the dataset description.
- Added description of more metadata fields relevant to DWI fieldmap correction.
- PhaseEncodingDirection is now expressed in “x”, “y” etc. instead of “PA” “RL” for DWI scans (so it’s the same as BOLD scans)

- Added rec-<label> flag to BOLD files to distinguish between different reconstruction algorithms (analogous to anatomical scans).
- Added recommendation to use `_physio` suffix for continuous recordings of motion parameters obtained by the scanner side reconstruction algorithms.

## 1.0.0-rc1

- Initial release

\* *This Change Log was automatically generated by [github\\_changelog\\_generator](#)*

### 1.1.1

- Improved the MEG landmark coordinates description.
- Replaced `ManufacturersCapModelName` in `meg.json` with `CapManufacturer` and `CapManufacturersModelName`.
- Remove `EEGSamplingFrequency` and `ManufacturersAmplifierModelName` from the `meg.json`.
- Improved the behavioral data description.

### 1.1.0

- Added support for MEG data (merged BEP008)
- Added `SequenceName` field.
- Added support for describing events with Hierarchical Event Descriptors [4.3 Task events].
- Added `VolumeTiming` and `AcquisitionDuration` fields [4.1 Task (including resting state) imaging data].
- Added `DwellTime` field.

### 1.0.2

- Added support for high resolution (anatomical) T2star images [4.1 Anatomy imaging data].
- Added support for multiple defacing masks [4.1 Anatomy imaging data].
- Added optional key and metadata field for contrast enhanced structural scans [4.1 Anatomy imaging data].
- Added `DelayTime` field [4.1 Task (including resting state) imaging data].
- Added support for multi echo BOLD data [4.1 Task (including resting state) imaging data].

### 1.0.1

- Added `InstitutionName` field [4.1 Task (including resting state) imaging data].
- Added `InstitutionAddress` field [4.1 Task (including resting state) imaging data].
- Added `DeviceSerialNumber` field [4.1 Task (including resting state) imaging data].
- Added `NumberOfVolumesDiscardedByUser` and `NumberOfVolumesDiscardedByScanner` field [4.1 Task (including resting state) imaging data].
- Added `TotalReadoutTime` to functional images metadata list [4.1 Task (including resting state) imaging data].

### 1.0.1-rc1

- Added T1 Rho maps [4.1 Anatomy imaging data].
- Added support for phenotypic information split into multiple files [3.2 Participant key file].
- Added recommendations for multi site datasets
- Added `SoftwareVersions`

- Added run-<run\_index> to the phase encoding maps. Improved the description.
- Added InversionTime metadata key.
- Clarification on the source vs raw language.
- Added trial\_type column to the event files.
- Added missing sub-<participant\_label> in behavioral data file names
- Added ability to store stimuli files.
- Clarified the language describing allowed subject labels.
- Added quantitative proton density maps.

## 1.0.0

- Added ability to specify fieldmaps acquired with multiple parameter sets.
- Added ability to have multiple runs of the same fieldmap.
- Added FLASH anatomical images.

## 1.0.0-rc4

- Replaced links to neurolex with explicit DICOM Tags.
- Added sourcedata.
- Added data dictionaries.
- Be more explicit about contents of JSON files for structural (anatomical) scans.

## 1.0.0-rc3

- Renamed PhaseEncodingDirection values from “x”, “y”, “z” to “i”, “j”, “k” to avoid confusion with FSL parameters
- Renamed SliceEncodingDirection values from “x”, “y”, “z” to “i”, “j”, “k”

## 1.0.0-rc2

- Removed the requirement that TSV files cannot include more than two consecutive spaces.
- Refactor of the definitions sections (copied from the manuscript)
- Make support for uncompressed .nii files more explicit.
- Added BIDSVersion to dataset.json
- Remove the statement that SliceEncodingDirection is necessary for slice time correction
- Change dicom converter recommendation from dcmstack to dcm2nii and dicm2nii following interactions with the community (see <https://github.com/moloney/dcmstack/issues/39> and <https://github.com/neurolabusc/dcm2niix/issues/4>).
- Added section on behavioral experiments with no accompanying MRI acquisition
- Add \_magnitude.nii.gz image for GE type fieldmaps.
- Replaced EchoTimeDifference with EchoTime1 and EchoTime2 (SPM toolbox requires this input).
- Added support for single band reference image for DWI.
- Added DatasetDOI field in the dataset description.
- Added description of more metadata fields relevant to DWI fieldmap correction.
- PhaseEncodingDirection is now expressed in “x”, “y” etc. instead of “PA” “RL” for DWI scans (so it’s the same as BOLD scans)

- Added rec-<label> flag to BOLD files to distinguish between different reconstruction algorithms (analogous to anatomical scans).
- Added recommendation to use \_\_physio suffix for continuous recordings of motion parameters obtained by the scanner side reconstruction algorithms.

## **1.0.0-rc1**

- Initial release